

# HSL110 Deluxe Universal Sling



## Lifestyle Slings

**MANUFACTURED IN  
GREAT BRITAIN BY CARE-  
ABILITY HEALTHCARE.**



More information [www.harvesthealthcare.co.uk](http://www.harvesthealthcare.co.uk)

**Harvest Healthcare**  
**Sheaf House**  
**Bradmarsh Way**  
**Bradmarsh Business Park S601BW**  
**T: 01709377172**  
**F: 01709377173**  
**E: [enquiries@harvesthealthcare.co.uk](mailto:enquiries@harvesthealthcare.co.uk)**  
**W: [www.harvesthealthcare.co.uk](http://www.harvesthealthcare.co.uk)**




harvesthealthcare®

# HSL110 Deluxe Universal Sling


- 1** Ensure that the guiding handles and identification label are on the outside.




- 2** Place the sling down the back of the client using the rear pocket. Ensure the aperture is tucked under the client's coccyx.

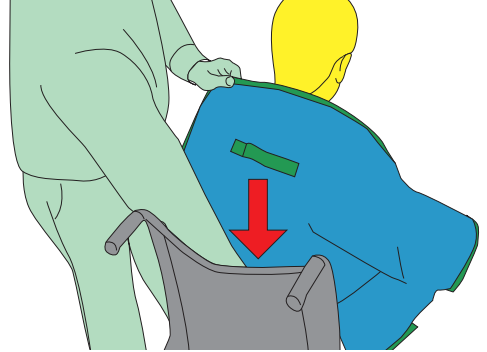

- 3** Straighten the leg section and push the side of the aperture underneath the client's hip. Repeat on the opposite side. The long outer leg straps should be pulled so they are positioned mid-thigh.

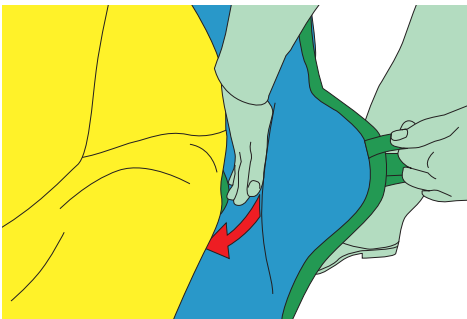

- 4** Carefully feed the leg of the sling underneath pushing the sling section down into the seating surface to avoid friction and shear against the clients legs. Repeat the procedure on the other leg.

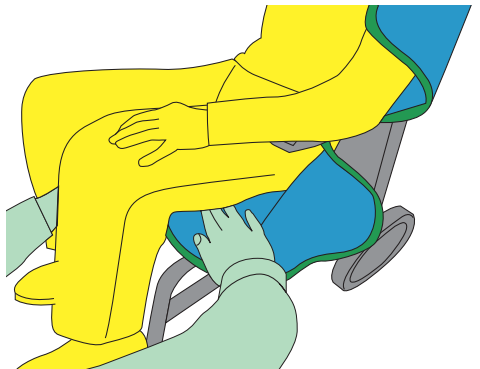


5

Ensure that the sling is not twisted or creased under the thighs and check both leg straps are equal in length. Feed the right leg through the loop on the back of the opposite leg.



6

Feed the clip on the left hand leg through the small dignity loop on the right hand leg.

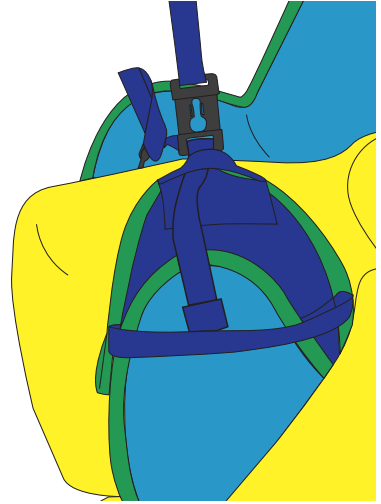


Feed the clip on the left leg through the short outer loop of the left leg and the right leg clip through the short outer loop of the right leg. Attach the leg clips to the spreader bar studs, ensuring that the clips clicks in to place and are secure, followed by the shoulder clips.



8

To ensure safe use, place the large hole of the clip over the hoist bar stud, pull down on the clip until it clicks into the narrower neck area. Ensure that the clip cannot be pulled from the bar. To remove the clip, pull the small strap up fully until the large hole is located over the stud. Failure to fully lift the clip over the stud when removing can damage the clip retainer prongs

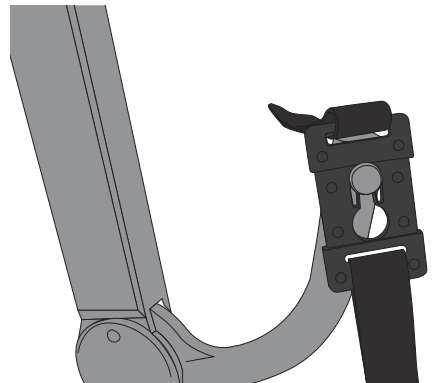

## EXAMINATION AND INSPECTION OF SLING

### CAUTION

Before every lift visually examine the slings structural integrity. If there is any sign of excessive damage to the fabric or straps **DO NOT USE THE SLING.**

### GENERAL EXTERNAL WEAR

In normal use this is unavoidable and is shown by fluffiness to the surface fibres. This may be caused by sharp edges on straps or Velcro. This is harmless unless it becomes too extensive. Normal material condition is smooth or flat.

### LOCAL ABRASION IN THE WEBBING STRAPS

This may be caused by the passage of the webbing over the sharp edges whilst under tension or Velcro. Slight damage to the outer fibres may be considered safe but in severe cases this should lead to instant rejection of the product.

### CUTS, HOLES OR BURNS


These are dangerous and must lead to instant rejection of the product.

### EXAMINATION OF STITCH PATTERNS


Examine the stitch patterns for broken, worn, pulled or unravelled stitches. Reject any product where the thread has become subject to excessive abrasion or broken stitches.



**DO NOT USE**

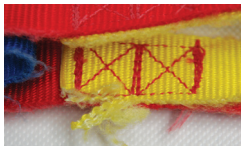

**DO NOT USE**

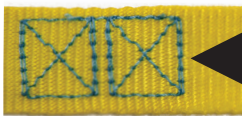

Critical stitch area on all lifting straps  
**THIS MUST BE INTACT**

## WASHING INSTRUCTIONS

Each sling can be laundered up to 80°C. Use a mild soap solution and then rinse. Slings will air dry quickly in a warm room or in the normal manner. Tumble dry on low setting only. Do not place over a convector heater or steam pipes.

Do not use biological detergents. Do not use bleach. Do not use strong chemicals. A faded label is an early indicator that the sling has been subject to harsh/abrasive cleaning methods: do not continue to use the sling. All slings fitted with a Velcro body belt should have the Velcro square removed during laundering and refitted afterwards. This will significantly extend the lifespan of the product.

### SAFE WORKING LOAD

Supersoft (Spacer) - 300kg | Polyester - 400kg | Mesh 300kg | Parachute 200kg

### SLING INSPECTION

Each product should be inspected prior to each use and with your six month or annual hoist inspection. This is in line with LOLER (Lifting Operations and Lifting Equipment Regulations)

### SAFETY AND QUALITY STANDARDS

Care-Ability Healthcare (manufacturer), are registered with the Medicines and Health Care Regulatory Agency of the UK Government to manufacture Class 1 Medical Devices. Our products are categorised as a Class 1 Medical Device and carry the CE mark under the European Medical Device Directive EU93-42. The products are designed and manufactured under the externally accredited BS EN ISO 9001:2008 quality system well beyond the requirements of BS EN 10535:2006, The Hoist and Sling Standard.

All of our products are thoroughly inspected by two staff members prior to the product being dispatched. Our Lifestyle slings are labeled with a unique serial number to assist with traceability and for recording LOLER examinations.

Harvest Healthcare, Sheaf House, Bradmarsh Way, Bradmarsh Business Park, Rotherham, S601BW

t: 01709377172 | f: 01709377172 | e: enquiries@harvesthealthcare.co.uk | www.harvesthealthcare.co.uk