




KINETIC AIR MATTRESS

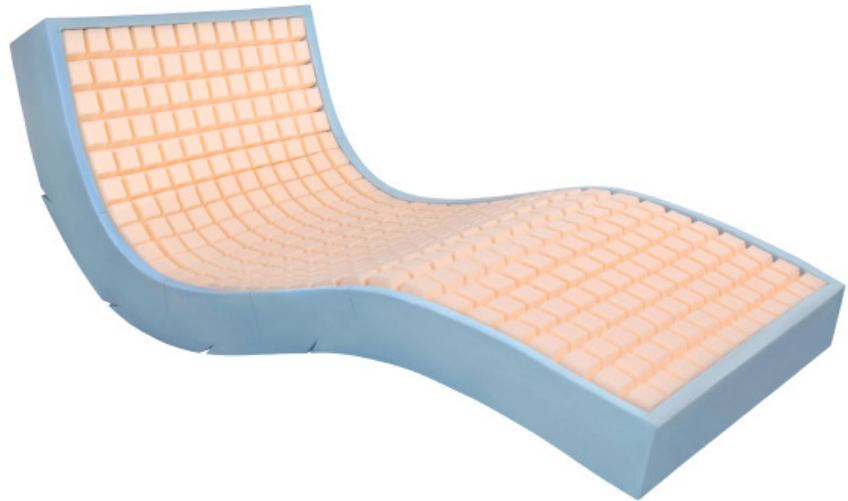

Product Codes: HFK100

- Description:**
- The Kinetic Air replacement mattress consists of a firm foam outer enclosing 11 interconnected foam filled air cells lying below a soft castellated foam topper
 - It combines air displacement technologies with the benefits of a static foam mattress.
 - As the patient moves, internal air pressure is redistributed in the cells through a system of valves.
 - Distribute patient weight and provide the comfort and support. The foam is key cut to allow the mattress to profile to the bed format.

Mattress Dimensions: 2000 x 900 x 150mm

Max. User Weight: 39 Stone (247kg)

Risk Level: High

Cover: (Available Separately) Welded, 2 way-stretch, water proof, PU coated cover in plum, zipped to a durable PVC coated base section. Silver-based antibacterial properties in the PU coating are wash proof and environmentally friendly. Cover is latex free.

Cleaning Guidance

The mattress foam may be hand washed with mild detergent. Do not immerse in water. Allow to dry naturally. Cover should be wiped down with hot water and mild detergent or machine washed up to 95°C and tumble dried. For a comprehensive list of tested cleaning products see our website.

Warranty:

Subject to the terms of Harvest Healthcare Guarantee and Warranties the Kinetic Air mattress is guaranteed for three years from date of purchase.

Exclusions:

1. Removable covers are guaranteed for one year
2. Damage to fabric covers through incorrect use and penetration by sharp instruments is excluded

Supplier Contact Details:

t: +44(0)1709 377172 e: enquiries@harvesthealthcare.co.uk; w: www.harvesthealthcare.co.uk
Sheaf house, Bradmarsh Way, Bradmarsh Business Park, Rotherham UK, S60 1BW

Important Safety Notice

Harvest Healthcare Ltd supplies equipment and appliances to cover a range of medical conditions. Specific choice of equipment is the sole responsibility of the clinical practitioner or carer overseeing individual patient care. Harvest Healthcare Ltd supplies a range of products such as foam mattresses, active air mattresses, fall protection mats and bed rail protection bumpers. These products are usually covered in vapour permeable PU fabric that could cause risk of suffocation in the event of prolonged facial contact.

(Mattress risk categories are to be used as a guide only) Risk assessments specific to user must be carried out prior to the use of all products.

Products and data continually undergo development and may therefore deviate from the details given.