




harvesthealthcare®

equipped to care

Static Mattresses

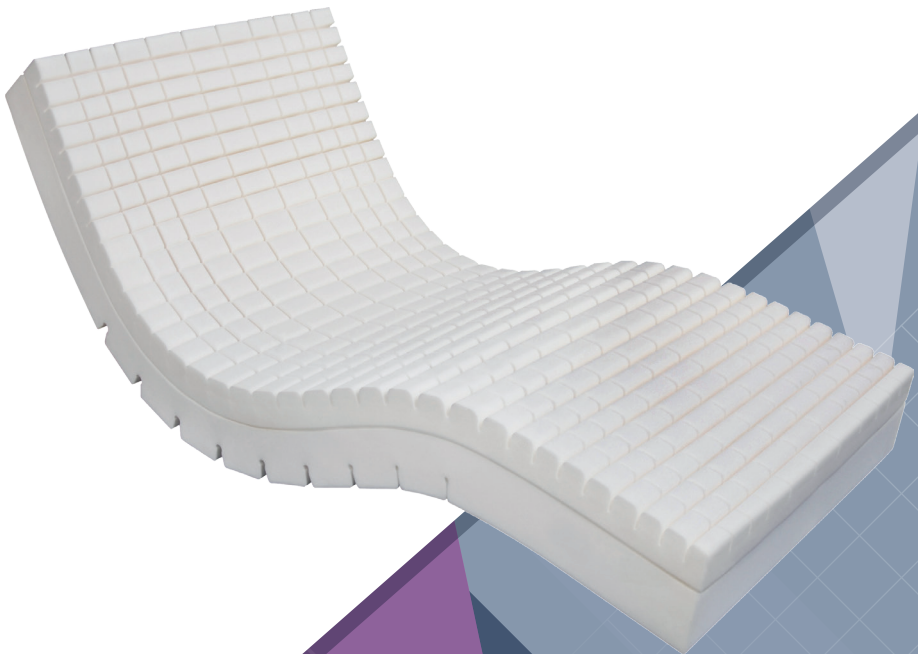

www.harvesthealthcare.co.uk

POWER TO
DELIVER

HARVEST HEALTHCARE
YOUR FIRST CHOICE...



TOTAL SUPPORT TEAM

We offer a full package to support you from installation to service and repair.



HARVEST FAST TRACK NATIONWIDE DELIVERY

We use our own transport for many deliveries which gives us the flexibility to meet your needs.



SUPERIOR BUILD QUALITY

Our products are built to provide years of service by using the latest materials and fabrication techniques.



KNOWLEDGEABLE SALES TEAM

Our experienced sales team is always available to help you find the best possible solution.



INNOVATIVE WELL DESIGNED PRODUCTS

Our products are designed to meet the ever changing needs of our customers. We continually provide new solutions.



EXCEPTIONAL VALUE FOR MONEY

Through careful design and production processes you can be assured that our products will provide value for money over their working life.



EXTENSIVE PRODUCT RANGES

To give you the best solution we maintain an extensive range of products for complete flexibility.

CONTENTS

STATIC MATTRESSES

3

STATIC MATTRESSES

Static Mattresses


04

Hybrids

15

Accessories

20

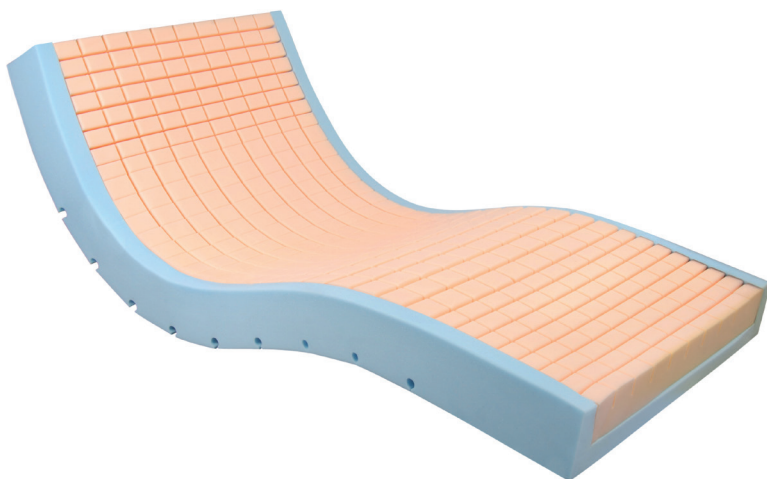

COMMUNITY

HSMCOM - COMMUNITY

The Community Static Mattress is a cost-effective comfort mattress for use on profiling beds.

The foam is highly resilient and fire retardant.

The Community mattress provides pressure distribution and can be turned to extend service life.

KEY FEATURES

Cost-effective

Highly Resilient foam

Four-way stretch, waterproof, anti-bacterial zipped cover with a crib 5 fire rating

Shaped to allow profiling

Two-way turn to extend life of the mattress

Additional sizes available upon request

Medium risk

Comes with a 3-year warranty

Max. patient weight: 25 Stone / 159kg



LOW MED HIGH V HIGH

PRESSURE RELIEF RISK

CODE	DIMENSIONS	MATTRESS WEIGHT
HSMCOM	2000 x 900 x 125mm	8kg
HSMCOM.12	2000 x 1200 x 125mm	10kg

COMMUNITY PLUS

HSMCOMP - COMMUNITY PLUS

The Community Plus Static Mattress is built on the Community mattress with higher quality but remaining cost-effective.

The foam is thicker offering better pressure distribution and support but keeping the same level of comfort.

KEY FEATURES

Cost-effective

Highly Resilient foam

Four-way stretch, waterproof, anti-bacterial zipped cover with a crib 5 fire rating

Shaped to allow profiling

Two-way turn to extend life of the mattress

Additional sizes available upon request

High risk

Comes with a 3-year warranty

Max. patient weight: 30 Stone / 190kg



LOW MED HIGH V HIGH

PRESSURE RELIEF RISK

CODE	DIMENSIONS	MATTRESS WEIGHT
HSMCOMP	2000 x 900 x 150mm	8.5kg
HSMCOMP.12	2000 x 1200 x 150mm	11kg

PRIME COMFORT

HSM PRIM - PRIME COMFORT

The Prime Comfort mattress is made of high resiliency foam offering excellent pressure relief. The castellated foam aids comfort and facilitates air movement inside the mattress, reducing heat build up.

KEY FEATURES

Excellent comfort and support

Four-way stretch, waterproof, anti-bacterial PU cover with welded seams and a crib 5 fire rating

Additional sizes available upon request

High risk

Comes with a 3-year warranty

Max. patient weight: 28 Stone / 177kg



LOW MED HIGH V HIGH

PRESSURE RELIEF RISK

CODE	DIMENSIONS	MATTRESS WEIGHT
HSM PRIM	2000 x 900 x 150mm	9 kg
HSM PRIM.12	2000 x 1200 x 150mm	12 kg

PRIME COMFORT PLUS

STATIC MATTRESSES

7

STATIC MATTRESSES

HSMPRIMP - PRIME COMFORT PLUS

The Prime Comfort Plus provides excellent patient comfort and support.

A castellated foam core facilitates pressure redistribution and air movement inside the mattress, reducing heat build up. It has a firm foam side wall encasing a resilient central insert, enhancing stability and assisting patient transfer.

Slots in the base enable the mattress to profile.

KEY FEATURES

Excellent comfort and support

Slots in the base enable profiling

Four-way stretch, waterproof, anti-bacterial PU cover with welded seams and a crib 5 fire rating

Additional sizes available upon request

High risk

Comes with a 3-year warranty

Max. patient weight: 39 Stone / 247kg



LOW MED HIGH V HIGH

PRESSURE RELIEF RISK

CODE	DIMENSIONS	MATTRESS WEIGHT
HSMPRIMP	2000 x 900 x 150mm	13kg
HSMPRIMP.1200	2000 x 1200 x 150mm	17kg

FIRE-TECH 7 COMFORT

HSMFIRET - FIRE-TECH 7 COMFORT

The Fire-Tech 7 Comfort is a pressure care mattress for environments at increased risk of fire. The mattress has a Crib 7 fire rating which is the most demanding fire test in the industry. It provides excellent pressure re-distribution and comfort for patients at high risk of developing pressure related injury.

KEY FEATURES

Excellent comfort and support

Four-way stretch, waterproof, anti-bacterial PU cover


Available welded or zipped

Crib 7 fire rating

Additional sizes available upon request
High risk

Comes with a 3-year warranty

Max. patient weight: 30 Stone / 190kg

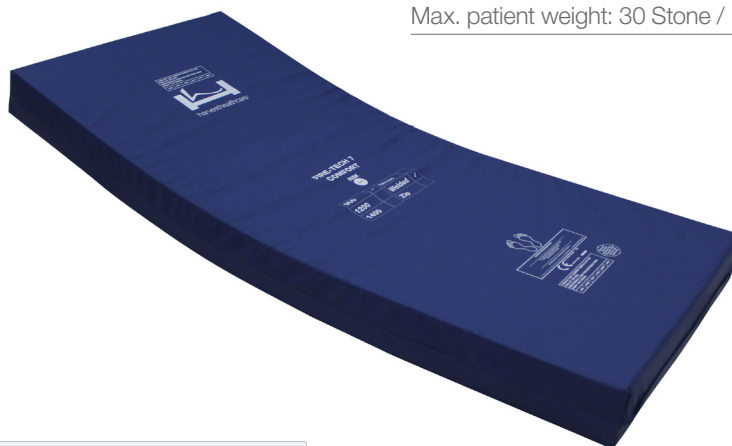

LOW MED HIGH V HIGH

PRESSURE RELIEF RISK

CODE	DIMENSIONS	MATTRESS WEIGHT
HSMFIRET.CFT	2000 x 900 x 150mm	8.5 kg
HSMFIRET.CFT.W	2000 x 1200 x 150mm	11.5 kg

VANDAL PROOF MATTRESSES

VANDAL PROOF MATTRESSES

Vandal proof mattresses are designed for challenging environments where you need a mattress that is fully scure and fire resistant.

KEY FEATURES

Slots in the base enable profiling

Four-way stretch, waterproof, anti-bacterial PU cover with welded seams

Available in both Crib 5 & Crib 7 fire rating

Additional sizes available upon request

High risk

Comes with a 3-year warranty

Max. patient weight: 30 Stone / 190kg

HSMCOM.VP: 25 Stone / 159kg

HSMCOMP.VP: 30 Stone / 190kg

HSMFIRET.VP: 30 Stone / 190kg



LOW MED **HIGH** V HIGH

PRESSURE RELIEF RISK

CODE	DIMENSIONS	MATTRESS WEIGHT	FIRE RATING
HSMCOM.VP	2000 x 900 x 125mm	13kg	5
HSMCOMP.VP	2000 x 900 x 150mm	8.5kg	5
HSMFIRET.VP	2000 x 900 x 150mm	8.5kg	7

REFLECT 2

HFR147 - REFLECT 2

The Reflect 2 mattress is a profiling mattress with a castellated viscoelastic top and a firm base.

The castellated viscoelastic foam offers improved immersion and pressure distribution. The firm base gives support and helps with even distribution of weight. There are keyhole slots in the base to aid profiling.

KEY FEATURES

Castellated viscoelastic top and firm foam base

Slots in the base enable profiling

Four-way stretch, waterproof, anti-bacterial PU cover with welded seams and a crib 5 fire rating

Additional sizes including plus sizes available upon request

High risk

Comes with a 3-year warranty

Max. patient weights:

HFR147: 39 Stone / 247kg

HFR227: 45 Stone / 286kg

HFR227.60: 60 Stone / 381kg



LOW

MED

HIGH

V HIGH

PRESSURE RELIEF RISK

CODE	DIMENSIONS	MATTRESS WEIGHT
HFR147	2000 x 900 x 165mm	16kg
HFR147.12	2000 x 1200 x 165mm	22kg
HFR227	2000 x 1200 x 180mm	30kg
HFR227.60	2000 x 1200 x 180mm	36kg

INTRODUCING SYSTEM

STATIC MATTRESSES

11


STATIC MATTRESSES

The System range is made of high quality moulded foam offering high maximum patient weights and risk levels. Extra features create superior levels of comfort and pressure relief.

Anatomical Shape Insert


Regular Viscoelastic Foam Mattress

The viscoelastic foam is strongly compressed in the area at risk, especially when profiled.

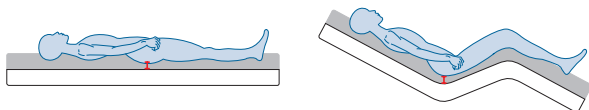

System Mattresses

As the insert follows the body shape in the sacrum area, the thickness of the viscoelastic foam is maximised to foster better prevention in areas at risk.

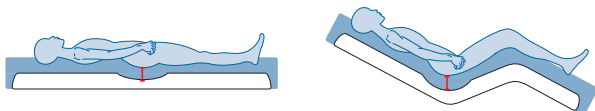

Multibearing Evolutive Concept


The multibearing surface enables a better pressure distribution through the transfer of support from areas with low bearing capacity toward areas with high bearing capacity.

White Area:


Areas with low bearing capacity (softer foam) giving increased prevention to areas at risk of pressure ulcerations (heels).

Blue Areas


Areas with high bearing capacity (firmer foam) for areas of the body with a lower level of pressure ulceration risk.

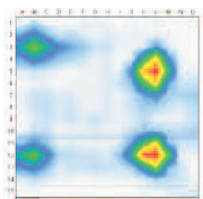

Viscoelastic Foam
Mattress

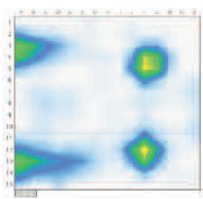

Viscoflex
Evolutive

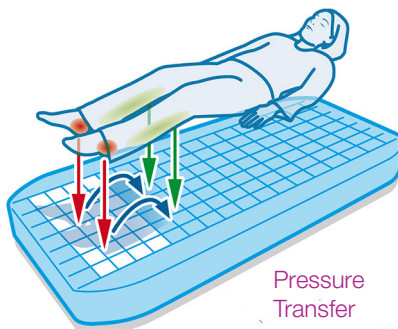

Pressure
Transfer

VISCOFLEX

P161MFP - VISCOFLEX

The Viscoflex mattress is a high-resilience mattress, moulded to ensure bonding. The mattress has an anatomically shaped base with a memory foam top.

The very high density top increases the contact area and improves pressure distribution with increased depth in sacral areas. The high density base works with the top to enhance support and comfort.

KEY FEATURES

Very high density viscoelastic foam

Uniquely shaped to accommodate bed accessories

Two-way stretch, waterproof, anti-bacterial PU coated cover with welded seams and a crib 5 fire rating


Comes with a 3-year warranty

High to very high risk

Max. patient weights:

P161MFPMHF: 20 stone / 130kg

P161MFPMHF120: 42 stone / 270kg

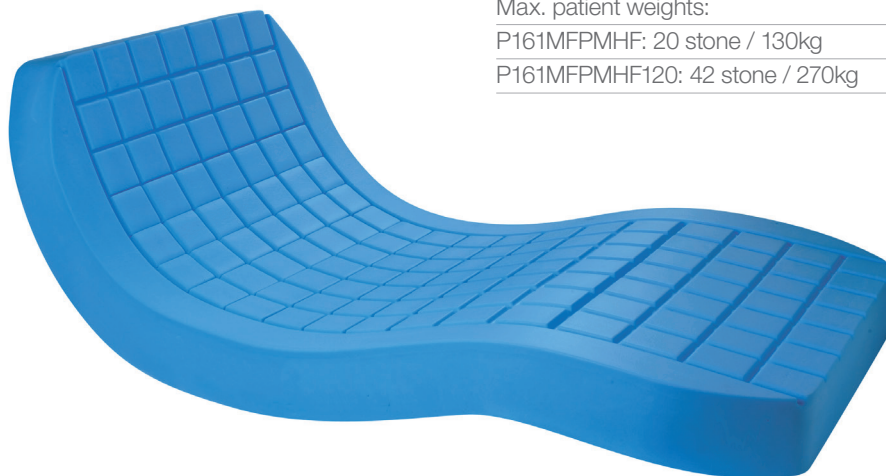

LOW MED HIGH **V HIGH**

PRESSURE RELIEF RISK

CODE

DIMENSIONS

P161MFPMHF	1990 x 860 x 140mm
P161MFPMHF120	1970 x 1150 x 140mm

VISCOFLEX WITH MULTIBEARING SURFACE

P161MOP VISCOFLEX WITH MULTIBEARING SURFACE

The Viscoflex Mattress with Multibearing Surface has all the great features of the Viscoflex mattress with an added multibearing surface.

The multibearing surface enables better pressure distribution through the transfer of supports from the areas of low bearing capacity toward areas with high bearing capacity.

KEY FEATURES

Very high density viscoelastic foam

Added multibearing surface


Uniquely shaped to accommodate bed accessories

Two-way stretch, waterproof, anti-bacterial PU coated cover with welded seams and a crib 5 fire rating

Comes with a 3-year warranty

High to very high risk

Max. patient weight: 20stone / 130kg

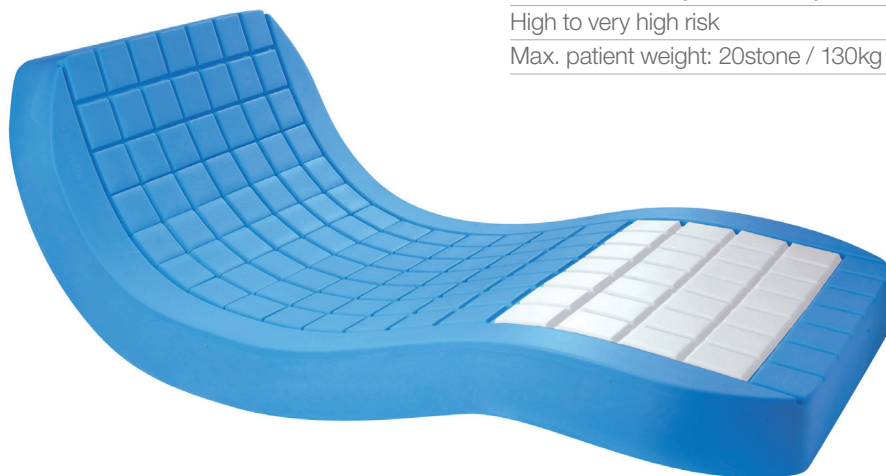

LOW MED HIGH **V HIGH**

PRESSURE RELIEF RISK

CODE

DIMENSIONS

P161MOPPMHF

1990 x 860 x 140mm

VISCOFLEX EVOLUTIVE

P161MFE VISCOFLEX EVOLUTIVE

The Viscoflex Evolutive Static Mattress is our enhanced mattress in the Viscoflex range.

Positioning specialist foam inserts across key areas of the mattress allows for more pressure relief and comfort. The varying firmness in foam at zonal areas contributes to the additional support.

The visco-elastic memory foam moulds to the users' body shape and the high resilience foam base increases stability.

The multibearing surface enables better pressure distribution through the transfer of supports from the areas of low bearing capacity toward areas with high bearing capacity.

KEY FEATURES

Very high density viscoelastic foam

Added multibearing surfaces

Uniquely shaped to accommodate bed accessories

Supplied a with two-way stretch, fully welded, waterproof, PU coated, polyester cover

Comes with a 3-year warranty

Very high risk

Max. patient weight: 20 stone / 130kg



LOW

MED

HIGH

V HIGH

PRESSURE RELIEF RISK

CODE

P161MFEPMHF

DIMENSIONS

1970 x 860 x 140mm

VISCOFLEX AIR

P161MOR VISCOFLEX AIR

The Viscoflex Air Mattress is built on the same base as the Viscoflex mattress, but with an enhanced performance.

Polyair cells have been added to the heel area, and a multi-bearing surface to the lower back area. This helps to improve the contact surface and pressure distribution through customised adjustments of the air sections.

Viscoelastic foam moulds to the shape of the body, which not only improves comfort, but also pressure distribution.

KEY FEATURES

Very high density viscoelastic foam

Different firmness in foam to zonal areas relative to each part of the body

Laprexene cells with high elasticity and tearing resistance


Uniquely shaped to accommodate bed accessories

Supplied with a two-way stretch, fully welded, waterproof, PU coated, polyester cover

Comes with a 3-year warranty

Very high risk

Max.patient weight: 20 stone / 130kg

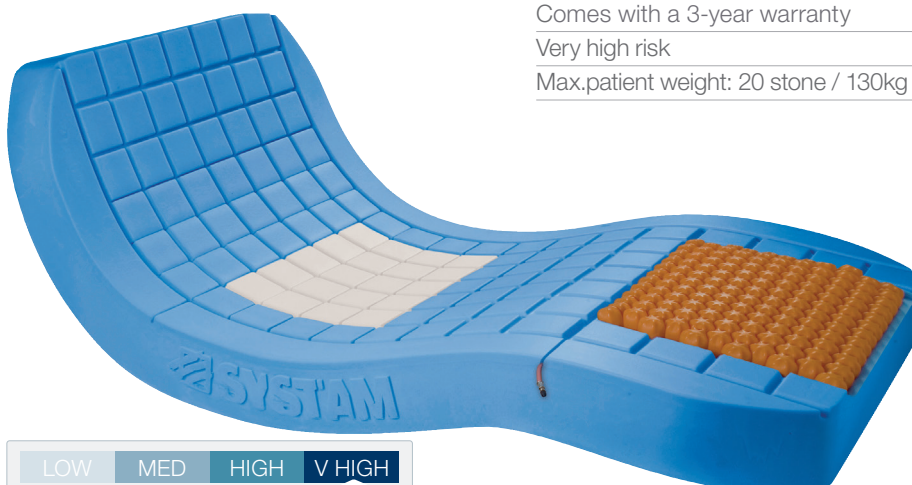

LOW

MED

HIGH

V HIGH

PRESSURE RELIEF RISK

CODE

P161MORPMHI

DIMENSIONS

1970 x 860 x 140mm

VISCOFLEX AIR PLUS AIR

P161MO2 VISCOFLEX AIR PLUS AIR

The Viscoflex Air Plus Air Mattress is built on the Viscoflex mattress base, but with an enhanced performance.

Polyair cells have been added to the heel and lower back areas, improving contact surface and pressure distribution through the customised adjustments of the air sections.

Viscoelastic foam moulds to the shape of the patient's body, which not only improves comfort, but also improves pressure distribution.

KEY FEATURES

Very high density viscoelastic foam

Laprexene cells with high elasticity and tearing resistance

Uniquely shaped to accommodate bed accessories

Supplied with two-way stretch, fully welded, waterproof, PU coated, polyester cover

Comes with a 3-year warranty

Very high risk

Max patient weight: 20 stone / 130kg



LOW

MED

HIGH

V HIGH

PRESSURE RELIEF RISK

CODE

DIMENSIONS

P161MO2RPMHI

1970 x 860 x 140mm

INTRODUCING HYBRIDS

HYBRID MATTRESSES

17

HYBRID MATTRESSES

Harvest Healthcare offer hybrid mattresses to suit our customers' needs. Hybrids enable carers to step up or step down the treatment of patients by either adding or removing a pump. This offers flexibility when caring for patients that could become at risk of a pressure ulcer. It means that users can benefit from the high levels of comfort that foam offers and the pressure relieving qualities of an active mattress when needed. Using one mattress for both functions helps reduce logistic costs and avoids disruption for the user as there is no need to change the mattress.



KINETIC AIR

HFK100 - KINETIC AIR

The Kinetic Air replacement mattress consists of a firm foam outer enclosing 11 interconnected foam filled air cells lying below a soft castellated foam topper

FEATURES

It combines air displacement technologies with the benefits of a static foam mattress.

As the patient moves, internal air pressure is redistributed in the cells through a system of valves.

Distributes patient weight, providing comfort and support.

The foam is key cut to allow the mattress to profile to the bed format.

Comes with a 3 year warranty

Max. patient weight: 39 stone / 247kg




LOW MED HIGH V HIGH

PRESSURE RELIEF RISK

TECHNICAL SPECIFICATIONS

Dimensions: 2000 x 900 x 150mm

High Risk

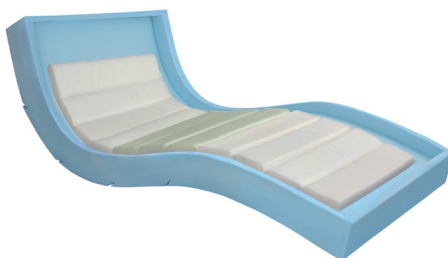

PRIME COMFORT ACTIVE

HHMPRIMA - PRIME COMFORT ACTIVE

The Prime Comfort Active mattress combines our Prime Comfort Plus with an active air cell component.

FEATURES

Firm 'U' shaped high density foam base for roll prevention and to aid transfers

Castellated soft and resilient central foam insert

10 alternating cells powered by the low noise level pump

Fully welded, antibacterial, waterproof, two way stretch cover

Durable PVC base

Step up step down with/without pump

Comes with a 3 year warranty

Compatible with all Harvest pumps

Max. patient weight: 39 Stone / 247kg



TECHNICAL SPECIFICATIONS

Dimensions: 2000 x 900 x 140mm

Weight: 14kg

High Risk (without pump)

Very High Risk (with pump)

Maximum Patient Weight: 39 Stone/247kg

LOW

MED

HIGH

V HIGH

PRESSURE RELIEF RISK

AIR-TECH

HHMAIRT001 - AIR-TECH REPLACEMENT

The Air-Tech hybrid mattress combines the comfort of a static mattress with the pressure relief of alternating cells. Foam within cells allows for the mattress to be stepped up or down giving flexibility to cover all nursing requirements.

KEY FEATURES


Firm 'U' shaped foam base for stability, roll prevention and to aid transfers.

Castellated soft and resilient central insert with zonal areas designed for support and maximum pressure relief.

Fully welded, antibacterial, waterproof, two-way stretch cover with a durable PVC base.

Powered by the quiet and efficient Harvest 4 pump.

Designed and manufactured in the UK

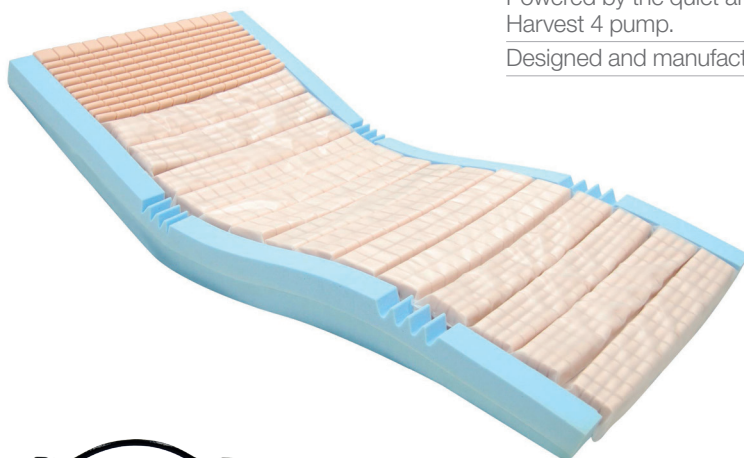

HARVEST 4 PUMP SPECIFICATIONS

Code: HPU4

Dimensions: 265 x 175 x 90mm

Cycle time: 12 minutes

LOW MED HIGH V HIGH

PRESSURE RELIEF RISK

MATTRESS TECHNICAL SPECIFICATIONS

Dimensions: 1990 x 900 x 150mm

Very High Risk

Max. Patient Weight With Pump: 39 Stone / 247kg

AIR-TECH ENHANCED HEEL PROTECTION

HHMAIRT002 - AIR-TECH ENHANCED HEEL PROTECTION SYSTEM

The Air-Tech hybrid mattress with heel protection offers the same great features of the Air-Tech mattress but with added heel protection to help prevent pressure ulcers in the heels. The tapered heel section is designed to reduce heel pressure in both static and dynamic modes. In dynamic mode air circulates in the cells to off load pressure in the heels area. It is a **very high risk** mattress.

KEY FEATURES

Zero heel pressure

20mm 12 cell MicroCell heel zone


Firm 'U' shaped foam base for stability, roll prevention and to aid transfers.

Castellated soft and resilient central insert with zonal areas designed for support and maximum pressure relief.

Fully welded, antibacterial, waterproof, two-way stretch cover with a durable PVC base.

Powered by the quiet and efficient Harvest 4 pump.

Designed and manufactured in the UK

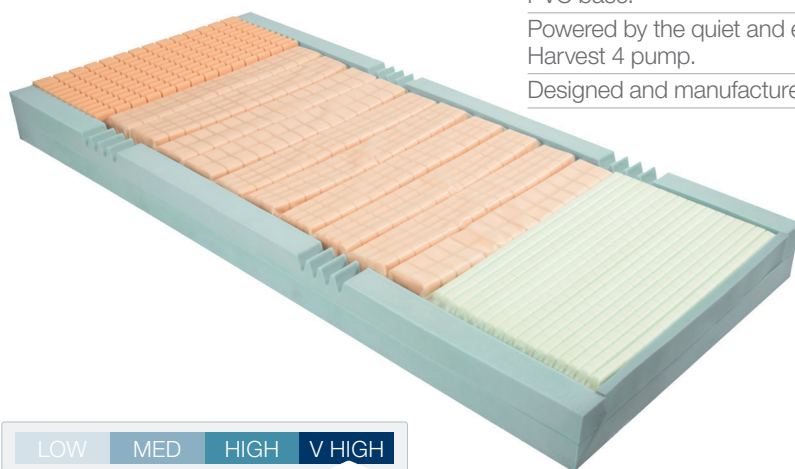

LOW MED HIGH V HIGH

PRESSURE RELIEF RISK

MATTRESS TECHNICAL SPECIFICATIONS

Dimensions: 1990 x 900 x 150mm

Very High Risk

Maximum Patient Weight With Pump: 39 Stone / 247kg



NOTES

NOTES

NOTES

23

NOTES

NOTES


25

ALSO AVAILABLE

ALSO AVAILABLE

26

ALSO AVAILABLE

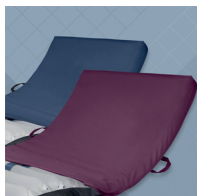

ACTIVE MATTRESSES




PROFILING BEDS AND ACCESSORIES



ACTIVE CUSHIONS



STATIC CUSHIONS

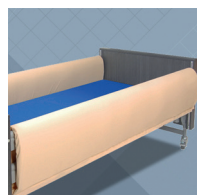


BED RAIL BUMPERS

ALSO AVAILABLE

ALSO AVAILABLE

27

ALSO AVAILABLE

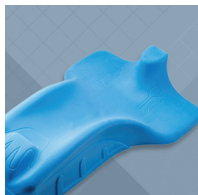

HEELS PADS AND POSITIONING AIDS



SLINGS, HOISTS AND STAND AIDS



EVACUATION PRODUCTS



FALL REDUCTION



SHOWERING AND TOILET AIDS



harvesthealthcare®

Sheaf House, Bradmarsh Business Park,
Bradmarsh Way, Rotherham S60 1BW

T +44 (0)1709 377172

F +44 (0)1709 377173

E sales@harvesthealthcare.co.uk

www.harvesthealthcare.co.uk

© Copyright Harvest Healthcare

Harvest Healthcare Limited. Company No: 07210261

Our Full Terms & Conditions are available by request or can be found on our Website.
Harvest Healthcare reserves the right to alter or amend this document without notice.