


MIXTE CUSHION

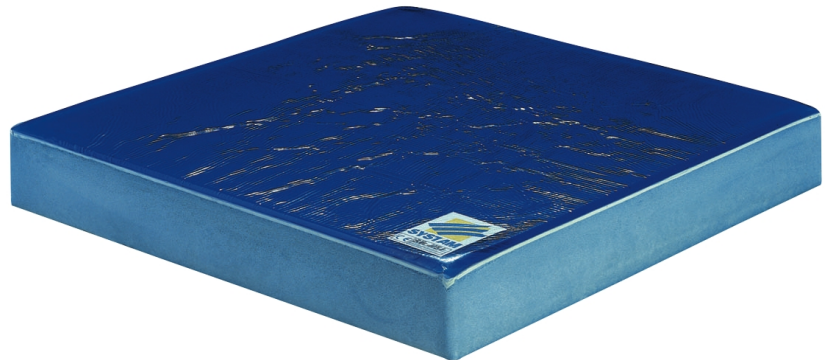

Product Code:

P314C40401HF

Description:

- Cushion of viscoelastic polyurethane gel and light resilience foam moulded together.
- While the gel resists penetration of the layer, it also assists pressure distribution.
- The thick high resilience foam base increases contact area for better pressure distribution, whilst improving stability and comfort.
- Supplied complete with a welded navy blue cover.
- Available in range of size and weight capacities to order.

Dimensions:

400 x 400 x 60mm

Weight of Cushion:

2.4 kg

User Weight:

30-90 kg

Pressure Relief Risk:

Medium to High Risk

Foam Density:

Base:

40 kg/m³

Cover:

Anti-allergenic, antibacterial, antimicrobial, antifungal, 2 way stretch, vapour permeable, washable, welded seams, 2 side zip. Printed user instructions on cover.

Fire Retardancy:

NF EN ISO 597 - 1 & 2

Cleaning Guidance:

The cushion foam may be hand washed with mild detergent. Do not immerse in water. Allow to dry naturally. Cover should be wiped down with hot water and mild detergent or machine washed up to 95°C and tumble dried. For a comprehensive list of tested cleaning products see our website.

Warranty

In terms of Harvest Healthcare Guarantee and Warranties the Mixte cushion is guaranteed for three years from date of purchase.

Exclusions:

1. Removable cushion covers are guaranteed for one year
2. Damage to fabric covers through incorrect use and penetration by sharp instruments is excluded

Supplier Contact Details:

t: +44(0)1709 377172 e: enquiries@harvesthealthcare.co.uk w: www.harvesthealthcare.co.uk
Sheaf house, Bradmarsh Way, Bradmarsh Business Park, Rotherham UK, S60

Important Safety Notice

Harvest Healthcare Ltd supplies equipment and appliances to cover a range of medical conditions. Specific choice of equipment is the sole responsibility of the clinical practitioner or carer overseeing individual patient care. Harvest Healthcare Ltd supplies a range of products such as foam mattresses, active air mattresses, fall protection mats and bed rail protection bumpers. These products are usually covered in vapour permeable PU fabric that could cause risk of suffocation in the event of prolonged facial contact.
(Mattress risk categories are to be used as a guide only)

Risk assessments specific to user must be carried out prior to the use of all products. Products and data continually undergo development and may therefore deviate from the details given.

Version 2: Oct 2017