


Sling Assessment

Selecting the appropriate sling size for your resident

Name _____ Company: _____
 Address: _____
 Telephone: _____ Email: _____

HSL200 TOILET/DRESS - LOOPS

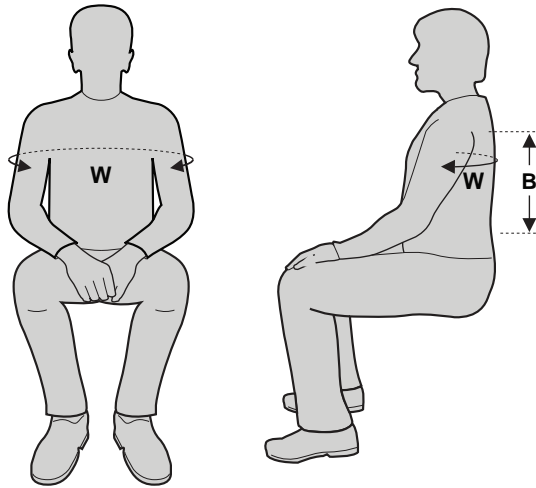

All measurements are in cm.

Using the Sling Sizing Chart

- Step 1:**
Measure the width of your resident between the outer midpoints of each upper arm as per the diagram (W).
- Step 2:**
Input your measurement into the corresponding box under 'User Width', next to the width range it fits within.
- Step 3:**
Measure the back length of your resident as per the diagram (B).
- Step 4:**
Input your measurement into the corresponding box under 'User Back Length', next to the width range it fits within. Please note some ranges overlap - in these cases select the most appropriate size.
- Step 5:**
If these two measurements correlate your resident fits a standard sling. If not, consider a bespoke sling or select the most appropriate standard sling after a risk assessment.

All sling selections should be accompanied by a risk assessment by a trained person. People of atypical anatomy may not be suited to a standard sling.

Sling Sizing Chart

HSL200 Size	User Width (W)	User Back Length (B)
Extra Small - Black	40-50	26-30
Small - Red	50-60	29-34
Medium - Yellow	62-68	34-40
Large - Green	68-78	36-42
Extra Large - Blue	78-85	38-44

Sling size selected for resident: _____

Loop Sling

Clip Sling

Material **Polyester** **Super Soft**

Comments / Additional info: _____

Date: _____ Signature: _____



t +44 (0)1709 377172 f +44 (0)1709 377173
 e sales@harvesthealthcare.co.uk
 www.harvesthealthcare.co.uk