

Sling Assessment


Selecting the appropriate sling size for your resident

Name _____ Company: _____

Address: _____

Telephone: _____ Email: _____

HSL500/510 STAND AID - LOOPS/CLIPS

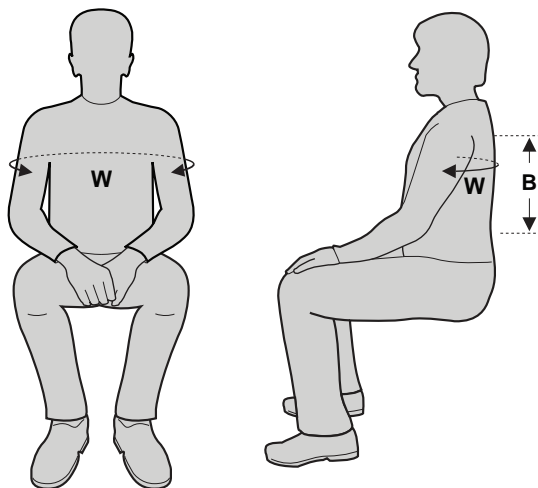

All measurements are in cm.

Using the Sling Sizing Chart

Step 1:
Measure the width of your resident between the outer midpoints of each upper arm as per the diagram (W).

Step 2:
Input your measurement into the corresponding box under 'User Width', next to the width range it fits within.

Step 3:
Measure the back length of your resident as per the diagram (B).

Step 4:
Input your measurement into the corresponding box under 'User Back Length', next to the width range it fits within. Please note some ranges overlap - in these cases select the most appropriate size.

Step 5:
If these two measurements correlate your resident fits a standard sling. If not, consider a bespoke sling or select the most appropriate standard sling after a risk assessment.

All sling selections should be accompanied by a risk assessment by a trained person. People of atypical anatomy may not be suited to a standard sling.

Sling Sizing Chart

HSL500/510 Size	User Width (W)	User Back Length (B)
Extra Small - Black	40-50	28-32
Small - Red	50-60	30-36
Medium - Yellow	60-70	33-39
Large - Green	70-80	37-42
Extra Large - Blue	80-90	40-46

Sling size selected for resident: _____

Loop Sling

Clip Sling

Material **Polyester** **Super Soft**

Comments / Additional info: _____

Date: _____

Signature: _____