


Sling Assessment

Selecting the appropriate sling size for your resident

| | |
|------------|----------|
| Name | Company: |
| Address: | |
| Telephone: | Email: |

HSL400/410 IN CHAIR HAMMOCK - LOOPS/CLIPS

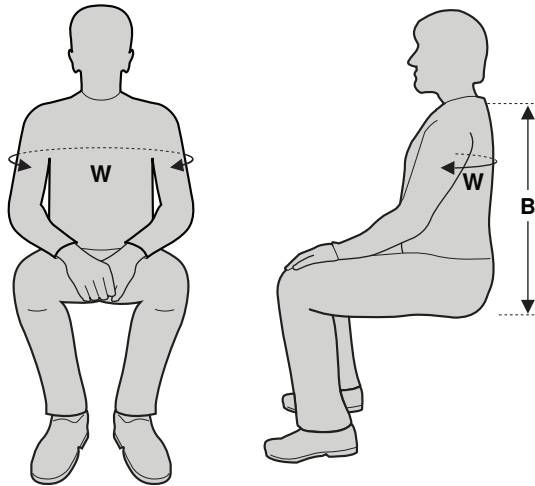

All measurements are in cm.



Using the Sling Sizing Chart

Step 1:

Measure the width of your resident between the outer midpoints of each upper arm as per the diagram (W).

Step 2:

Input your measurement into the corresponding box under 'User Width', next to the width range it fits within.

Step 3:

Measure the back length of your resident as per the diagram (B).

Step 4:

Input your measurement into the corresponding box under 'User Back Length', next to the width range it fits within. Please note some ranges overlap - in these cases select the most appropriate size.

Step 5:

If these two measurements correlate your resident fits a standard sling. If not, consider a bespoke sling or select the most appropriate standard sling after a risk assessment.

All sling selections should be accompanied by a risk assessment by a trained person. People of atypical anatomy may not be suited to a standard sling.

Sling Sizing Chart

| HSL400/410 Size | User Width (W) | User Back Length (B) |
|---------------------|---------------------------------|--------------------------------|
| Extra Small - Black | 45-52 <input type="checkbox"/> | 48-56 <input type="checkbox"/> |
| Small - Red | 52-60 <input type="checkbox"/> | 53-63 <input type="checkbox"/> |
| Medium - Yellow | 60-76 <input type="checkbox"/> | 57-67 <input type="checkbox"/> |
| Large - Green | 76-92 <input type="checkbox"/> | 62-72 <input type="checkbox"/> |
| Extra Large - Blue | 92-100 <input type="checkbox"/> | 65-75 <input type="checkbox"/> |

Sling size selected for resident: _____

Loop Sling ☐

Clip Sling ☐

Material

Super Soft ☐

Comments / Additional info:

Date:

Signature: