


REFLECT 2 BARIATRIC MATTRESS

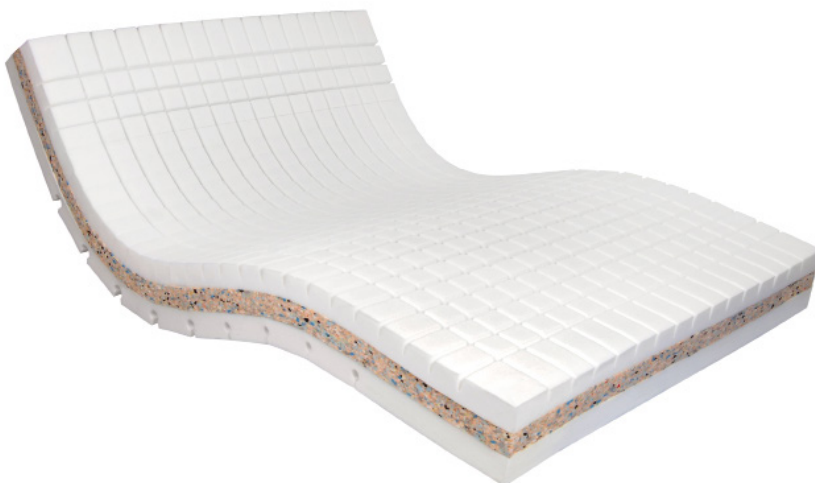

Product Codes:

HFR227
HFR227.14
HFR227.60
HFR227.1400.60

Description:

- Deep layer of castellated viscoelastic foam and a firm foam base combined with a sandwich layer of high load bearing foam allow for greater prevention of pressure ulcers for bariatric patients.
- Distribute patient weight and provide the comfort and support along with the ability to profile.
- Latex free navy blue cover is PU coated, water resistant and antibacterial with a durable PVC base section
- does not require monthly turning, reducing hazardous lifting
- Available in 1200mm and 1400mm widths which make it suitable for our Woburn Ultimate 1200 and 1400 profiling beds.

Mattress Dimensions and Weight, per size:

2000 x 1200 x 180mm	30kg
2000 x 1400 x 180mm	36kg
2000 x 1200 x 180mm	30kg
2000 x 1400 x 180mm	36kg

Max. User Weight:

45 Stone (285kg)
45 Stone (285kg)
60 Stone (381kg)
60 Stone (381kg)

Risk Level:

High

Cover: (Available Separately)

Four way-stretch, water resistant PU coated cover in navy blue, zipped to a durable PVC coated base section. Silver-based antibacterial properties in the PU coating are washable and environmentally friendly. Cover does not contain latex.

Cleaning Guidance

The cushion foam may be hand washed with mild detergent. Do not immerse in water. Allow to dry naturally. Cover should be wiped down with hot water and mild detergent or machine washed up to 95°C and tumble dried. For a comprehensive list of tested cleaning products see our website.

Warranty:

Subject to the terms of Harvest Healthcare Guarantee and Warranties the Reflect 2 Bariatric mattress is guaranteed for three years from date of purchase.

Exclusions:

1. Removable covers are guaranteed for one year
2. Damage to fabric covers through incorrect use and penetration by sharp instruments is excluded

Supplier Contact Details:

t: +44(0)1709 377172 e: enquiries@harvesthealthcare.co.uk; w: www.harvesthealthcare.co.uk
Sheaf house, Bradmarsh Way, Bradmarsh Business Park, Rotherham UK, S60 1BW

Important Safety Notice

Harvest Healthcare Ltd supplies equipment and appliances to cover a range of medical conditions. Specific choice of equipment is the sole responsibility of the clinical practitioner or carer overseeing individual patient care. Harvest Healthcare Ltd supplies a range of products such as foam mattresses, active air mattresses, fall protection mats and bed rail protection bumpers. These products are usually covered in vapour permeable PU fabric that could cause risk of suffocation in the event of prolonged facial contact.

(Mattress risk categories are to be used as a guide only) Risk assessments specific to user must be carried out prior to the use of all products.

Products and data continually undergo development and may therefore deviate from the details given.