

VISCOFLEX MATTRESS



Product Codes:

P161MFPMHF
P161MFPMHF120

Description:

- High resilience and viscoelastic foams moulded in one to ensure bond.
- Anatomically shaped base with memory foam top.
- Very high density viscoelastic top increases contact area and improves pressure redistribution with increased depth in sacral areas.
- High density base works with top to enhance support and comfort.
- Supplied with a fully welded water resistant PU coated stretch polyester cover.
- Uniquely shaped to accommodate bed accessories.
- Delivered compressed and rolled.

Dimensions:

P161MFPMHF: 860 x 1990 mm to suit 900 x 2000 mm (Bed)
P161MFPMHF120: 1150 x 1970 mm to suit 1200 x 2000 mm (Bed)

Max. User Weight:

P161MFPMHF: 20 Stone (130 kg)
P161MFPMHF120: 42 Stone (270kg)

Risk Level:

High to Very High

**Cover:
(Available Separately)
Cover Code:**

Hypoallergenic, waterproof, four way stretch cover with non-slip base and three sided zip P452PF

Foam Density


80kg/m³ foam (visco upper layer)
40kg/m³ foam (firm base layer)

Fire Retardancy:

BS 7177: 2008 crib 5

Cleaning Guidance

The cushion foam may be hand washed with mild detergent. Do not immerse in water. Allow to dry naturally. Cover should be wiped down with hot water and mild detergent or machine washed up to 95°C and tumble dried. For a comprehensive list of tested cleaning products see our website.

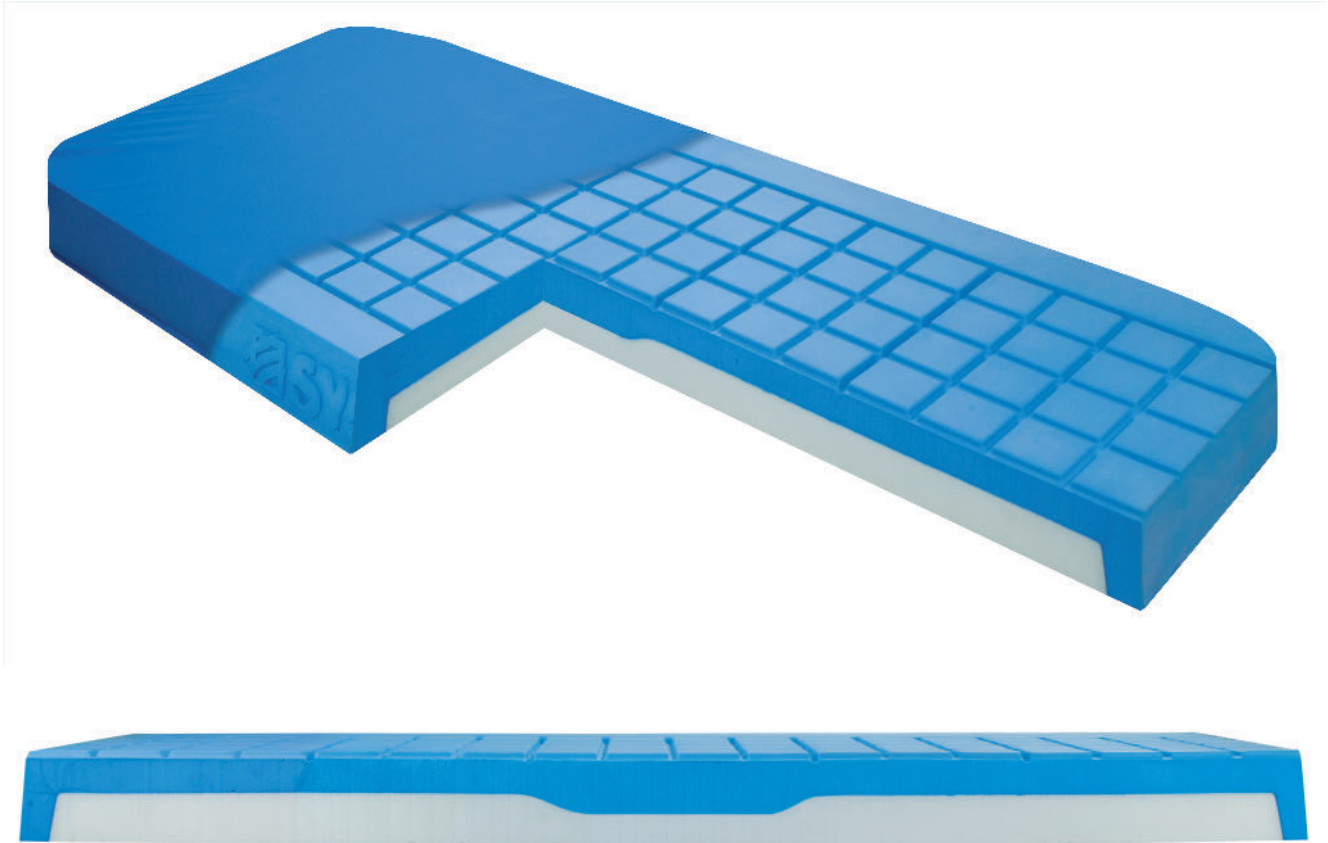

Warranty:

Subject to the terms of Harvest Healthcare Guarantee and Warranties the Viscoflex Mattress is guaranteed for three years from date of purchase.

Exclusions:

1. Removable covers are guaranteed for one year
2. Damage to fabric covers through incorrect use and penetration by sharp instruments is excluded

Supplier Contact Details:

t: +44(0)1709 377172 e: enquiries@harvesthealthcare.co.uk; w: www.harvesthealthcare.co.uk
Sheaf house, Bradmarsh Way, Bradmarsh Business Park, Rotherham UK, S60 1BW

Important Safety Notice

Harvest Healthcare Ltd supplies equipment and appliances to cover a range of medical conditions. Specific choice of equipment is the sole responsibility of the clinical practitioner or carer overseeing individual patient care. Harvest Healthcare Ltd supplies a range of products such as foam mattresses, active air mattresses, fall protection mats and bed rail protection bumpers. These products are usually covered in vapour permeable PU fabric that could cause risk of suffocation in the event of prolonged facial contact. (Mattress risk categories are to be used as a guide only)

Risk assessments specific to user must be carried out prior to the use of all products.

Products and data continually undergo development and may therefore deviate from the details given.

Version 4: March 2018