

harvesthealthcare®
a prismhealthcare company

Harvest Healthcare Equipment Guide

August 2023

www.harvesthealthcare.co.uk

POWER TO DELIVER

HARVEST HEALTHCARE YOUR FIRST CHOICE...

WHO WE ARE & WHAT WE DO

Harvest Healthcare has more than 20 years experience in the design, manufacture and service of pressure care equipment, beds, moving, handling and showering aids.

We offer a variety of support services from mattress repair and decontamination, service and repair provision and client assessment through to project management for installed solutions.

Since 2019 we have been proud to be part of the Prism Healthcare Group, which has enabled us to invest significantly in R&D, expand facilities, and serve our customers even better. Prism Healthcare has an unrivalled range of products and services across all healthcare sectors - see pages 144-145 for further details of the Group, its brands and subsidiaries.



SUPERIOR PRODUCT SOLUTIONS

We provide a complete pressure care equipment, beds and moving & handling product portfolio consisting of static and active (dynamic) mattresses and cushions, profiling beds, hoists, ceiling track hoists, mobile hoists, stand aids, gantry systems, slings, shower chairs, shower stretchers and much more as well as the supply and service of specialist seating.



TOTAL SUPPORT TEAM

We offer a full package to support you from installation to service and maintenance.



HARVEST FAST TRACK NATIONWIDE DELIVERY

We use our own transport for many deliveries which gives us the flexibility to meet your needs.



SUPERIOR BUILD QUALITY

Our products are built to provide years of service by using modern materials and fabrication techniques.



KNOWLEDGEABLE SALES TEAM

Our experienced sales team is always available to help you find the best possible solution.



INNOVATIVE WELL DESIGNED PRODUCTS

Our products are designed to meet the ever-changing needs of our customers. We continually provide new solutions.



EXCEPTIONAL VALUE FOR MONEY

Through careful design and production processes, you can be assured that our products will provide value for money throughout their working life.



EXTENSIVE PRODUCT RANGES

To give you the best solution we maintain an extensive range of products for complete flexibility.



www.harvesthealthcare.co.uk

Sheaf House, Bradmarsh Business Park, Bradmarsh Way, Rotherham S60 1BW
T +44 (0)1709 377172 **F** +44 (0)1709 377173 **E** sales@harvesthealthcare.co.uk

CONTENTS

Your Contacts	03
1. Beds	04
1.1. Woburn Community Profiling Beds	05
1.2. Woburn Profiling Beds	09
1.3. Woburn Ultimate Profiling Beds	14
1.4. Cadence Profiling Beds	16
1.5. Woburn Bed Accessories	20
1.6. Bed Bumpers	26
2. Mattresses and Cushions	33
2.1. Static Mattresses	34
2.2. Hybrid Mattresses	46
2.3. Active Mattresses	48
2.4. Pumps	59
2.5. Fall Mats, Underlays and Overlays	60
2.6. Cushions	61
3. Positioning Aids	74
4. Nurse Call Alarm System	89
5. Moving and Handling	93
5.1. Stand Aids, Hoists and Transfer Aids	94
5.2. Slings	103
5.3. Bed Management Systems	112
5.4. Wheelchairs	116
6. Showering and Toilet Aids	118
7. Evacuation Products	128
8. Respose Specialist Seating	132
9. What Our Service Offers You	140
10. Mattress Management Guide	141

HERE TO HELP

Harvest Healthcare Contact Number: 01709 377172


Option 1 - PRODUCT SUPPLY / CUSTOMER SERVICES – sales@harvesthealthcare.co.uk

Option 2 - SERVICE & MAINTENANCE / REPAIR SCHEDULING – servicing@harvesthealthcare.co.uk

Option 3 - MATTRESS REPAIR DEPARTMENT – mattress.repairs@harvesthealthcare.co.uk

Option 4 - ACCOUNTS / PURCHASING – accounts@harvesthealthcare.co.uk

FOR MORE INFORMATION VISIT WWW.HARVESTHEALTHCARE.CO.UK

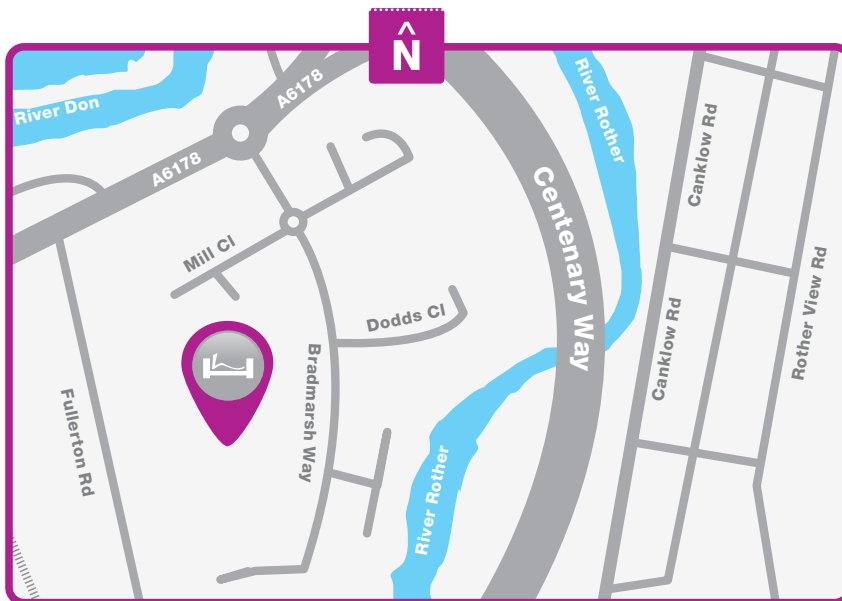

FIND US TODAY

Harvest Healthcare Ltd
Sheaf House, Bradmarsh Business
Park,
Bradmarsh Way,
Rotherham S60 1BW

T +44 (0)1709 377172

F +44 (0)1709 377173

E sales@harvesthealthcare.co.uk

www.harvesthealthcare.co.uk



www.harvesthealthcare.co.uk

Sheaf House, Bradmarsh Business Park, Bradmarsh Way, Rotherham S60 1BW
T +44 (0)1709 377172 **F** +44 (0)1709 377173 **E** sales@harvesthealthcare.co.uk

1. BEDS

Harvest Healthcare manufacture a wide range of profiling beds designed to provide support to a variety of service users. The exclusive range of Woburn beds, along with the Cadence range, offer ease of movement to both service user and carer, providing a higher quality of life. The profiling platform has an adjustable knee break, head section, and the entire bed frame is height adjustable. Safety and durability are key aspects of our design ethos – each bed has something different to offer.

In this section, you will find:

Woburn Community Profiling Beds (Community, Community Low, Community 1200)

Woburn Profiling Beds (Woburn, Woburn Low, Ultra Low and Ultra Low 800)

Woburn Ultimate Profiling Beds (Ultimate 1200, Ultimate 1400)

Cadence Profiling Beds (Comfort, Select, Premium)


Woburn Bed Accessories

Bumpers

We offer simple to use bedrail height measurement tools so that you can quickly and easily check if you are compliant with Side Rail Regulations 60601/2/52.

The regulations state that the height of the top edge of the side rail must be more than 220mm above the non-compressed mattress, and the gap between bedrails must be smaller than 120mm.

If you would like to be sent copies of the rail height measuring guide (**Product Code: RAILGUIDE**), please get in touch.

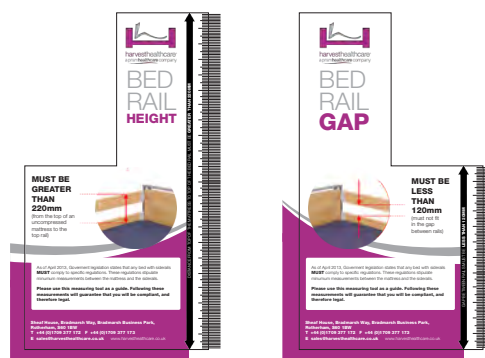

WOBURN COMMUNITY PROFILING BEDS



Harvest Healthcare Woburn Community Profiling beds are designed to fit within the community markets. These versatile beds are intended to enhance comfort, safety, and functionality, making them invaluable for individuals with diverse needs.

In a community setting, the Woburn Community Profiling beds provide essential support for elderly or disabled individuals, facilitating easier mobility, reducing the risk of pressure ulcers, and assisting caregivers in delivering proper care. The beds have integrated safety features like side rails and braking systems, ensuring patient safety and reducing the risk of falls.

All the Woburn beds are built to comply with Side Rail Regulations. This allows rails to be used on beds without compromising safety for the user. Identifying the needs of the service user before choosing a bed is of great importance. Harvest Healthcare have a knowledgeable team able to provide you with enough information to make such a decision.

WOBURN COMMUNITY BED RANGE

	Max. Height (mm)	Min. Height (mm)	Safe Working Load (kg/stone)	Max. Patient Weight (kg/stone)	Overall Length (mm)	Overall Width (mm)	Max. Mattress Thickness With Rails (mm)	Warranty (Years)	200mm Extension Available	In-Built Extension (Optional)	Autoregression	Trendelenburg
Woburn Community	800	380	220/34	185/29	2115	1015	205	5	✓	✓	✓	✓
Woburn Community Low	640	220	220/34	185/29	2115	1015	205	5	✓	✓	✓	✓
Woburn Community 1200	815	385	193/30	158/25	2120	1320	205	5	✓	∅	✓	✓

WOBURN COMMUNITY BED



4-section profiling bed designed to fit within the community markets, the Woburn Community bed is built with the same quality across the range - at a more cost-effective price.

Vinyl coated, the bed aids with cleaning and infection control. It is also suited for multiple dismantles across the community, part of the more robust design.



Part numbers and bed information labels are provided on each bed for ease of servicing and repair.

KEY FEATURES & BENEFITS

- Low voltage power supply (no mains voltage on the bed chassis to reduce the risk of electric shock)
- Electric backrest and knee break with mechanical lower leg support
- Locking handset
- Integral 115 mm side rails (will take up to 8 inch mattress) (optional)
- Inbuilt emergency battery backup
- Extendable by 200mm
- Autoregression (optional)
- Quick release actuators
- Trendelenburg as standard
- Replaceable PCB
- Replaceable actuator pistons
- Easy to install
- Padded headboards & other wood finish options are available on request

5 year warranty

185kg/29st
max user weight

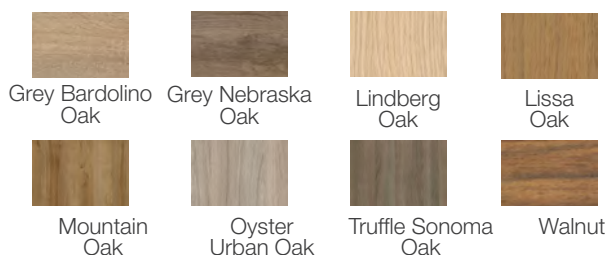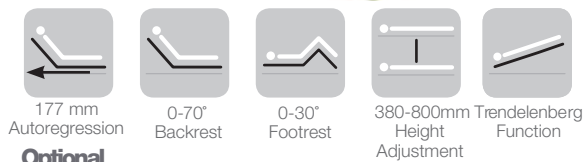

Technical Information

Mattress Platform Dimensions	2000(L)x900(W)mm
Adjustable Height	380-800mm
Safe Working Load (SWL)	220kg/34 Stone

Product Information

Product Code	Product Description
HLBC791.02.TVC	Woburn Community Bed
HLBC791.02.TVC.SR	Woburn Community Bed with Side Rails

Optional Extras

- Length Extension: **HLDWC123** (Includes Base & Side Rails)
 - Integral length extension is also available
 - Grab Handle: **HLB685**
 - Mattress Extension: **HFIM2016**
 - Standard Length Rails: **HLDWC110**
- If not purchased with the bed you may also need:
- Side Rail Fingers x 4: **HLDWC112**
 - End Caps x 8: **HLDWC001**
 - Bottom Stops x 4: **HLDWC122**
 - Handset Clip: **HLDW008**
 - Reusable Cable Retainer: **HLB603**

If bumpers are required please see pages (26-32)

WOBURN COMMUNITY LOW BED



In addition to the Woburn Community, the Woburn Community Low profiling bed features a lower entry and exit height. Ideal for many community environments, the bed will lower to 220mm from ground level, allowing for easy access to service users.

Vinyl coated, the bed aids with cleaning and infection control. The handset is lockable to prevent tampering by those other than carers.

Part numbers and bed information labels are provided on each bed for ease of servicing and repair.

KEY FEATURES & BENEFITS

- Low voltage power supply (no mains voltage on the bed chassis to reduce the risk of electric shock)
- Electric backrest and knee break with mechanical lower leg support
- Locking handset
- Integral 115 mm side rails (will take up to 8 inch mattress) (optional)
- Variable height (220 - 640mm)
- Inbuilt emergency battery backup
- Extendable by 200mm
- Trendelenburg as standard
- Quick release actuators
- Replaceable PCB
- Replaceable actuator pistons
- Easy to install
- Padded headboards & other wood finish options are available on request

5 year warranty

185kg/29st
max user weight



- 177 mm Autoregression
 - 0-70° Backrest
 - 0-30° Footrest
 - 220-640mm Height Adjustment
 - Trendelenburg Function
- Optional**



Technical Information

Mattress Platform Dimensions	2000(L)x900(W)mm
Adjustable Height	220-640mm
Safe Working Load (SWL)	220kg/34 Stone

Product Information

Product Code	Product Description
HLBC791L.TVC	Woburn Community Low Bed
HLBC791L.TVC.SR	Woburn Community Low Bed with Side Rails

Optional Extras

- Length Extension: **HLDWC123** (Includes Base & Side Rails)
- Integral length extension is also available
- Grab Handle: **HLB685**
- Mattress Extension: **HFIM2016**
- Standard Length Rails: **HLDWC110**

If not purchased with the bed you may also need:

- Side Rail Fingers x 4: **HLDWC112**
- End Caps x 8: **HLDWC001**
- Bottom Stops x 4: **HLDWC122**
- Handset Clip: **HLDW008**
- Reusable Cable Retainer: **HLB603**

If bumpers are required please see pages (26-32)

WOBURN COMMUNITY 1200 BED



The Woburn Community 1200 is Harvest's first wide profiling bed, built to the same high-quality as other beds in the range. Although not a bariatric bed, it has a safe working load of 30 Stone.

The generous width allows patients to move freely. The lockable handset allows for quick adjustments of height, backrest and knee break, and the attractive vinyl-clad head and footboards facilitate infection control. This bed is for one-person use only and is not for plus size patients.



The bed is designed for easy cleaning to aid infection control. Part numbers and bed information labels are provided on each bed for ease of servicing and repair.

KEY FEATURES & BENEFITS

- Low voltage power supply (no mains voltage on the bed chassis to reduce the risk of electric shock)
- Electric backrest and knee break with mechanical lower leg support
- Locking handset
- Integral 115 mm side rails (will take up to 8 inch mattress) (Included)
- Inbuilt emergency battery backup
- In-built autoregression
- Trendelenburg as standard
- Quick release actuators
- Replaceable PCB
- Replaceable actuator pistons
- Easy to install
- Also available in vinyl coated wooden head and footboards
- Padded headboards & other wood finish options are available on request

5 year warranty

158kg/25st max user weight

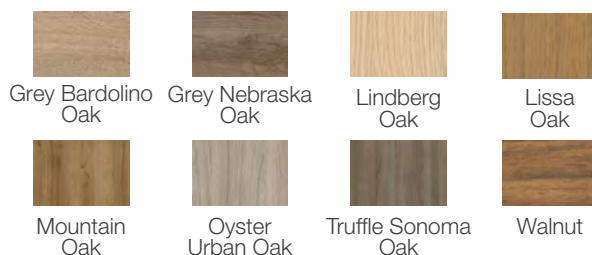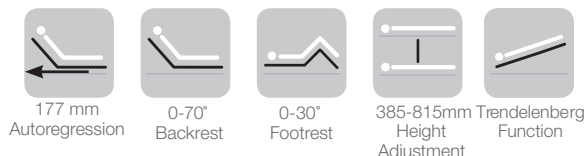

Technical Information

Mattress Platform Dimensions	2000(L)x1200(W)mm
Adjustable Height	385-815mm
Safe Working Load (SWL)	193kg/30st

Product Information

Product Code	Product Description
HLBC791.12AR.TVC	Woburn Community 1200
HLBC791.12AR.TVC.SR	Woburn Community 1200 with Side Rails

Optional Extras

- Length Extension Metal Base: [ENWIDE002](#)
- Length Extended Rails: [HLDWC111](#)
- Integral length extension is also available
- Grab Handle: [HLB685](#)
- Mattress Extension: [HFIM2016.12](#)
- Standard Length Rails: [HLDWC110](#)

If not purchased with the bed you may also need:

- Side Rail Fingers x 4: [HLDWC112](#)
- End Caps x 8: [HLDWC001](#)
- Bottom Stops x 4: [HLDWC122](#)
- Handset Clip: [HLDW008](#)
- Reusable Cable Retainer: [HLB603](#)

If bumpers are required please see pages (26-32)

WOBURN PROFILING BEDS



Harvest Healthcare Woburn Profiling beds are designed to fit within the long-term care markets. The beds significantly improve the quality of care and comfort for long-term care patients while supporting caregivers in delivering effective, safe, and efficient care.

The Woburn Profiling beds allow for easy adjustment of the bed's height, backrest, and leg sections, enabling personalised positioning that minimises the risk of pressure ulcers and promotes circulation. The ability to change positions also aids in respiratory and digestive functions. The beds have integrated safety features like side rails and braking systems, ensuring patient safety and reducing the risk of falls. For caregivers, these beds make patient care more manageable by enabling easier transfers, reducing strain on caregivers during repositioning, and facilitating daily care tasks.

A nationwide team of engineers are on hand to install, provide training and service profiling beds. At Harvest Healthcare, we provide a cost-effective, quality service with an honest approach – preferring to repair rather than replace wherever possible.

Our compliance checks ensure all equipment serviced by Harvest Healthcare is reliable and efficient, enabling staff to fully carry out their role effectively. Health and safety of service users and staff is important to us, so with ongoing support through regular training sessions and equipment audits, we will ensure that the quality of patient care increases and professional standards are maintained.

WOBURN BED RANGE

	Max. Height (mm)	Min. Height (mm)	Safe Working Load (kg/stone)	Max. Patient Weight (kg/stone)	Overall Length (mm)	Overall Width (mm)	Max. Mattress Thickness With Rails (mm)	Warranty (Years)	200mm Extension Available	In-Built Extension (Optional)	Autoregression	Trendelenburg
Woburn	795	365	220/34	185/29	2135	1025	205	5	✓	✓	✓	✓
Woburn Low	650	220	220/34	185/29	2135	1025	205	5	✓	✓	✓	✓
Woburn Ultra-Low	640	67	185/29	175/27	2130	900	205	5	✓	∅	✓	✓
Woburn Ultra-Low 800	800	90	185/29	175/27	2399	900	205	5	✓	∅	✓	✓

WOBURN BED



The quality features and reliability of the Woburn profiling bed provide outstanding comfort and ease of care for both occupant and carer.

The robust wooden side rails are collapsible below mattress level and extend higher to allow a safe mattress depth of 205mm, giving scope for improved pressure relief and comfort.

Backrest, knee break and platform height can all be adjusted at the touch of a button.

Part numbers and bed information labels are provided on each bed for ease of servicing and repair.



KEY FEATURES & BENEFITS

- Low voltage power supply (no mains voltage on the bed chassis to reduce the risk of electric shock)
- Electric backrest and knee break with mechanical lower leg support
- Locking handset
- Integral 115 mm side rails (will take up to 8 inch mattress) (Included)
- Inbuilt emergency battery backup
- Autoregression (optional)
- Trendelenburg (optional)
- Wraparound wooden head and footboards
- Quick release actuators
- Replaceable PCB
- Replaceable actuator pistons
- Easy to install
- Padded headboards & other wood finish options are available on request

5 year warranty



185kg/29st
max user weight

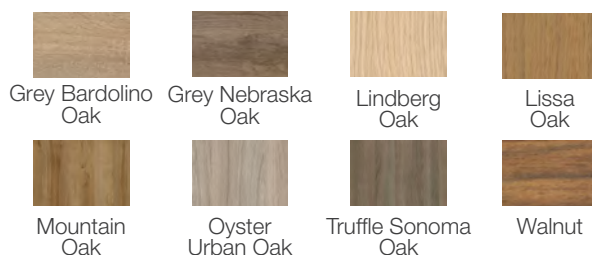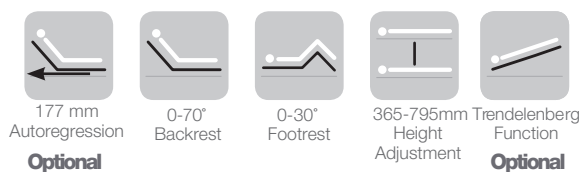

Technical Information

Mattress Platform Dimensions	2000(L)x900(W)mm
Adjustable Height	365-795mm
Safe Working Load (SWL)	220kg/34 Stone

Product Information

Product Code	Product Description
HLB791.03	Woburn Bed
HLB791.03.SR	Woburn Bed with Side Rails

Optional Extras

- Length Extension: **HLDWC123** (Includes Base & Side Rails)
- Integral length extension is also available
- Grab Handle: **HLB685**
- Mattress Extension: **HFIM2016**
- Standard Length Rails: **HLDWC110**

If not purchased with the bed you may also need:

- Side Rail Fingers x 4: **HLDWC112**
- End Caps x 8: **HLDWC001**
- Bottom Stops x 4: **HLDWC122.02**
- Reusable Cable Retainer: **HLB603**

If bumpers are required please see pages (26-32)

WOBURN LOW BED



harvesthealthcare
a prismhealthcare company

This dependable 4-section profiling bed features a high-quality control system, a reduced minimum platform height and a higher load capacity. By lowering the platform to only 220mm from floor level, the risk of injury for those with a tendency to roll out of bed is significantly reduced.

Part numbers and bed information labels are provided on each bed for ease of servicing and repair.

KEY FEATURES & BENEFITS

- Low voltage power supply (no mains voltage on the bed chassis to reduce the risk of electric shock)
- Electric backrest and knee break with mechanical lower leg support
- Locking handset
- Integral 115 mm side rails (will take up to 8 inch mattress) (optional)
- Inbuilt emergency battery backup
- Autoregression (optional)
- Trendelenburg (optional)
- Wraparound wooden head and footboards
- Quick release actuators
- Replaceable PCB
- Replaceable actuator pistons
- Easy to install
- Padded headboards & other wood finish options are available on request

5 year warranty



185kg/29st
max user weight



177 mm
Autoregression

Optional



0-70°
Backrest



0-30°
Footrest



220-650mm
Height
Adjustment



Trendelenburg
Function

Optional



Grey Bardolino
Oak



Grey Nebraska
Oak



Lindberg
Oak



Lissa
Oak



Mountain
Oak



Oyster
Urban Oak



Truffle Sonoma
Oak



Walnut

Technical Information

Mattress Platform Dimensions	2000(L)x900(W)mm
Adjustable Height	220-650mm
Safe Working Load (SWL)	220kg/34st

Product Information

Product Code	Product Description
HLB797.03	Woburn Low Bed
HLB797.03.SR	Woburn Low Bed with Side Rails

Optional Extras

- Length Extension: [HLDWC123](#) (Includes Base & Side Rails)
- Integral length extension is also available
- Grab Handle: [HLB685](#)
- Mattress Extension: [HFIM2016](#)
- Standard Length Rails: [HLDWC110](#)

If not purchased with the bed you may also need:

- Side Rail Fingers x 4: [HLDWC112](#)
- End Caps x 8: [HLDWC001](#)
- Bottom Stops x4: [HLDWC122.02](#)
- Reusable Cable Retainer: [HLB603](#)

If bumpers are required please see pages (26-32)

WOBURN ULTRA LOW BED



harvesthealthcare
a prismhealthcare company

The Woburn Ultra Low profiling bed is a unique system, which allows the platform to be lowered to as little as 67mm from ground level.

Intended for use without side rails, the service user can safely roll onto a crash mat, should they be inclined to roll whilst in bed.

Built with the Trendelenburg function as standard, service users are able to raise themselves to an almost stood or seated position, allowing residents to be more comfortable whilst in bed.

The setting of height, backrest or knee break is simple, using the lockable handset. Part numbers and bed information labels are provided on each bed for ease of servicing and repair.

KEY FEATURES & BENEFITS

- Low voltage power supply (no mains voltage on the bed chassis to reduce the risk of electric shock)
- Electric backrest and knee break with mechanical lower leg support
- Locking handset
- In-built emergency battery backup
- In-built autoregression
- Trendelenburg as standard
- Available with an enclosed headboard and footboard
- Wooden head and footboards
- Easy Cleaning - complies with MHRA guidelines
- Quick release actuators
- Replaceable PCB
- Replaceable actuator pistons
- Padded headboards & other wood finish options are available on request

5 year warranty



175kg/27st
max user weight



177 mm
Autoregression



0-70°
Backrest



0-30°
Footrest



67-640mm
Height
Adjustment



Trendelenburg
Function



Grey Bardolino
Oak



Grey Nebraska
Oak




Lindberg
Oak



Lissa
Oak



Mountain
Oak

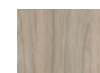

Oyster
Urban Oak



Truffle Sonoma
Oak



Walnut

Technical Information

Mattress Platform Dimensions	2000(L)x900(W)mm
Adjustable Height	67-640mm
Safe Working Load (SWL)	185kg/29st

Product Information

Product Code	Product Description
HLB799.03AR	Woburn Ultra Low Bed

Optional Extras

- Length Extension Metal Base: [HLDWUL113.02](#)
- Grab Handle: [HLDWUL685](#)
- Mattress Extension: [HFIM2016](#)
- Mesh Side Rails for Woburn Ultra Low: [HLDWUL110](#)
- Reusable Cable Retainer: [HLB603](#)

If fall mats are required see page 60

WOBURN ULTRA LOW 800 BED



A fully functioning four-section Ultra-Low profiling bed. Platform can be lowered to 90mm off the floor for patient safety and raised to 800mm to make it easier for the carer when providing care for the patient. Electronically setting of height, back rest or knee brake is simple, using the lockable handset. The Woburn Ultra-Low 800 has a Trendelenburg function as standard. 275-690mm


Optional bed levels, lighting and accessories available.

KEY FEATURES & BENEFITS

- Low voltage power supply (no mains voltage on the bed chassis to reduce the risk of electric shock)
- Electric backrest and knee break with mechanical lower leg support
- Locking handset
- Emergency battery backup
- In-built autoregression
- Trendelenburg as standard
- Wooden head and footboards
- Easy Cleaning - complies with MHRA guidelines
- Compatible with environmental control system
- Padded headboards & other wood finish options are available on request

5 year warranty

175kg/27st
max user weight

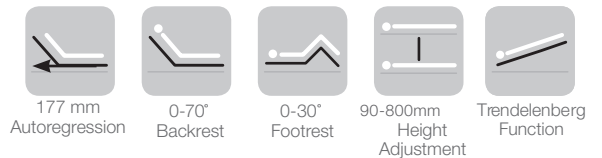

Technical Information

Mattress Platform Dimensions	2000(L)x900(W)mm
Adjustable Height	90-800mm
Safe Working Load (SWL)	185kg/29st

Product Information

Product Code	Product Description
HLB800.01AR	Woburn Ultra Low 800 Bed

Optional Extras

- Length Extension Metal Base: [HLDWUL113.03](#)
- Grab Handle: [HLDWUL685](#)
- Mattress Extension:
- Mesh Side Rails for Woburn Ultra Low 800: [HLDWUL8110](#)
- Reusable Cable Retainer: [HLB603](#)

If fall mats are required see page 60

WOBURN ULTIMATE PROFILING BEDS



Harvest Healthcare Woburn Ultimate Profiling beds are designed to fit for the plus size residents and service users. At Harvest Healthcare we aim to relieve the stress of sourcing and fitting plus sized equipment, leaving you to focus on providing the best quality care and support.

Our plus size range includes beds, mattresses, hoists and bathroom equipment, all designed to make life easier for residents and for their carers and give everyone access to the love and support they deserve.

Our knowledgeable customer services and sales teams can talk you through any tricky situations. Our sister company Prism Medical can also support with ceiling track hoists, an excellent way to make use of smaller spaces.

WOBURN ULTIMATE BED RANGE

	Max. Height (mm)	Min. Height (mm)	Safe Working Load (kg/stone)	Max. Patient Weight (kg/stone)	Overall Length (mm)	Overall Width (mm)	Max. Mattress Thickness With Rails (mm)	Warranty (Years)	200mm Extension Available	In-Built Extension (Optional)	Autoregression	Trendelenburg
Woburn Ultimate 1200	690	275	381/60	300/47	2195	1340	260	5	✓	∅	✓	✓
Woburn Ultimate 1400	690	275	381/60	300/47	2195	1540	260	5	✓	∅	✓	✓

WOBURN ULTIMATE BED



harvesthealthcare
a prismhealthcare company

The Woburn Ultimate Profiling Bed is the highest capacity bed within the Woburn range. Featuring triple aluminium side rails and under-bed strengthening, the Ultimate excels in patient support and ease of infection control.

Available in 1200mm and 1400mm wide, with a 60stone weight capacity, this bed is able to support most bariatric patients. Comfort adjustments are made with ease due to the powerful lifting components and the high-quality electronic control system.

Part numbers and bed information labels are provided on each bed for ease of servicing and repair.

KEY FEATURES & BENEFITS

- Low voltage power supply (no mains voltage on the bed chassis to reduce the risk of electric shock)
- Electric backrest and knee break with mechanical lower leg support
- Locking handset
- Aluminium side rails
- Wooden head and footboards
- Inbuilt emergency battery backup
- In-built autoregression
- Trendelenburg as standard
- In-built stability system
- Quick release actuators
- Replaceable PCB
- Replaceable actuator pistons
- Padded headboards & other wood finish options are available on request

5 year warranty



300kg/47st
max user weight

Technical Information

	HLB798.12	HLB798.14
Mattress Platform Dimensions	2000(L)x1200(W)mm	2000(L)x1400(W)mm
Adjustable Height	275-690mm	275-690mm
Safe Working Load (SWL)	381kg/60st	381kg/60st

Product Information

Product Code	Product Description
HLB798.12	Woburn Ultimate 1200 with Side Rails
HLB798.14	Woburn Ultimate 1400 with Side Rails



177 mm
Autoregression



0-70°
Backrest



0-30°
Footrest



275-690mm
Height
Adjustment



Trendelenburg
Function



Grey Bardolino
Oak



Grey Nebraska
Oak



Lindberg
Oak



Lissa
Oak



Mountain
Oak



Oyster
Urban Oak



Truffle Sonoma
Oak



Walnut

Optional Extras

- Length Extension Kit for Ultimate 1200 Wide Bed: [HLDWC123.12](#) (includes metal base & rails extensions)
- Mattress Extension 1200 Wide: [HFIM2016.12](#)
- Length Extension Kit for Ultimate 1400 Wide Bed: [HLDWC123.14](#) (includes metal base & rails extensions)
- Mattress Extension 1400 Wide: [HFIM2016.14](#)
- Grab Handle: [HLB686](#)
- Reusable Cable Retainer: [HLB603](#)

If bumpers are required please see pages (26-32)

CADENCE PROFILING BEDS



With cutting-edge technology and comfort at the forefront, our Cadence Beds are the ultimate solution for resident care.

Experience unparalleled comfort and style with our advanced features and state-of-the-art technology. But that's not all - the Cadence Bed Range is also fully customisable to your exact specifications. Choose from a wide range of bold colours and calming fabrics to make your bed the perfect addition to any care environment.

Design the perfect bed and elevate your resident's well-being and safety with the Harvest Healthcare Cadence Bed Range.



CADENCE BED RANGE

	Max. Height (mm)	Min. Height (mm)	Safe Working Load (kg/stone)	Max. Patient Weight (kg/stone)	Overall Length (mm)	Overall Width (mm)	Max. Mattress Thickness With Split Side Rails (mm)	Max. Mattress Thickness With Continuous Side Rails (mm)	Warranty (Years)	Alternative Finish Available	200mm Extension Available	Autoregression	Trendelenburg
Cadence Beds	800	240	235/37	200/31	2115	1050	150	205	5	✓	✓	✓	✓

CADENCE COMFORT BED



harvesthealthcare
a prismhealthcare company

The Cadence Comfort Bed is a stylish, versatile and cost-effective option that features a standard head and footboard with a curved design, complete with side skirts. The bed can be used with either continuous side rails or split side rails, which are available as optional extras, providing the flexibility to meet a variety of needs and preferences.

With its elegant design and adaptable features, the Cadence Comfort Bed is a comfortable and stylish choice for any bedroom.

A full Cadence Beds brochure is available on request for more details.

KEY FEATURES & BENEFITS

- Lockable handset
- Electric backrest and knee break with mechanical lower leg support
- Autoregression
- Trendelenburg
- Stylish wooden continuous and split side rails (optional)
- 4 braked castors as standard
- Central brake available on request
- Alternative wood finish available
- Full upholstered headboards are available on request

5 year warranty



200kg/31st
max user weight



Grey Bardolino
Oak



Grey Nebraska
Oak



Lindberg
Oak



Lissa
Oak



Mountain
Oak



Oyster
Urban Oak



Truffle Sonoma
Oak



Walnut

Product Information

Product Code	Mattress Platform Dimensions	Adjustable Height	Safe Working Load (SWL)
CADENCE.COMFORT.C	2000(L)x900(W)mm	240-800mm	235kg/37st

CADENCE SELECT BED



harvesthealthcare
a prismhealthcare company

The Cadence Select bed is a sleek and functional mid-range option that boasts a standard headboard and a low footboard, both adorned with side skirts for added style. The bed is designed to be used with split side rails, which are available as optional extras, providing additional support and customisation options to suit individual needs.

With its sophisticated design and versatile features, the Cadence Select bed is an ideal choice for those looking for a comfortable and practical bed frame that doesn't sacrifice style.

A full Cadence Beds brochure is available on request for more details.

KEY FEATURES & BENEFITS

- Lockable handset
- Electric backrest and knee break with mechanical lower leg support
- Autoregression
- Trendelenburg
- Stylish wooden continuous and split side rails (optional)
- 4 braked castors as standard
- Central brake available on request
- Alternative wood finish available
- Full upholstered headboards are available on request

5 year warranty



200kg/31st
max user weight



Grey Bardolino
Oak



Grey Nebraska
Oak



Lindberg
Oak



Lissa
Oak



Mountain
Oak



Oyster
Urban Oak



Truffle Sonoma
Oak



Walnut

Product Information

Product Code	Mattress Platform Dimensions	Adjustable Height	Safe Working Load (SWL)
CADENCE.SELECT.LFB	2000(L)x900(W)mm	240-800mm	235kg/37st

CADENCE PREMIUM BED



harvesthealthcare
a prismhealthcare company

The Cadence Premium bed is a luxurious and elegant bed frame that boasts a curved headboard and low curved footboard, both accented with stylish side skirts. The bed is designed to be used with split side rails, which can be added as optional extras, providing enhanced support and customisation options to meet individual needs.

With its refined design and customisable features, the Cadence Premium bed is a perfect choice for those seeking a high-end and personalised sleeping experience.

A full Cadence Beds brochure is available on request for more details.

KEY FEATURES & BENEFITS

- Lockable handset
- Electric backrest and knee break with mechanical lower leg support
- Autoregression
- Trendelenburg
- Stylish wooden continuous and split side rails (optional)
- 4 braked castors as standard
- Central brake available on request
- Alternative wood finish available
- Full upholstered headboards are available on request

5 year warranty



200kg/31st
max user weight



Grey Bardolino
Oak



Grey Nebraska
Oak



Lindberg
Oak



Lissa
Oak



Mountain
Oak



Oyster
Urban Oak



Truffle Sonoma
Oak




Walnut

Product Information

Product Code	Mattress Platform Dimensions	Adjustable Height	Safe Working Load (SWL)
CADENCE.PREMIUM	2000(L)x900(W)mm	240-800mm	235kg/37st

BED ACCESSORIES

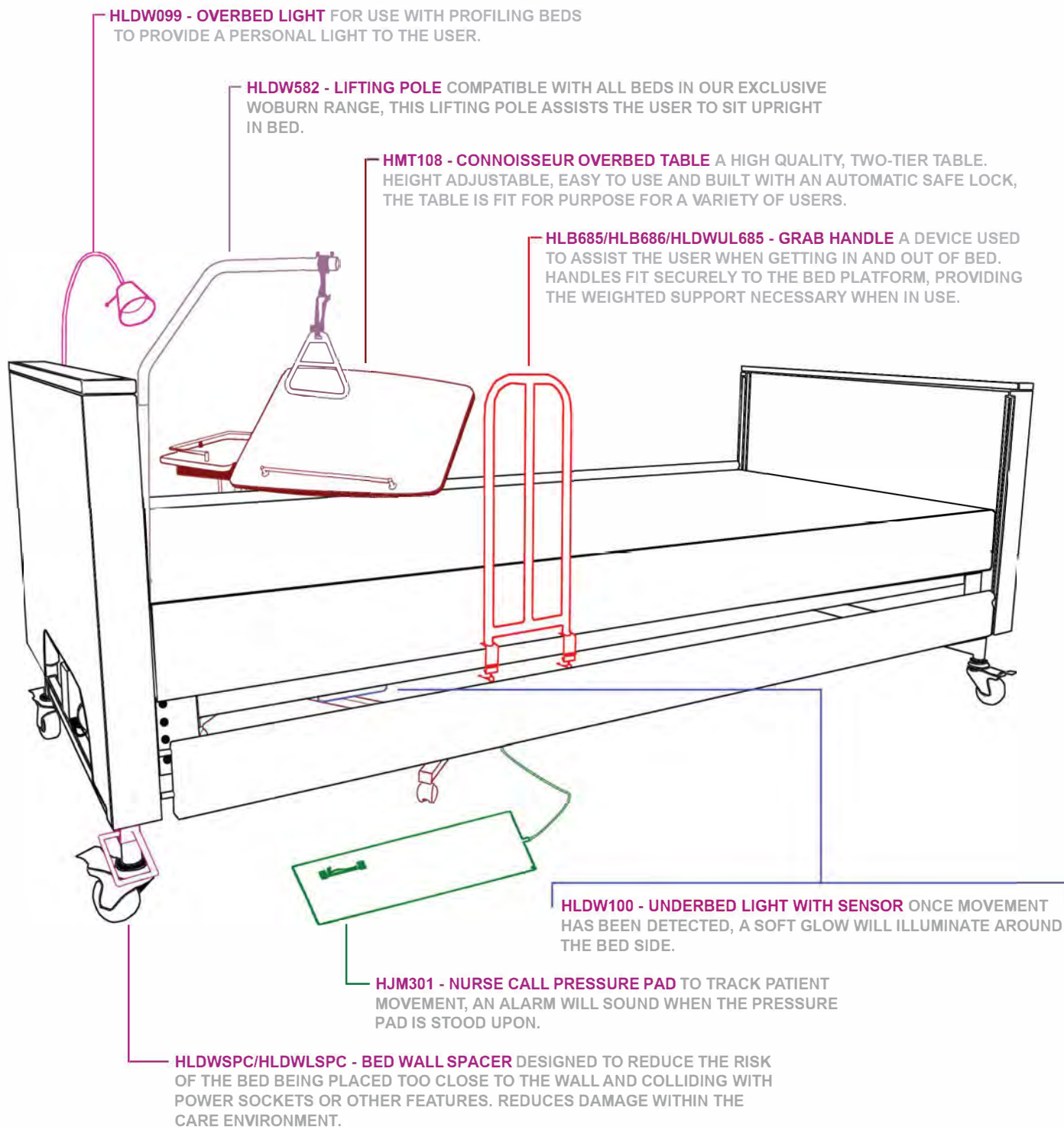

HLDU582/HLDU582.14 - ULTIMATE LIFTING POLE COMPATIBLE WITH OUR EXCLUSIVE WOBURN ULTIMATE BED, THIS ULTIMATE LIFTING POLE ASSISTS THE USER TO SIT UPRIGHT IN BED.

HLDW300/HLDU300 - BED HOOKS FOR USE WITH ANY PUMPS WHICH POWER AN ACTIVE MATTRESS, THE HOOKS ATTACH TO THE BED END, AND CAN THEN HOLD THE PUMPS.

HJR075 - BED ALARM SYSTEM CONTROL BOX IS DESIGNED FOR USE WITHIN A BED, NOTIFYING CARERS OF PATIENT MOVEMENT. THIS HELPS REDUCE THE LIKELIHOOD OF FALLS.

HLDW / HLDWC / HLD - BED LENGTH EXTENSION KIT IS USED TO EXTEND YOUR PROFILING BED BY AN ADDITIONAL 200MM. SOME USERS ARE NOT COMFORTABLE IN A TYPICAL 2000MM LONG PROFILING BED, WHICH CAN CAUSE ENTRAPMENT ISSUES



IN-BUILT EXTENSION

To remove the need for removable platforms, Harvest has introduced the In-Built Platform Extension.

These add an additional 200mm to the Woburn beds only when needed. The modern beds are now more adaptable for a wider range of residents.



AUTOREGRESSION

AutoRegression Backrests provide patients with a more comfortable position.


The mechanisms steadily raise or lower the patient whilst controlling posture - this reduces the 'push-down' effect on the body, keeping the resident more central in the bed and removing shear and friction forces upon the body.

HOW CAN AUTOREGRESSION HELP PREVENT PRESSURE INJURY?

A standard backrest pushes the user down the bed as it raises, creating shear forces and increasing the risk of tissue damage.

Autoregression backrests move slowly backwards at the sacral area as they raise. This allows the bed to move fully around the user, rather than the user being moved to accommodate the bed.

This means that the sacral area and heels do not move at all as the bed profiles. This removes the risk of shear and friction, protecting against tissue damage, prevents the need for repositioning and maintains user comfort.

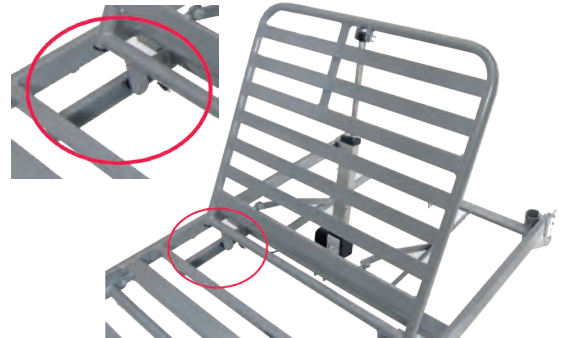

TRENDELENBURG FUNCTION

The Trendelenburg Function is available across the Woburn Bed Range. This function enables the carer to raise the foot end of the profiling bed higher than the head end, beneficial for use with respiratory patients. Anti-Trendelenburg enables the head end to be elevated higher than the feet, improving posture, allowing patients to feel as though they are standing or are seated whilst in bed and can increase blood flow.



Trendelenburg Handset



BED ACCESSORIES



harvesthealthcare
a prismhealthcare company

GRAB HANDLE FOR WOBURN BED

The Grab Handles assist the patient to sit up and get in and out of bed. They fit securely to the bed frame - no additional fixings required. They have a durable epoxy coating and can be used without removing the side rails. They are available to be used with Woburn, Woburn Ultimate and Woburn Ultra-Low beds.

Product Information

Product Code	Description	Safe Working Load (SWL)
HLB685	Grab Handle for Woburn Standard	60kg/9st
HLB686	Grab Handle for Woburn Ultimate	60kg/9st
HLDWUL685	Grab Handle for Woburn Ultra Low	60kg/9st

Dimensions: 240mm(W)x610mm(H) from platform



STANDARD LIFTING POLE

The Lifting Pole is made to be compatible with the Woburn Standard bed. It slots securely into sockets provided in the bed frame and has a durable powder coated finish. It assists the occupant to sit up in bed by holding onto the triangular handle.

Product Information

Product Code	Description	Safe Working Load (SWL)
HLDW582	Standard Lifting Pole with Handle	80kg/12st
HLB680	Handle	N/A



LIFTING POLE - ULTIMATE

The Ultimate 1200 Lifting Pole is made to be compatible with the Woburn Ultimate 1200 bed. It slots securely into sockets provided in the bed frame and has a durable powder coated finish. It assists the occupant to sit up in bed by holding onto the triangular handle.

Product Information

Product Code	Description	Safe Working Load (SWL)
HLDU582	Ultimate 1200 Lifting Pole with Handle	80kg/12st
HLDU582.14	Ultimate 1400 Lifting Pole with Handle	80kg/12st
HLB680	Handle	N/A





OVER-BED READING LIGHT

The Over-Bed Light is compatible with the full Woburn bed range. It is safe and simply connects to the low voltage power supply of the bed. It comes complete with an adjustable switch position and manoeuvrable head section. The light provides a valuable touch of home comfort, providing convenient directed lighting.



Product Information

Product Code	Description
HLDW099	Overbed Light

UNDER-BED LIGHT WITH MOVEMENT SENSOR

The Under-Bed Light uses an energy efficient LED light, compatible with the full Woburn bed range. It simply connects to the low voltage bed supply. It is safe and maintains a low temperature. The built-in sensor detects movement when the user leaves the bed and illuminates underneath and around the bed, helping to prevent trips and falls.



Product Information

Product Code	Description
HLDW100	Under Bed Light with Sensor for Woburn Beds

BED WALL SPACER

These spacers are designed to reduce the risk of the bed being placed too close to the wall and colliding with power sockets or other features. This can be costly to fix and often requires an emergency callout; the Bed Wall Spacer can eliminate these costs.

It increases safety by reducing the risk of electric shock due to socket damage. The Bed Wall Spacer also protects the bed chassis from damage if being operated against window sills/skirting boards or other features, which can be expensive to fix.

Product Information

Product Code	Description
HLDWLSPC	Bed Wall Spacer (Pair) for the Woburn & Woburn Low Beds



BED ACCESSORIES



BED LENGTH EXTENSION

The Bed Length Extension Platform is available to extend your profiling bed by an additional 200mm. Some users are not comfortable in a typical 2000mm long Profiling Bed, which can cause entrapment issues.

Product Information

Product Code	Dimensions
HLDWC110	Full Length Wooden Side Rails for Woburn Beds (set of 4)
HLDWC111	Extended Wooden Side Rails (set of 4)
HLDW113	Bed Length Extension Section 900mm(W)x200mm(L)
HLDW123	Woburn Side Rails Length Extension Pack 900mm Wide Bed
HLDWC123	Side Rail Length Extension Pack 115mm Height Rails
HLDU123.12	Ultimate 1200 Bed Length Extension Kit
HLDU123.14	Ultimate 1400 Bed Length Extension Kit



MESH SIDE RAILS FOR WOBURN ULTRA LOW AND ULTRA LOW 800

The Mesh Side Rail System is intended to be used to provide security to the occupant of the bed whilst the bed is raised to a correct height for giving care. This is to reduce the risk of a fall from the bed or potential entrapment from a side rail system. It can be used with the enclosed head and footboards with some adjustments.

Product Information

Product Code	Description
HLDWUL110	Full Length Mesh Side Rails for Woburn Ultra Low Bed
HLDWUL8110	Full Length Mesh Side Rails for Woburn Ultra Low 800 Bed



MATTRESS LENGTH EXTENSION

The Static Mattress Length Extensions have been designed for use with our Static Mattress range.

Ideal for use with existing mattresses, the piece will add a further 200mm to the length of the mattress and close any gaps between the mattress and bed frame. It is available in both 900mm and 1200mm widths.

Once in place, the kit continues to provide pressure relief around the heels and zonal areas.

Product Information

Product Code	Dimensions
HFIM2013	200(L) x 900(W) x 130(H)mm
HFIM2016	200(L) x 900(W) x 160(H)mm
HFIM2016.12	200(L) x 1200(W) x 160(H)mm
HFIM2020	200(L) x 900(W) x 200(H)mm



CONNOISSEUR OVER BED TABLE



harvesthealthcare
a prismhealthcare company

The Connoisseur Over-bed Table is a high-quality, two tier table. Height adjustable, easy to use and built with an automatic safety lock, the table is fit for purpose for a variety of users.

2 vinyl wrapped trays allow for more usability and provide more space. Available in Oak and Walnut options.

It can be used for both Woburn and Cadence range profiling beds.

KEY FEATURES & BENEFITS

- Height adjustable, hand-operated steel column
- Easy-to-use height control handle
- Automatic safe lock
- Available in oak and walnut and more
- U-shaped base
- 2 vinyl wrapped trays



Oak



Walnut

Product Information

Product Code	Dimensions	Adjustable Height	Castors Size
HMT107W (Walnut) HMT107OK (Oak)	440mm x 890mm (Overall Surface Area) 830mm x 410mm (U-shaped base) 40mm x 40mm (Tube Steel)	760mm-1150mm	4x50mm

STANDARD BUMPER ANTI-ENTRAPMENT BUMPER



STANDARD BUMPER WITH NET/MESH INSERT

The Standard Bumper with net/mesh inserts is designed to fit over side rails, providing safe protection for users in bed. The net/mesh inserts are entwined with the bumpers, to improve visibility between each side rail, reducing the feeling of entrapment. The bumper fits all Woburn beds. All the bumpers come with covers that are removable and machine washable, for better infection control.



Product Information

Product Code	Description
HCSBUMP112	Standard Bumper with Net/Mesh Inserts for Woburn, Community & 1200 Bed (L: 2000 mm)
HCSBUMP212	Extended Bumper with net/mesh inserts for Woburn, Community & 1200 Bed (L: 2200 mm)
HCSBUMP130	Standard Bumper with Net/Mesh Inserts for Woburn Ultimate Bed (L: 2000 mm)

ANTI-ENTRAPMENT BUMPER WITH NET/MESH INSERT

The Anti-Entrapment Bumper is a 2 bar bumper with net/mesh insert and anti-entrapment system with link sheet. Designed for users who move excessively in bed, to prevent limbs from becoming trapped. The Anti-Entrapment Bumper fits all standard width Woburn beds, whilst also being suitable for most other manufacturer's beds with two full length wooden side rails at each side of the bed. We also offer the Anti-Entrapment Bumper for beds with 3 bar side rails (Ultimate/Safeguard).

All the bumpers come with covers that are removable and machine washable, for better infection control.



Product Information

Product Code	Description
HCSAEBUMP791	Anti-Entrapment Bumper with Net/Mesh Inserts for Woburn & Community Bed
HCSAEBUMP140	Anti-Entrapment Bumper with Net/Mesh Inserts for Beds with 3 Bar Side Rails (Ultimate/Safeguard)

EXTENDED STANDARD BUMPER

EXTENDED ANTI-ENTRAPMENT BUMPER



EXTENDED LENGTH STANDARD BUMPER WITH NET/MESH INSERT

The 200mm Length Extension suitable for the Standard Bumper with net/mesh inserts is designed to fit over side rails, providing safe protection for users in bed. Created for use when the bed has been extended, ensuring full protection for the user. All the bumpers come with covers that are removable and machine washable, for better infection control.



Product Information

Product Code	Description
HCSBUMP112 + HCSBUMP112.EX20	Extended Standard Bumper with Net/Mesh Inserts for Woburn, Community & 1200 Bed
HCSBUMP130 + HCSBUMP130.EX20	Extended Standard Bumper with Net/Mesh Inserts for Woburn Ultimate Bed

EXTENDED LENGTH ANTI-ENTRAPMENT BUMPER WITH NET/MESH INSERT

The 200mm Length Extension is suitable for the Anti-Entrapment Bumper with net/mesh inserts and anti-entrapment system with link sheet. Designed for users who move excessively in bed, to prevent limbs from becoming trapped. Created for use when the bed has been extended, ensuring full protection for the user. All the bumpers come with covers that are removable and machine washable, for better infection control.



Product Information

Product Code	Description
HCSAEBUMP791 + HCSAEBUMP791.EX20	Extended Anti-Entrapment Bumper with Net/Mesh Inserts for Woburn & Community Bed
HCSAEBUMP140 + HCSAEBUMP140.EX20	Extended Anti-Entrapment Bumper with Net/Mesh Inserts for Woburn Ultimate Bed

CONTACT PAD HEAD AND FOOT BUMPER



CONTACT PAD FOR ANTI-ENTRAPMENT BUMPER

Designed for use with the Harvest Anti-Entrapment Bumper, the contact pads protect users with involuntary movement from skin damage and prevents the risk of friction injury from contact with the mesh inserts. The pads can be positioned anywhere inside the bed with a Hook and Loop fastening. All the bumpers come with covers that are removable and machine washable, for better infection control.



Product Information

Product Code	Description
HCS140CONTP	Contact Pad for Anti-Entrapment Bumper

HEAD AND FOOTBOARD BUMPER

The Head and Footboard bumper fits all Woburn and Community Beds. Providing extra protection for the user at the head and foot end of the bed. All the bumpers come with covers that are removable and machine washable, for better infection control.



Product Information

Product Code	Description
HBCWOB	Head and Footboard Bumper for Woburn, Woburn Low and Woburn Community Beds
HBCCOM1.12	Head and Footboard Bumper for Woburn Community 1200 Bed
HBCULT	Head and Footboard Bumper for Woburn Ultimate Beds

SPECIALIST ANTI-ENTRAPMENT BUMPER



SPECIALIST ANTI-ENTRAPMENT BUMPER


Adults of small stature or juniors are at increased risk of entrapment in a standard adult bed. This innovative solution quickly transforms a Community Woburn bed, making it suitable for smaller occupants, and meets junior bed regulations when paired with a locking handset box. Please find the installation video on our website. All the bumpers come with covers that are removable and machine washable, for better infection control.



Product Information

Product Code	Description
HJB002	Specialist Anti-Entrapment Bumper for Community Woburn Bed

Please click on the image below to see how to install our Specialist Anti-Entrapment Bumper. Alternatively, please go to Harvest Healthcare's YouTube page and find the Specialist Anti-Entrapment Bumper video. <https://tinyurl.com/Specialist-Bumper>

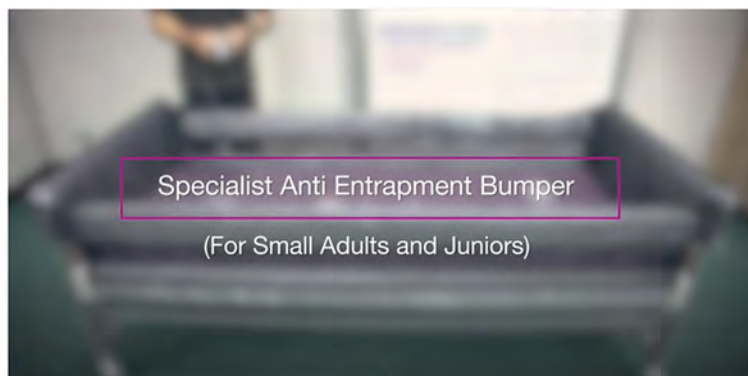

WRAPOVER BUMPER STANDARD AND HEIGHT EXTENSION



WRAPOVER BUMPER

The Wrapover Bumper is used to surround the bed side rails, providing a safe and cushioned barrier for patients. Designed with users in mind who are inclined to roll whilst in bed, meaning extra protection may be necessary. The three quarter length bumper wraps over the side rails and secures with Hook and Loop at the bottom, ensuring ease of use. All the bumpers come with covers that are removable and machine washable, for better infection control.



WRAPOVER BUMPER STANDARD & HEIGHT EXTENSION (200MM)

The 200mm height extension suitable for the Wrapover Bumper is used to surround the bed side rails, providing a safe and cushioned barrier for patients. Designed for use when the height attachment side rails are placed onto the bed, ensuring full protection for the user.

WRAPOVER BUMPER STANDARD & LENGTH EXTENSION (500MM)

The 500mm length extension suitable for the Wrapover Bumper is used to surround the bed side rails, providing a safe and cushioned barrier for patients. Designed for use when the bed has been extended, ensuring full protection for the user.



Product Information

Product Code	Description
HCSWRAP113	Harvest Wrapover Bumper Standard
HCSWRAP113 + HCSWRAP113.EX20HEIGHT	Harvest Wrapover Bumper Standard & Height Extension (200mm)
HCSWRAP113 + HCSWRAP113.EX50	Harvest Wrapover Bumper Standard & Length Extension (500mm)

UNIVERSAL SAFE ROLL BUMPER



UNIVERSAL SAFE ROLL

The Universal Safe Roll system is a state-of-the-art crash mat, designed to dramatically change the landscape of falls protection and support the reduction of the use of side rails. Traditional static floor mats are often a trip hazard, have to be moved to accommodate moving and handling equipment, and are frequently too thin to provide true protection from falls. This radial new design eliminates all these issues, providing safety for the resident and exceptional ease of use for carers.

The folding mat attaches to the side of the bed, and Hook and Loop up when the bed is raised for care (image 1), enabling access for carers and for equipment. As the bed is lowered the deep foam mat folds out to lie alongside the bed, as a true safety feature (image 2). Designed to fit our Woburn Ultra Low, Low and standard beds, the Universal Safe Roll is suitable for a wide range of users and is ideal for use within nursing homes and the domestic environment.

All the bumpers come with covers that are removable and machine washable, for better infection control.



Image 1



Image 2


Product Information

Product Code	Description
HFB003	Universal Safe Roll Bumper

Please click on the image below to see how to install our Universal Safe Roll.

Alternatively, please go to Harvest Healthcare's YouTube page and find the Universal Safe Roll video.

<https://tinyurl.com/Universal-Safe-Roll>

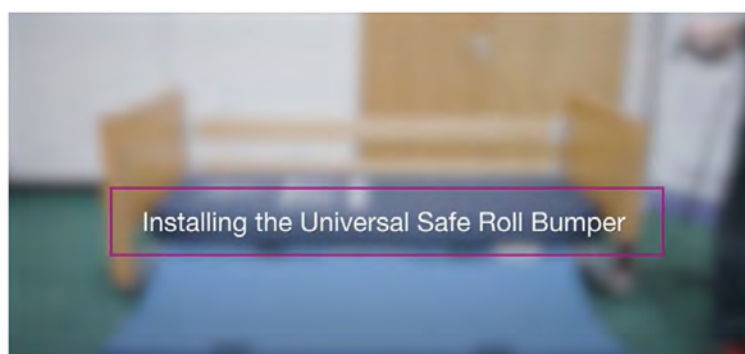

SNUG WALL BUMPER / BED SIDE PLATFORM BUMPER / SAFETY WEDGES



SNUG WALL BUMPER

The Snug Wall Bumper has been designed to fit onto our **Woburn bed range** or similar types of profiling beds produced by other manufacturers. Beds often need to be placed slightly away from the wall due to room fittings and fixtures, which can leave a gap. This can be an entrapment risk, and the Snug Bumper unobtrusively and comfortably fills this space. This removes the entrapment risk without the need for side rails to be fitted. The easy to fit bumper helps to prevent accidental contact with the wall and is created using soft to the touch material and spacer mesh to prevent the risk of suffocation and increase patient safety. All the bumpers come with covers that are removable and machine washable, for better infection control.

Please click on the image to see how to install the bumper.

Alternatively, please go to Harvest Healthcare's YouTube page and find the Snug Wall Bumper video.
<https://tinyurl.com/Snug-Wall>



Product Information

Product Code	Description
HWB003	Snug Wall Bumper

BED SIDE PLATFORM BUMPER

The Bed Side Platform Bumper has been designed to easily fit onto the outer edge of the bed frame to reduce the risk of injury for patients on the frame when getting in and out of bed. Suitable for use on all our Woburn bed range and similar types of profiling beds produced by other manufacturers. Created using soft to the touch material, protecting the user and increasing patient safety.

All the bumpers come with covers that are removable and machine washable, for better infection control.



Product Information

Product Code	Description
HWB004	Bed Side Platform Bumper

SAFETY WEDGE

The Safety Bed Wedges provide a safe, firm, sloping barrier to keep the user comfortably positioned whilst in bed. The wedges are connected to a waterproof base sheet meaning the product is easy to fit and will be secure once fitted onto the bed. All the bumpers come with covers that are removable and machine washable, for better infection control.



Product Information

Product Code	Description
HCS096	Safety Wedges with Double Duty Hook and Loop Base Sheet
HCS096B	Safety Wedges 900x900mm with Hook and Loop Loop

2. MATTRESSES AND CUSHIONS

Harvest Healthcare have a wide range of active and static mattresses and cushions, designed with comfort, pressure relief and weight distribution in mind to meet your needs.

Varying levels of protection are available for static mattresses; it all depends on the foam chosen and where the profile sections are cut.

To make the choice easier for active mattresses, we have created a mattress selector based on maximum user weight, mattress width and user risk. This selector can be seen below.

We also offer hybrid mattresses to give you the ability to step up or step down the pressure relieving abilities of the mattress.

In this section you will find:

Static Mattresses

Hybrid Mattresses

Active Mattresses

Pumps

Fall Mats, Underlays and Overlays

Cushions

Active Mattress Selector

Width	900mm				1200mm	1400mm	
Maximum User Weight	25 Stone 158kg	26 Stone 165kg	29 Stone 184kg	39 Stone 247kg	40 Stone 254kg	50 Stone 317kg	50 Stone 317kg
High Risk	Kensington Overlay with Harvest 1 Pump Sandringham Replacement with Harvest 1 Pump		Kensington Overlay with Harvest 2 Pump Sandringham Replacement with Harvest 2 Pump	<i>Prime Comfort Active Replacement without pump</i>			
Very High Risk	Hampton Replacement with Harvest 1 Pump	Duke Replacement with Harvest 2 Pump	Hampton Replacement with Harvest 2 Pump Windsor Replacement with Harvest 2 Pump Balmoral Replacement with Harvest 2 Pump	<i>Prime Comfort Active Replacement with pump</i> Duke Replacement with Harvest 3 Pump	Hampton Replacement with Harvest 3 Pump Windsor Replacement with Harvest 3 Pump Balmoral Replacement with Harvest 3 Pump	Duke Extra 2.0 Replacement with Harvest 3 Pump Hampton Extra Replacement with Harvest 3 Pump	

COMMUNITY STATIC MATTRESS



The Community static mattress is a cost-effective comfort mattress for use on profiling beds. The foam is highly resilient and fire retardant.

The Community mattress provides pressure distribution and can be turned to extend service life.

The mattress is available with castellated foam on request.

KEY FEATURES & BENEFITS

- Cost-effective
- Highly Resilient foam
- Four-way stretch, waterproof, anti-bacterial zipped cover with a crib 5 fire rating
- Shaped to allow profiling
- Two-way turn to extend life of the mattress
- Additional sizes available upon request
- Medium risk

3 year warranty

159kg/25st
max user weight



Product Information

Product Code	Description	Mattress Dimensions	Mattress Weight
HSMCOM	Community 900mm Wide	2000(L)x900(W)x125(H)mm	8kg
HSMCOM.12	Community 1200mm Wide	2000(L)x1200(W)x125(H)mm	10kg

COMMUNITY PLUS STATIC MATTRESS



The Community Plus static mattress is built on the Community mattress with higher quality but remaining cost-effective. The foam is thicker offering better pressure distribution and support but keeping the same level of comfort.

The mattress is available with castellated foam on request.

KEY FEATURES & BENEFITS

- Cost-effective
- Highly Resilient foam
- Four-way stretch, water resistant, anti-bacterial zipped cover with a crib 5 fire rating
- Shaped to allow profiling
- Two-way turn to extend life of the mattress
- Additional sizes available upon request
- High risk

3 year warranty



190kg/30st
max user weight



LOW MED HIGH V HIGH

PRESSURE RELIEF RISK

Product Information

Product Code	Description	Mattress Dimensions	Mattress Weight
HSMCOMP	Community Plus 900mm Wide	2000(L)x900(W)x150(H)mm	8.5kg
HSMCOMP.12	Community Plus 1200mm Wide	2000(L)x1200(W)x150(H)mm	11kg

PRIME COMFORT STATIC MATTRESS



The Prime Comfort static mattress is made of high resiliency foam offering excellent pressure relief.

The castellated foam aids comfort and facilitates air movement inside the mattress, reducing heat build up.

KEY FEATURES & BENEFITS

- Excellent comfort and support
- Four-way stretch, water resistant, anti-bacterial PU cover with welded seams and a crib 5 fire rating
- Additional sizes available upon request
- High risk

3 year warranty

177kg/28st
max user weight



Product Information

Product Code	Description	Mattress Dimensions	Mattress Weight
HSMPRIM	Prime Comfort 900mm Wide	2000(L)x900(W)x150(H)mm	9kg
HSMPRIM.12	Prime Comfort 1200mm Wide	2000(L)x1200(W)x150(H)mm	12kg

PRIME COMFORT PLUS STATIC MATTRESS



The Prime Comfort Plus static mattress provides excellent patient comfort and support.

A castellated foam core facilitates pressure redistribution and air movement inside the mattress, reducing heat build up. It has a firm foam side wall encasing a resilient central insert, enhancing stability and assisting patient transfer.

Slots in the base enable the mattress to profile.

KEY FEATURES & BENEFITS

- Excellent comfort and support
- Slots in the base enable profiling
- Four-way stretch, water resistant, anti-bacterial PU cover with welded seams and a crib 5 fire rating
- Additional sizes available upon request
- High risk

3 year warranty

247kg/39st
max user weight



Product Information

Product Code	Description	Mattress Dimensions	Mattress Weight
HSMPRIMP	Prime Comfort Plus 900mm Wide	2000(L)x900(W)x140(H)mm	13kg
HSMPRIMP.12	Prime Comfort Plus 1200mm Wide	2000(L)x1200(W)x140(H)mm	17kg

REFLECT 2 STATIC MATTRESS



The Reflect 2 mattress is a profiling mattress with a castellated viscoelastic top and a firm base. The zonal castellations give tailored pressure relief, and allow airflow through the mattress, improving the skin microclimate.

The castellated viscoelastic foam offers improved immersion and pressure distribution. The firm base gives support and helps with even distribution of weight. There are keyhole slots in the base to aid profiling.

This comfortable pressure care solution is a dependable staple for plus sized pressure care.

KEY FEATURES & BENEFITS

- Castellated viscoelastic top and firm foam base
- Slots in the base enable profiling
- Four-way stretch, water resistant, anti-bacterial PU cover with welded seams and a crib 5 fire rating
- Available in 2 different weight limits
- High risk

3 year warranty



Product Information

Product Code	Description	Mattress Dimensions	Max. Patient Weight
HFR147	Reflect 2, 900mm Wide	2000(L)x900(W)x165(H)mm	47kg/39st
HFR147.12	Reflect 2, 1200mm Wide	2000(L)x1200(W)x165(H)mm	247kg/39st
HFR147.1370	Reflect 2, 1370mm Wide	2000(L)x1370(W)x165(H)mm	247kg/39st
HFR227	Reflect 2 Plus Size 45st 1200mm Wide	2000(L)x1200(W)x180(H)mm	285kg/45st
HFR227.14	Reflect 2 Plus Size 45st 1400mm Wide	2000(L)x1400(W)x180(H)mm	285kg/45st
HFR227.60	Reflect 2 Plus Size 60st 1200mm Wide	2000(L)x1200(W)x180(H)mm	381kg/60st
HFR227.1400.60	Reflect 2 Plus Size 60st 1400mm Wide	2000(L)x1400(W)x180(H)mm	381kg/60st

FIRE-TECH 7 COMFORT STATIC MATTRESS



The Fire-Tech 7 Comfort is a pressure care static mattress for environments at increased risk of fire. The mattress has a Crib 7 fire rating which is the most demanding fire test in the industry.

It provides excellent pressure redistribution and comfort for patients at high risk of developing pressure related injury.

KEY FEATURES & BENEFITS

- Excellent comfort and support
- Four-way stretch, waterproof, anti-bacterial PU cover
- Available welded or zipped
- Crib 7 fire rating
- Additional sizes available upon request
- High risk

3 year warranty



190kg/30st
max user weight



LOW MED HIGH V HIGH

PRESSURE RELIEF RISK

Product Information

Product Code	Description	Mattress Dimensions	Mattress Weight
HSMFIRET.CFT	Fire-Tech 7 Comfort 900mm Wide	2000(L)x900(W)x150(H)mm	8.5kg
HSMFIRET.CFT.W	Fire-Tech 7 Comfort Welded Cover	2000(L)x900(W)x150(H)mm	8.5kg
HSMFIRET.CFT.12	Fire-Tech 7 Comfort 1200mm Wide	2000(L)x1200(W)x150(H)mm	11kg

FIRE-TECH 7 STATIC MATTRESS



The Fire-Tech 7 Static Mattress is a vandal proof mattress for challenging environments at increased risk of fire.

KEY FEATURES & BENEFITS

- Block foam in the base
- Anti-bacterial PU cover with welded seams
- Crib 7 fire rating
- Additional sizes available upon request



Product Information

Product Code	Description	Mattress Dimensions	Mattress Weight
HSMFIRE7.VP	Fire-Tech 7 900mm Wide Vandal Proof	2000(L)x900(W)x150(H)mm	8.5kg

DIVAN MATTRESS



The Divan mattress is an affordable, comfortable option for those in need of a non-profiling mattress.

Covered with a cotton, woven material, the mattress also has a PVC, water resistant coating which allows the mattress to breathe. This is important, as the foam within will last longer and regain shape after use.

The Divan mattress includes:

- 13.5 gauge
- 3 inch diameter springs

Please note, the mattress has no pressure relieving properties and is not suitable for a profiling bed. The vandal proof and Fire-Tech versions of this mattress are also available on request for challenging environments, or environments at increased risk of fire.

TURNING

The National Bed Federation recommend that all mattresses should be turned once a month to allow for an even wear on the fillings, and to prevent mattress settlement. Turning a mattress is a condition of our manufacturer's warranty and may invalidate the warranty if this is not done.

CLEANING

The mattress should be wiped down with a solution of mild disinfectant and water. Do not use neat disinfectant as this will eventually affect the waterproof membrane. The mattress should not be left for a long time with incontinent patients on them as in extreme cases, fluid has been known to seep to the edges of the mattress and through the tape edging.



127kg/20st
max user weight



Our Divan bed is a great affordable option for your nursing home or community residents who are not in need of a profiling bed.

Product Information

Product Code	Description	Dimensions
HLB709	Divan Bed Base	1900(L)x930(W)x370(H)mm
HLB708	Divan Mattress	1900(L)x930(W)x160(H)mm
HLB7092	Divan Bed Headboard	900mm(W)

SYSTEM RANGE STATIC MATTRESSES




harvesthealthcare
a prismhealthcare company

The System range is made of high-quality moulded foam offering high maximum patient weights and risk levels. Extra features create superior levels of comfort and pressure relief.

Anatomical Shape Insert


Regular Viscoelastic Foam Mattress

The viscoelastic foam is strongly compressed in the area at risk, especially when profiled.

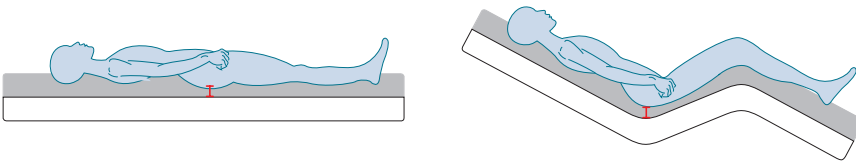

System Mattresses

As the insert follows the body shape in the sacrum area, the thickness of the viscoelastic foam is maximised to foster better prevention in areas at risk.

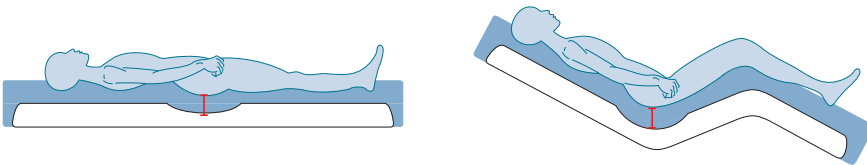

Multibearing Evolutive Concept


The multibearing surface enables a better pressure distribution through the transfer of support from areas with low bearing capacity toward areas with high bearing capacity.

White Area:


Areas with low bearing capacity (softer foam) giving increased prevention to areas at risk of pressure ulcerations (heels).

Blue Areas

Areas with high bearing capacity (firmer foam) for areas of the body with a lower level of pressure ulceration risk.

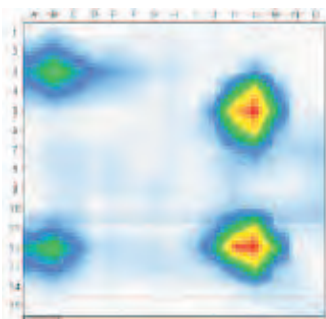

Viscoelastic Foam
Mattress

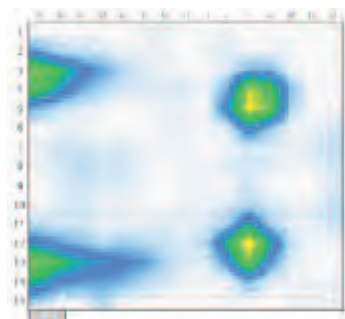

Viscoflex
Evolutive



Pressure
Transfer

VISCOFLEX STATIC MATTRESS




The Viscoflex mattress is a high-resilience mattress, moulded to ensure bonding. The mattress has an anatomically shaped base with a memory foam top.

The very high density top increases the contact area and improves pressure distribution with increased depth in sacral areas. The high density base works with the top to enhance support and comfort.

KEY FEATURES & BENEFITS

- Very high density viscoelastic foam
- Uniquely shaped to accommodate bed accessories
- Four-way stretch, waterproof, anti-bacterial PU coated cover with welded seams and a crib 5 fire rating
- Comes with a 3-year warranty
- High to very high risk

3 year warranty

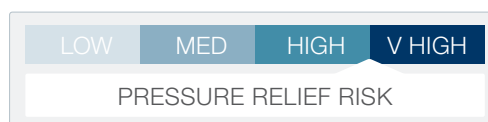

Product Information

Product Code	Description	Mattress Dimensions	Max Patient Weight
P161MFPMHF	Viscoflex 860mm Wide	1990(L)x860(W)x140(H)mm	130kg/20st
P161MFPMHF120	Viscoflex 1200mm Wide	2000(L)x1200(W)x140(H)mm	270kg/42st

VISCOFLEX MULTI BEARING STATIC MATTRESS



The Viscoflex Mattress with Multibearing Surface has all the great features of the Viscoflex mattress with an added multibearing surface.

The multibearing surface enables better pressure distribution through the transfer of supports from the areas of low bearing capacity toward areas with high bearing capacity. White areas with low bearing capacity (softer foam) give increased prevention to areas at risk of pressure ulcerations (heels). Blue areas with high bearing capacity (firmer foam) are for areas of the body with a lower level of pressure ulceration risk.

KEY FEATURES & BENEFITS

- Very high density viscoelastic foam
- Added multibearing surface
- Uniquely shaped to accommodate bed accessories
- Four-way stretch, water resistant, anti-bacterial PU coated cover with welded seams and a crib 5 fire rating
- Very high risk

3 year warranty

 130kg/20st
max user weight



Product Information

Product Code	Description	Mattress Dimensions
P161MFPPMHF	Viscoflex Multi Bearing 860mm Wide	1990(L)x860(W)x140(H)mm

VISCOFLEX EVOLUTIVE STATIC MATTRESS



The Viscoflex Evolutive static mattress is our enhanced mattress in the Viscoflex range.

Positioning specialist foam inserts across key areas of the mattress allows for more pressure relief and comfort. The varying firmness in foam at zonal areas contributes to the additional support.

The visco-elastic memory foam moulds to the users' body shape and the high resilience foam base increases stability.

The multibearing surface enables better pressure distribution through the transfer of support from the areas of low bearing capacity toward areas with high bearing capacity.

KEY FEATURES & BENEFITS

- Very high density viscoelastic foam
- Added multibearing surfaces
- Uniquely shaped to accommodate bed accessories
- Supplied with a four-way stretch, fully welded, waterproof, PU coated, polyester cover
- Very high risk

3 year warranty

130kg/20st
max user weight



Product Information

Product Code	Description	Mattress Dimensions
P161MFEPMHF	Viscoflex Evolutive 860mm Wide	1970(L)x860(W)x140(H)mm

KINETIC AIR HYBRID MATTRESS



The Kinetic Air hybrid mattress is a replacement mattress system, combining air displacement technologies with the benefits of a static foam mattress.

A castellated foam top and high resilience base aid in comfort and support for the patient. The walls of the mattress provide additional support by containing the foam core.

As the patient moves, internal air pressure is redistributed through a system of valves. This redistribution prevents the patient from overheating whilst laying on the mattress, as the foam can breathe.

KEY FEATURES & BENEFITS

- Suitable for fully profiling beds
- Castellated visco foam top and high resilience foam base and walls
- Combination of softer and firmer cells at key locations for full and comfortable support
- Fully welded, antibacterial, water resistant, four way stretch cover
- High risk

3 year warranty



247kg/39st
max user weight



LOW MED HIGH V HIGH

PRESSURE RELIEF RISK

Product Information

Product Code	Description	Mattress Dimensions
HFK100	Kinetic Air Hybrid Mattress	2000L)x900(W)x150(H)mm

PRIME COMFORT ACTIVE HYBRID MATTRESS



The Prime Comfort Active mattress combines all the benefits of our Prime Comfort Plus mattress, plus so much more. The addition of cells gives the ability to step up or step down patient care with or without a pump.

Built using both static foam and alternating air cells, the mattress can be used with or without the pump. This combination means patients at a very high risk of developing pressure ulcers, will still be protected when the mattress is deflated.

Harvest 1 pump can be purchased separately.

KEY FEATURES & BENEFITS

- Firm 'U' shaped high density foam base for roll prevention and to aid transfers
- Castellated soft and resilient central foam insert
- 10 alternating cells powered by the low noise level pump
- Fully welded, antibacterial, water resistant, four way stretch cover
- Durable PVC base
- Step up step down with/without pump
- Compatible with Harvest pumps 1 and 2, allowing for a range of maximum patient weights
- Very high risk with pump
- High risk without pump

3 year warranty



247kg/39st
max user weight



Compatible with:



The Harvest 1
pump



The Harvest 2
pump



Pump Grip (Small - Compatible with Harvest 1 Pump) - PU115

Pump Grip (Large - Compatible with Harvest 2 Pump) - PU215

Product Information

Product Code	Description	Mattress Dimensions	Mattress Weight
HHMPRIMA.P1	Prime Comfort Active only for Harvest 1 Pump	2000(L)x900(W)x140(H)mm	9kg
HHMPRIMA.P2	Prime Comfort Active only for Harvest 2 Pump	2000(L)x900(W)x140(H)mm	9kg

KENSINGTON OVERLAY ACTIVE MATTRESS



The Kensington Active Mattress is Harvest Healthcare's alternating overlay mattress, designed for patients at high risk of developing pressure ulcers. This active overlay is designed to be used on top of the patient's existing static mattress. It should not be used on its own.

The Kensington is light and easily manoeuvrable. If you have many patients in need of pressure reduction and relief, this overlay mattress is perfect for community use and provides a cost-effective, high-quality solution.

KEY FEATURES & BENEFITS

- 18 cells: 16 alternating cells made of TPU material to reduce noise when cells inflate and deflate and 2 static cells for patient stability and comfort
- Pillow section for patient stability
- All cells interchangeable for ease of repair and cleaning
- Low weight to allow the carer to easily move the mattress
- 4-way stretch welded PU cover which is vapour permeable, water resistant and has a Crib 5 fire rating
- High risk
- Compatible with Harvest pumps 1 and 2, allowing for a range of maximum patient weights
- Mattress weight: 6kg

3 year warranty



LOW MED HIGH V HIGH

PRESSURE RELIEF RISK

Compatible with:



The Harvest 1 pump



The Harvest 2 pump



Pump Grip (Small - Compatible with Harvest 1 Pump) - PU115

Pump Grip (Large - Compatible with Harvest 2 Pump) - PU215

Product Information

Product Code	Description	Mattress Dimensions	Max. Patient Weight
HAMKENS001.P1	Kensington Overlay only for Harvest 1 Pump	2000(L)x900(W)x110(H)mm	158kg/25st
HAMKENS001.P2	Kensington Overlay only for Harvest 2 Pump	2000(L)x900(W)x110(H)mm	184kg/29st

SANDRINGHAM REPLACEMENT ACTIVE MATTRESS



The Sandringham Active Mattress is Harvest Healthcare's basic replacement mattress with 18 alternating cells and an integral 50mm layer of foam providing an overall depth of 160cm. A pillow section has 2 static cells for stability and comfort. The cells are made of TPU material to reduce noise disturbance when cells alternately inflate and deflate on a 10-minute cycle.

Being easy to maintain with interchangeable cells, this product provides the customer with excellent value for money when choosing an alternating mattress replacement to achieve pressure reduction for high risk patients in long-term care.

KEY FEATURES & BENEFITS

- 18 cells: 16 alternating cells made of TPU material to reduce noise when cells inflate and deflate and 2 static cells for patient stability and comfort
- 50mm insert of high-quality foam
- Pillow section for patient stability
- All cells interchangeable for ease of repair and cleaning
- Low weight to allow the carer to easily move the mattress
- 4-way stretch welded PU cover which is vapour permeable, water resistant and has a Crib 5 fire rating
- High risk
- Compatible with Harvest pumps 1 and 2, allowing for a range of maximum patient weights
- Mattress weight: 9kg

3 year warranty



LOW MED HIGH V HIGH

PRESSURE RELIEF RISK

Compatible with:



The Harvest 1 pump



The Harvest 2 pump



Pump Grip (Small - Compatible with Harvest 1 Pump) - PU115

Pump Grip (Large - Compatible with Harvest 2 Pump) - PU215

Product Information

Product Code	Description	Mattress Dimensions	Max. Patient Weight
HAMSAND001.P1	Sandringham Mattress only for Harvest 1 Pump	2000(L)x900(W)x160(H)mm	158kg/25st
HAMSAND001.P2	Sandringham Mattress only for Harvest 2 Pump	2000(L)x900(W)x160(H)mm	184kg/29st

HAMPTON REPLACEMENT ACTIVE MATTRESS



The Hampton is Harvest Healthcare's durable active replacement mattress, designed for the prevention and management of pressure ulcers in patients at very high risk. It is also recommended for the treatment of up to stage 4 of pressure ulcers.

The mattress is made up of 20 individual cells. 17 alternating air cells and 3 static pillow section cells. The cells are made of TPU material to reduce noise disturbance when cells alternately inflate and deflate on a 10-minute cycle.

The Hampton mattress is an all-rounder for comfort and pressure area care solutions.

KEY FEATURES & BENEFITS

- 20 cells made of TPU material to reduce noise when cells inflate and deflate
- Pillow section for patient stability and comfort
- All cells interchangeable for ease of repair and cleaning
- 4-way stretch welded PU cover which is vapour permeable, water resistant and has a Crib 5 fire rating
- Very high risk
- Compatible with all Harvest pumps, allowing for a range of maximum patient weights
- Mattress weight: 8kg

3 year warranty



Compatible with:



The Harvest 1 pump



The Harvest 2 pump



The Harvest 3 pump



Pump Grip (Small - Compatible with Harvest 1 Pump) - PU115

Pump Grip (Large - Compatible with Harvest 2 Pump) - PU215

Product Information

Product Code	Description	Mattress Dimensions	Max. Patient Weight
HAMHAMP001.P1	Hampton Mattress only for Harvest 1 Pump	2000(L)x900(W)x200(H)mm	158kg/25st
HAMHAMP001.P2	Hampton Mattress only for Harvest 2 Pump	2000(L)x900(W)x200(H)mm	184kg/29st
HAMHAMP001.P3	Hampton Mattress only for Harvest 3 Pump	2000(L)x900(W)x200(H)mm	254kg/40st

WINDSOR REPLACEMENT ACTIVE MATTRESS



The Windsor is Harvest Healthcare's active replacement mattress designed for the prevention and management of pressure ulcers in patients at very high risk.

The mattress has 17 alternating air cells and 3 static pillow section cells for patient stability.

The Windsor mattress is an all-rounder for pressure area care solutions. Being made of a very soft and quiet material, it has got the critical factors to gain patient cooperation and concordance with pressure reduction therapy. It provides high levels of comfort and pressure care for users who find other alternating mattresses uncomfortable.

KEY FEATURES & BENEFITS

- 20 cells made of pure PU, a softer and quieter material for patient comfort
- Base sheet is a super easy wipe down material and has a Crib 5 fire rating
- Pillow section for patient stability
- All cells interchangeable for ease of repair and cleaning
- 4-way stretch welded PU cover which is vapour permeable, welded, water resistant and has a Crib 5 fire rating
- The PU resists water to make the cells easier to clean
- Very high risk
- Compatible with Harvest 2 and Harvest 3 pumps
- Cable management system
- Mattress weight: 10kg

3 year warranty



LOW MED HIGH V HIGH

PRESSURE RELIEF RISK

Compatible with:



The Harvest 2 pump



The Harvest 3 pump



Pump Grip (Large - Compatible with Harvest 2 Pump) - PU215

Product Information

Product Code	Description	Mattress Dimensions	Max. Patient Weight
HAMWIND001.P2	Windsor Mattress only for Harvest 2 Pump	2000(L)x900(W)x200(H)mm	184kg/29st
HAMWIND001.P3	Windsor Mattress only for Harvest 3 Pump	2000(L)x900(W)x200(H)mm	254kg/40st

BALMORAL REPLACEMENT ACTIVE MATTRESS



The Balmoral Replacement is Harvest Healthcare's top of the range pressure reducing active mattress used in the prevention and management of pressure ulcers. It is recommended for use by a patient who is at very high risk of developing pressure ulcers, and for the treatment of up to stage 4 pressure ulcers.

It is made up of 27 individual cells: 23 alternating air cells and 4 static pillow section cells, on a 50mm thick layer of foam. Cells are made of pure PU and are all interchangeable for ease of repair and cleaning. In addition, the mattress is designed to ensure bedrail compliance on any UK manufactured profiling bed.

KEY FEATURES & BENEFITS

- 27 individual cells made of pure PU, a softer and quieter material for patient comfort
- The mattress has an integrated foam base
- Pillow section for patient stability
- All cells are interchangeable for ease of repair and cleaning
- 4-way stretch welded PU cover which is vapour permeable and water resistant and has a Crib 5 fire rating
- Compatible with Harvest pumps 2 and 3
- Very high risk
- Mattress weight: 8kg

3 year warranty



LOW MED HIGH V HIGH

PRESSURE RELIEF RISK

Compatible with:



The Harvest 2 pump



The Harvest 3 pump



Pump Grip (Large - Compatible with Harvest 2 Pump) - PU215

Product Information

Product Code	Description	Mattress Dimensions	Max. Patient Weight
HAMBALM001.P2	Balmoral Mattress only for Harvest 2 Pump	2000(L)x880(W)x150(H)mm	184kg/29st
HAMBALM001.P3	Balmoral Mattress only for Harvest 3 Pump	2000(L)x880(W)x150(H)mm	254kg/40st

DUKE REPLACEMENT ACTIVE MATTRESS



The Duke Active Mattress is a high-quality replacement mattress. Designed for a frailer user who could be vulnerable to developing pressure ulcers, the system provides a very high degree of pressure relief whilst protecting the user. It is recommended for use by a patient who is at very high risk of developing pressure ulcers, and for the treatment of up to stage 4 pressure ulcers.

The mattress is made up of 25 alternating air cells and 5 static pillow section cells. A 50mm integral foam base rests beneath the removable cell structure, providing more support and redistributing pressure where necessary.

KEY FEATURES & BENEFITS

- 30 independent cells provides more comfort and support
- Alternating pressure relief distributes the bodyweight of the service user to provide treatment for pressure ulcers
- Pillow section for patient stability
- All cells interchangeable for ease of repair and cleaning
- 4-way stretch welded PU cover which is vapour permeable, welded, water resistant and has a Crib 5 fire rating
- The PU resists water to make the cells easier to clean
- Very high risk
- Substantial pressure relief for service users
- Compatible with Harvest 2 and Harvest 3 pumps
- Cable management system
- Mattress weight: 9kg

3 year warranty



LOW MED HIGH V HIGH

PRESSURE RELIEF RISK

Compatible with:



The Harvest 2 pump



The Harvest 3 pump



Pump Grip (Large - Compatible with Harvest 2 Pump) - PU215

Product Information

Product Code	Description	Mattress Dimensions	Max. Patient Weight
HR440.P2	Duke Mattress only for Harvest 2 Pump	1980(L)x900(W)x140(H)mm	165kg/26st
HR440.P3	Duke Mattress only for Harvest 3 Pump	1980(L)x900(W)x140(H)mm	247kg/39st

PARTNERSHIP ACTIVE MATTRESS



The Partnership is a double sized active mattress, made up of both foam and alternating air cells – one on either side. Half is comprised of viscoelastic foam, which provides pressure reduction, comfort and support for the patient. The alternating air cells make up the other half, with a 10-minute cell cycle of inflation and deflation to provide effective pressure reduction.

It is a great solution for couples who share a bed and one of the partners is at high risk of developing pressure ulcers. The foam section provides comfort and support for the patient's partner, whilst the alternating cells provide pressure reduction and relief. This is the perfect mattress for high risk patients who have a double bed but require a solution for pressure reduction and relief to either prevent or promote healing of pressure ulcers.

KEY FEATURES & BENEFITS

- 15 independent cells
- Alternating pressure relief
- Visco topped static mattress and foam walls
- Maintains quality of life
- Removable cell structure aids decontamination
- Electronic pressure control
- Timed static mode
- Low running costs
- Quiet Operation
- CPR quick release facility
- 4-way stretch welded PU cover which is vapour permeable, welded, water resistant and has a Crib 5 fire rating
- High risk
- Compatible with Harvest 1 and Harvest 2 pumps

3 year warranty



LOW MED HIGH V HIGH

PRESSURE RELIEF RISK

Compatible with:



The Harvest 1
pump



The Harvest 2
pump



Pump Grip (Small - Compatible with Harvest 1 Pump) - PU115

Pump Grip (Large - Compatible with Harvest 2 Pump) - PU215

Product Information

Product Code	Description	Mattress Dimensions	Max. Patient Weight
HS500.P1	Partnership Mattress only for Harvest 1 Pump	1980(L)x1370(W)x180(H)mm	158kg/25st
HS500.P2	Partnership Mattress only for Harvest 2 Pump	1980(L)x1370(W)x180(H)mm	184kg/29st

HAMPTON EXTRA REPLACEMENT PLUS SIZE ACTIVE MATTRESS



The Hampton Extra replacement is Harvest Healthcare's plus size active mattress. Designed for plus size patients at very high risk of developing pressure ulcers.

It has 17 alternating cells made of Nylon with a TPU lining to reduce noise when the cells inflate and deflate and a pillow section with 3 static cells for patient stability.

The foam surround allows for a wider mattress to be used on a 1200 wide bed, like the Community 1200 and the Ultimate 1200 and provides extra stability for transfers. It also improves stability when transferring from bed to chair, and prevents bottoming out. This mattress is essential for those wanting to provide the right equipment for plus size patients.

KEY FEATURES & BENEFITS

- 20 cells made of Nylon with a TPU lining to reduce noise when cells inflate and deflate
- Pillow section for patient stability
- All cells interchangeable for ease of repair and cleaning
- 4-way stretch welded PU cover which is vapour permeable, water resistant and has a Crib 5 fire rating
- Foam surround
- Very high risk
- Compatible with Harvest 3 pump

3 year warranty



Compatible with:



The Harvest 3 pump

Product Information

Product Code	Description	Mattress Dimensions	Max. Patient Weight
HAMHAMP002.P3	Hampton Extra Mattress only for Harvest 3 Pump	2000(L)x1200(W)x230(H)mm	317kg/50st

DUKE EXTRA 2.0 REPLACEMENT PLUS SIZE ACTIVE MATTRESS



The Duke Extra 2.0 is an alternating pressure relieving mattress system used in the prevention and treatment of pressure ulcers. It is recommended for use by a patient who is at very high risk of developing pressure ulcers, and for the treatment of up to stage 4 pressure ulcers.

The mattress is made up of 25 alternating air cells and 4 static pillow section cells at the head end. With a maximum patient weight of 50 Stone, this mattress is suitable for larger patients.

Edge to edge alternating cells provide pressure reduction across the full width of the mattress offering alternating pressure care at all points. It has been designed to fit the wider beds in our product range such as the Community 1200 and the Ultimate 1200 wide profiling beds.

KEY FEATURES & BENEFITS

- 29 cells provide more comfort and support
- Alternating pressure relief
- CPR rapid deflation valve
- 4-way stretch welded PU cover which is vapour permeable, water resistant and has a Crib 5 fire rating
- Quiet operation
- Electronic pressure control
- Very high risk
- Compatible with Harvest 3 pump

3 year warranty



Compatible with:



The Harvest 3 pump

Product Information

Product Code	Description	Mattress Dimensions	Max. Patient Weight
HAMDUKE.12.P3	Duke Extra 2.0 Mattress only for Harvest 3 Pump	2000(L)x1200(W)x200(H)mm	317kg/50st

WONDERMAT ACTIVE MATTRESS



The Wondermat Active Mattress is designed to constantly maintain interface pressures, below capillary closure powered by the Wondermat pump for very high risk patients.

The low air loss Wondermat Replacement System is designed to constantly maintain low interface pressures to prevent capillary closure.

This mattress is aimed at patients who are at the highest risk of developing pressure ulcers and moisture damage and need equipment that is right for them in terms of pressure reduction and relief, control of moisture, and optimised comfort.

KEY FEATURES & BENEFITS

- Rapid auto firm
- Tri-zonal feature
- Fowler mode
- Electronic pressure control
- Lockout to avoid unauthorised changes to settings
- Low pressure / power failure visual and audible alarms
- Self diagnosis check
- CPR quick release facility
- Complete with quilted vapour permeable water resistant cover
- Transport facility
- Powered by the Wondermat pump
- Suitable for fully profiling beds
- Very high risk

3 year warranty



LOW MED HIGH V HIGH

PRESSURE RELIEF RISK

Product Information

Product Code	Description	Mattress Dimensions	Max. Patient Weight
HW600	Wondermat Active System with Wondermat Pump	2060(L)x900(W)x220(H)mm	190kg/30st

ROTATIONAL ACTIVE MATTRESS




The Rotational Active Mattress is a high-quality, replacement active mattress system which offers three-in-one therapy for the effective prevention of pressure ulcers. This includes lateral turning therapy, alternating therapy, and continuous low-pressure therapy. Designed to assist with pulmonary complications and treatment of existing pressure ulcers, the mattress is suitable for patients in need of regular re-positional changes where this cannot be achieved by the patient or carers.

During lateral turning therapy, the carer can achieve easier repositioning with support from the mattress. The mattress helps facilitate postural drainage, especially for immobile patients. Also, during lateral turning therapy, the alternating therapy can be activated simultaneously. This function improves circulation and skin reperfusion and provides effective pressure relief via selectable cycle time. The continuous low-pressure therapy can provide body immersion for optimal comfort as a larger skin area will be supported by lower air pressure. This active mattress can significantly lower the risk of pressure ulcers while helping carers facilitate repositioning with this automatic system. Full training on use is provided by the team at Harvest and this system is compatible with the Harvest profiling beds.

KEY FEATURES & BENEFITS

- Up to 30° of bilateral turning therapy to assist patient repositioning and enhance pressure relief
- Offering extra support in the sacral area when in and during sitting position
- Fast inflation time & electronic pressure control
- Heel pressure relief function
- Improving circulation & skin reperfusion
- Lockout to avoid unauthorised changes to settings
- Low pressure / power failure visual and audible alarms
- CPR quick release facility
- Complete with low friction and sheer vapour permeable, high flame retardant, water resistant cover
- Transport mode
- Cable management for a safe operational environment
- Suitable for fully profiling beds
- High to very high risk

3 year warranty

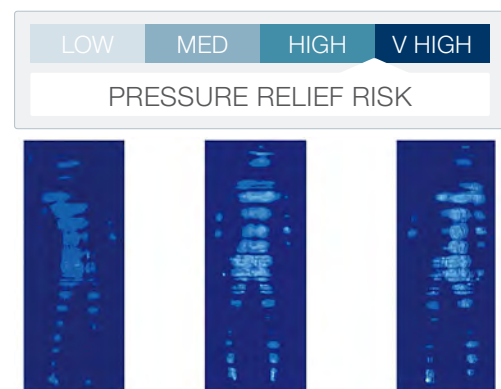

Pressure Mapping of 30° bilateral turning (weight 75kg)

Product Information

Product Code	Description	Mattress Dimensions	Max. Patient Weight
HR800	Rotational Active Mattress with Pump	2000(L)x900(W)x250(H)mm	180kg/28st
Pump Dimensions	Pump Weight	Pump Voltage	Pump Cycle Time
341 x 165 x 260mm	4.5kg	AC230v/50Hz	10, 15, 20, 25 min.
			Pressure Setting
			LED

HARVEST PUMPS



HARVEST 1 PUMP

The Harvest 1 Pump is Harvest Healthcare's manually operated pump. It is ideal for use in care home settings as it has a simple interface with a dial. It is easy to connect and disconnect and is lightweight, creating ease of use for the carer. It has visual and audible low pressure alarms for safety.



Product Information

Product Code	Dimensions	Pump Weight	Cycle Time	Air Flow Output	Pressure Setting
HPU1	100 x 130 x 280mm	1.7kg	10 minutes	7 lpm	Manual

PUMP POWER REQUIREMENTS

Power rating: 8W- Voltage: AC230v/50Hz- Fuse:1A
Electrical Safety Standard: EN 60601-1. EN 60601-1-2

HARVEST 2 PUMP

The Harvest 2 is Harvest Healthcare's digital pump. It is ideal for use in community and care home settings as it is easy to use.

It has a power failure alarm, L.E.D display, lock-out mode and static mode.



Product Information

Product Code	Dimensions	Pump Weight	Cycle Time	Air Flow Output	Pressure Setting
HPU2	130 x 140 x 270mm	1.7kg	10 minutes	9 lpm	LED

PUMP POWER REQUIREMENTS

Power rating: 9VA - Voltage: AC230v/50Hz - Fuse: T 1AH 250V
Electrical Safety Standard: EN 60601-1. EN 60601-1-2

HARVEST 3 PUMP

The Harvest 3 Pump is Harvest Healthcare's high-end digital pump. It is ideal for use in all settings including acute. It caters for plus size products and more complicated clinical needs. It is easy to use and has an L.E.D display, lock-out mode, static mode, maxi firm mode, 4 cycle times and visual and audible low pressure alarms for safety.



Product Information

Product Code	Dimensions	Pump Weight	Cycle Time	Air Flow Output	Pressure Setting
HPU3	340 x 125 x 220mm	3.5kg	10/15/20/ 25 min	12 lpm	LED

PUMP POWER REQUIREMENTS

Power rating: 8W- Voltage: AC230v/50Hz- Fuse:T1A 250v- Pressure range: 20-765mmHg
Electrical Safety Standard: EN 60601-1. EN 60601-1-2

FALL MATS UNDERLAYS AND OVERLAYS



FOLDING FALL MAT

The Harvest Folding Fall Mat offers the same level of safety as the standard fall mat, with additional advantages.


Manufactured using similar foam as our static mattresses, in the event of a patient fall, a soft landing is provided. This particular fall mat can be folded and stored away when not in use.

Fall mats are ideal for patients in Low or Ultra Low Beds; with the beds being closer to the floor itself, the patients can simply roll onto the fall mat.

The fall mats are complete with a vapour permeable cover and non-slip base.

Product Information

Product Code	Dimensions
HFM104NS	1800(L)x640(W)x64(H)

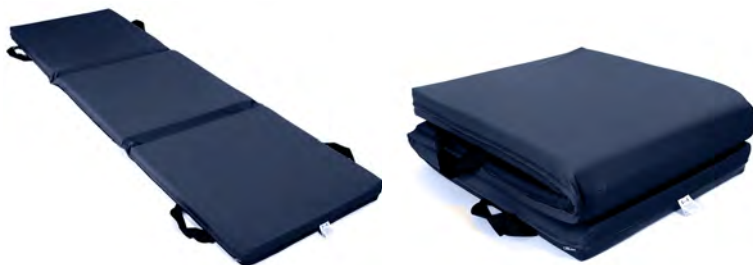

FALL MATS

Harvest's high density, foam fall mats provide additional safety to those who are at risk of rolling. Should users be in a Woburn Low or Woburn Ultra-Low Bed, without the use of side rails, the fall mat will provide a soft landing should they roll too far. Available in various sizes.



Product Information

Product Code	Description	Dimensions	Max. Weight
HFO122	Modular Foam Overlay	1880(L)x865(W)x95(H) mm	114kg/18st
HFO119	Modular Foam Overlay	1900(L)x1370(W)x95(H) mm	114kg/18st
HFU110	50mm White Foam Underlay	2000(L)x900(W)x50(H) mm	
HFC102	Standard Fall Mat	1830(L)x610(W)x100(H) mm	
HFC103	Standard Fall Mat	1420(L)x815(W)x64(H) mm	

COMMUNITY MODULAR FOAM CUSHION



The Community Modular Foam Cushion combines stability with a high level of comfort, offering the user more than expected.

Cut with specific castellations, the highly breathable foam can prevent the user from becoming too hot.

The cushion is complete with a zipped, vapour permeable, four-way stretch cover with an anti slip base.

KEY FEATURES & BENEFITS

- Castellated block foam
- Zipped, vapour permeable, water resistant, four-way stretch cover
- Non-slip base
- Can be paired with our Community mattress for 24 hour pressure care
- Available in 2 different sizes
- Medium risk

3 year warranty



Paired Mattress:
Community (HSMCOM)
Medium Risk - SWL:159kg

Product Information

Product Code	Dimensions	Paired Mattress	Safe Working Load
HFC121	430(L)x430(W)x95(H) mm	Community (HSMCOM)	95kg /15st
HFC116	500(L)x500(W)x95(H) mm	Community (HSMCOM)	108kg/17st

PRIME COMFORT VISCO TOP CUSHION




A 50mm layer of viscoelastic foam tops the standard foam base, transforming the cushion into a high risk product. The top layer of foam provides extra support and comfort to high risk patients.

The Prime Comfort Visco Top cushion has a zip with flap, four-way stretch cover with an anti slip base.

KEY FEATURES & BENEFITS

- Foam topped with viscofoam
- Zipped with flap, water resistant, four-way stretch cover
- Non-slip base
- Can be paired with our Prime Comfort mattress for 24 hour pressure care
- Available in 2 different sizes
- High risk

3 year warranty

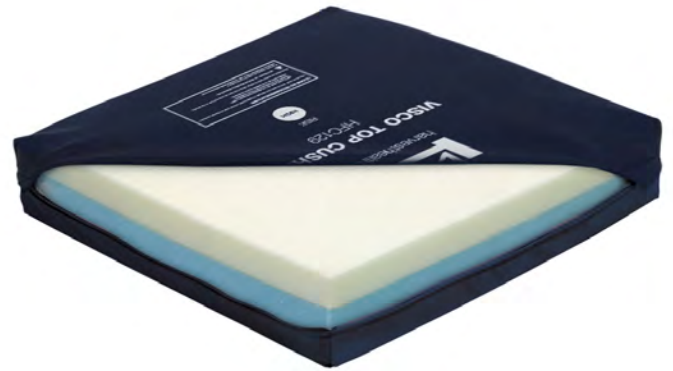

Paired Mattress:
Prime Comfort (HSMPRIM)
High Risk - SWL: 177kg

Product Information

Product Code	Dimensions	Paired Mattress	Safe Working Load
HFC129	430(L)x430(W)x79(H) mm	Prime Comfort (HSMPRIM)	127kg/20st
HFC123	510(L)x510(W)x76(H) mm	Prime Comfort (HSMPRIM)	139kg/22st

PRIME COMFORT PLUS VISCO TOP CUSHION



A high risk viscoelastic and foam bonded cushion with improved immersion of visco top and a combination of a second layer of viscoelastic foam and a firm foam base add support.

The foam combination enhances pressure distribution and combines stability and comfort for patients at high risk.

The Prime Comfort Plus Visco Top cushion has a navy blue, four-way stretch cover with an anti slip base.

KEY FEATURES & BENEFITS

- Foam topped and bottomed with viscofoam
- Welded flap, water resistant, four-way stretch cover
- Non-slip base
- Can be paired with our Prime Comfort Plus mattress for 24 hour pressure care
- High risk

3 year warranty



Paired Mattress:
Prime Comfort Plus
(HSMPRIMP)
High Risk - SWL: 247kg

Product Information

Product Code	Dimensions	Paired Mattress	Safe Working Load
HFC129.01	430(L)x430(W)x75(H) mm	Prime Comfort Plus (HSMPRIM)	158kg/25st

REFLECT CUSHION




A layer of castellated, viscoelastic foam rests on a high density foam base. The castellations allow the foam to breathe as well as providing more patient comfort. The base gives resilience and stability.

The Reflect cushion has a navy blue, four-way stretch, highly breathable cover with an anti slip base.

KEY FEATURES & BENEFITS

- Castellated viscofoam on foam
- Zipped, highly breathable, water resistant, four-way stretch cover
- Non-slip base
- Can be paired with our Reflect 2 mattress for 24 hour pressure care
- High risk

3 year warranty

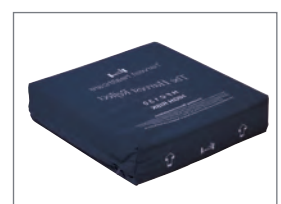

Paired Mattress:
Reflect 2 (HFR147)
High Risk - SWL: 247kg

Product Information

Product Code	Dimensions	Paired Mattress	Safe Working Load
HFC130	430(L)x430(W)x100(H) mm	Reflect 2 (HFR 147)	127kg/20st

VISCOFLEX PLUS CUSHION



The Viscoflex Plus cushion is a very high density viscoelastic, moulded cushion with an ergonomically shaped load bearing core.

The core is shaped to give optimal immersion in sacral/ischial areas, reducing pressure and facilitating blood circulation. The outer foam is moulded to promote stability and increase the contact area. The skin effect of the moulded foam also gives protection from tearing or crumbling.

This ultra-light cushion is specially designed to provide high levels of stability and pressure relief for the user. It is supplied with a vapour-permeable, water resistant cover.

KEY FEATURES & BENEFITS

- Injected foam
- Welded, water resistant, four-way stretch cover
- Can be paired with our Viscoflex Multibearing mattress for 24 hour pressure care
- Available in 3 different sizes
- High to very high risk

3 year warranty



Paired Mattress:
Viscoflex Multibearing
(P161MOPPMHF)

High-Very High Risk - SWL: 130kg



The inside of the
Viscoflex Plus Cushion

Product Information

Product Code	Dimensions	Paired Mattress	Safe Working Load
P361CAF42421HF	420(L)x420(W)x80(H) mm	Viscoflex Multibearing (P161MOPPMHF)	100kg/15st
P361CAF45421HF	420(L)x450(W)x80(H) mm	Viscoflex Multibearing (P161MOPPMHF)	110kg/17st
P361CAF51461HF	460(L)x510(W)x80(H) mm	Viscoflex Multibearing (P161MOPPMHF)	150kg/23st

MIXTE CUSHION




The Mixte cushion is a combination of gel and foam. Reticulated polyurethane gel avoids migration and thus prevents penetration (punching) of the cushion while assisting pressure distribution.

The thick high resilience foam base increases the contact area, providing better distribution, comfort and stability. The two layers are permanently bonded together, eliminating the risk of separation.

KEY FEATURES & BENEFITS

- Gel top on foam
- Welded, water resistant, four-way stretch cover
- Can be paired with our Prime Comfort mattress for 24 hour pressure care
- Medium to high risk

3 year warranty

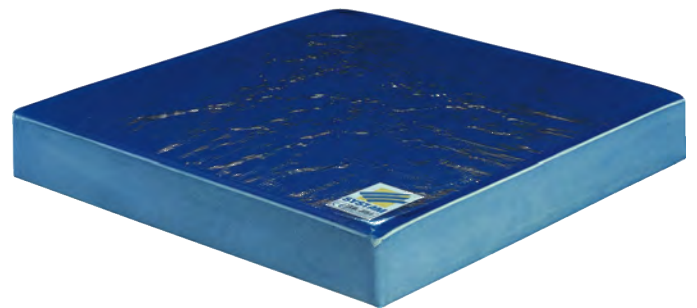

LOW MED HIGH V HIGH

PRESSURE RELIEF RISK



Paired Mattress:
Prime Comfort (HSMPRIM)
High Risk - SWL: 177kg

Product Information

Product Code	Dimensions	Paired Mattress	Safe Working Load
P314C40401HF	400(L)x400(W)x60(H) mm	Prime Comfort (HSMPRIM)	90kg/14st

DUOFORM CUSHION



The Duoform® cushion is an anatomically moulded base of high-resilience foam surmounted by a two-part viscoelastic gel. The white pad fluid gel ensures pressure reduction on high risk areas, while the blue viscoelastic gel resists fluid migration thus preventing penetration and providing added support.

The high-resilience foam base enables the cushion to better adapt to bony protuberances, spreads contact surfaces and gives better pressure distribution. It has a highly breathable cover.

KEY FEATURES & BENEFITS

- Anatomically moulded two-part gel on foam
- Welded, four-way stretch, highly breathable cover
- Can be paired with our Prime Comfort Plus mattress for 24 hour pressure care
- Available in 4 different sizes
- Medium to high risk

3 year warranty



LOW MED HIGH V HIGH

PRESSURE RELIEF RISK



Paired Mattress:
Prime Comfort Plus (HSMPRIMP)
High Risk - SWL: 247kg

Product Information

Product Code	Dimensions	Paired Mattress	Safe Working Load
P341C42421HF	420(L)x420(W)x80(H) mm	Prime Comfort Plus (HSMPRIM)	110kg/17st
P341C45421HF	420(L)x450(W)x80(H) mm	Prime Comfort Plus (HSMPRIM)	120kg/18st
P341C45461HF	460(L)x450(W)x80(H) mm	Prime Comfort Plus (HSMPRIM)	120kg/18st
P341C51461HF	460(L)x510(W)x80(H) mm	Prime Comfort Plus (HSMPRIM)	160kg/25st

AIRSET CUSHION



The Harvest Healthcare AirSet Cushion is a state-of-the-art pressure care solution that boasts an innovative hybrid design, combining the benefits of foam and air to provide optimal pressure relief.

Users can easily adjust the cushion to their desired level of pressure relief, thanks to its self-inflating technology, ensuring maximum comfort and support. The cushion is climate-controlled, helping maintain a comfortable seating temperature, reducing the risk of heat build-up and pressure ulcers. The AirSet Cushion's unique construction allows it to distribute pressure evenly across the seating surface, promoting healthy circulation and reducing the risk of injury.

KEY FEATURES & BENEFITS

- 36 Cells
- Innovative hybrid design, combining the benefits of foam and air
- Self-inflating technology for easy pressure adjustment
- Climate-controlled for a comfortable seating temperature
- Even pressure distribution across the seating surface
- Promotes healthy circulation
- Washable with fully enclosed foam

3 year warranty



Product Information

Product Code	Dimensions	Safe Working Load
HHC001	400(L)x400(W)x70(H) mm	120kg/19st

VISCOFLEX PLUS BARIATRIC CUSHION

Plus Size
Static Cushions

69



The innovative Viscoflex Plus Bariatric cushion consists of a moulded Viscoelastic foam with an ergonomic foam insert.


This plus size cushion is specially designed to provide high levels of stability and pressure relief for the user.

P361CAF56461HF is supplied with a welded, vapour-permeable, water resistant cover and P361CA61501HF is supplied with a welded, vapour-permeable, water resistant cover with an anti slip base.


KEY FEATURES & BENEFITS

- Injected foam
- Available in two sizes
- Welded vapour-permeable, water resistant cover for P361CAF56461HF (560(L)x460(W)x80(H) mm)
- Stitched vapour-permeable, water resistant cover with an anti slip base available for P361CA61501HF (610(L)x500(W)x80(H) mm)
- Can be paired with our Viscoflex Bariatric mattress for 24 hour pressure care
- High to very high risk


3 year warranty

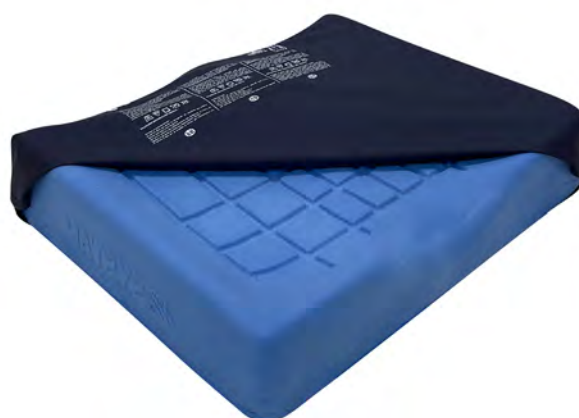

With Standard Cushion

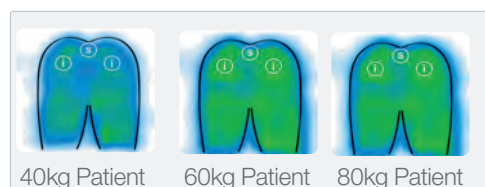

With Viscoflex Plus

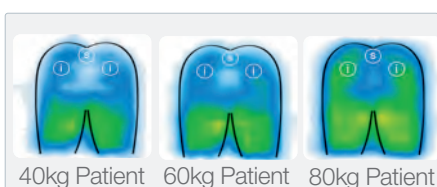

<p>ONE REDUCE THE LOAD ON RISK ZONES</p> <p>Protected risk zones (sacrum/ischium), even in the cases of pelvic retroversion, in order to provide partial relief of bony protuberances through support.</p>	<p>TWO ANTI-PUNCHING SYSTEM</p> <p>Specifically positioned in relation to the ischium. Distributes the pressure to the buttock peripheral zone, stabilising pelvis support and preventing it from reaching the bottom zone.</p>	<p>THREE PRE-ISCHIAL RIDGE</p> <p>Blocks the ischium, thus preventing forward sliding.</p>
---	--	---



Paired Mattress:
Viscoflex Bariatric (P161MFPMHF120)
High-Very High Risk - SWL: 285kg

*Please note information supplied here is a guide, if you would like more information or help with making the correct choice, please call our team today.

Product Information

Product Code	Dimensions	Paired Mattress	Safe Working Load
P361CA56461HF	460(L)x560(W)x80(H) mm	Viscoflex Bariatric (P161MFPMHF120)	230kg/36st
P361CA61501HF	500(L)x610(W)x80(H) mm	Viscoflex Bariatric (P161MFPMHF120)	270kg/42st

GEL CUSHION



The Gel Cushion is made of high resilience, impermeable, viscoelastic polyurethane gel with highly elastic properties.

It is designed to enhance existing support surfaces.

KEY FEATURES & BENEFITS

- Thin, gel only
- Welded, water resistant, four-way stretch cover
- Can be paired with our Community mattress for 24 hour pressure care
- Available in 2 different sizes
- Low to medium risk

3 year warranty



LOW MED HIGH V HIGH

PRESSURE RELIEF RISK



Paired Mattress:
Community (HSMCOM)
Medium Risk - SWL:159kg

Product Information

Product Code	Dimensions	Paired Mattress	Safe Working Load
P311C40401HF	400(L)x400(W)x25(H) mm	Community (HSMCOM)	120kg/18st
P311C45431HF	430(L)x450(W)x25(H) mm	Community (HSMCOM)	140kg/22st

ULTRA-THIN CUSHION



The Ultra-Thin Cushion has been developed for use with a range of chairs. The viscoelastic foam offers excellent comfort, support and pressure distribution.

Available separately, each 35mm thick cushion provides enhanced pressure distribution to the existing chair cushion, without compromising user posture. The thin nature of the cushion also ensures the patient remains safely in the chair as the arms can still be a barrier.

This lightweight cushion is also available as a leg and lumbar system and comes with a full strap system to secure it to the chair. Comes with a water resistant, vapour permeable cover with an anti slip base.

KEY FEATURES & BENEFITS

- Viscoelastic foam
- Zipped, vapour permeable, water resistant, four-way stretch cover
- Non-slip base
- Can be paired with our Reflect mattress for 24 hour pressure care
- Available in 2 different sizes
- Complete system consisting of seat, leg and lumbar pieces is also available (HFC143COM)
- High risk

3 year warranty



LOW MED **HIGH** V HIGH

PRESSURE RELIEF RISK



Paired Mattress:
Reflect 2 (HFR147)
High Risk - SWL: 247kg

Product Information

Product Code/Description	Dimensions	Paired Mattress	Safe Working Load
HFC143 - Seat Only	430(L)x430(W)x35(H) mm	Reflect 2 (HFR147)	152kg/23st
HFC143.20 - Seat Only (Larger)	508(L)x508(W)x35(H) mm	Reflect 2 (HFR147)	152kg/23st
HFC143LEG	440(L)x430(W)x25(H) mm	Reflect 2 (HFR147)	152kg/23st
HFC143LUM	700(L)x430(W)x25(H) mm	Reflect 2 (HFR147)	152kg/23st
HFC143COM (Seat, Leg and Lumbar Pieces)	3 Piece Set as Above	Reflect 2 (HFR147)	152kg/23st

Arm rests are available on request; Product Code: HFC143ARM

ACTIVE SEAT CUSHION



A one-piece alternating cushion that provides maximum patient comfort whilst achieving optimum pressure relief whilst seated.


Especially beneficial to patients who can no longer sleep in a hospital type bed.

The cushion has 6 alternating cells.

KEY FEATURES & BENEFITS

- 6 alternating cells
- Zipped, water resistant, four-way stretch cover
- Non-slip base
- Can be paired with our Hampton mattress for 24 hour pressure care
- Compatible with our Harvest 1 (analogue) and Harvest 2 (digital) pumps
- Very high risk

3 year warranty

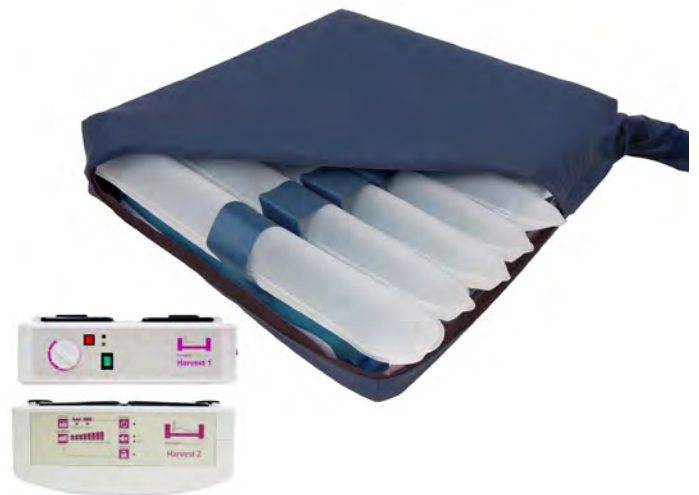

LOW MED HIGH **V HIGH**

PRESSURE RELIEF RISK



Paired Mattress:
Hampton (HAMHAMP001)
Very High Risk
SWL: HAMHAMP001.P1: 158kg
HAMHAMP001.P2: 184kg

Product Information

Product Code	Dimensions	Paired Mattress	Safe Working Load
HH20.P1	400(L)x430(W)x90(H) mm	Hampton (HAMHAMP001.P1)	158kg/25st
HH20.P2	400(L)x430(W)x90(H) mm	Hampton (HAMHAMP001.P2)	184kg/29st

BLENHEIM ACTIVE CUSHION



A one-piece alternating cushion that provides maximum patient comfort whilst achieving optimum pressure relief whilst seated.

Especially beneficial to patients who can no longer sleep in a hospital type bed.

The Blenheim Active cushion has 9 alternating cells.

KEY FEATURES & BENEFITS

- 9 alternating cells
- Welded, water resistant, four-way stretch cover
- Non-slip base
- Can be paired with our Balmoral mattress for 24 hour pressure care
- Compatible with our Harvest 2 (digital) pump
- Very high risk

3 year warranty



LOW MED HIGH **V HIGH**

PRESSURE RELIEF RISK



Paired Mattress:
Balmoral (HAMBALM001)
Very High Risk
SWL: HAMBALM001.P2: 184kg

Product Information

Product Code	Dimensions	Paired Mattress	Safe Working Load
HACBLEN001.P2	460(L)x460(W)x60(H) mm	Balmoral (HAMBALM001.P1)	152kg / 24st

3. POSITIONING AIDS

The Heel Pads and Positioning Aids from Harvest Healthcare, feature high-quality, ergonomically designed aids, which support and comfort patients. For best use when in beds and chairs, the positioning aids protect the relevant parts of the body from pressure damage and help to prevent skin damage at vulnerable locations.

INTEGRAL HEEL SUPPORT



The Integral Heel Support is made of viscoelastic foam with memory effect. The heel is cushioned by the moulded high density memory foam, which is shaped to protect wear on the achilles tendon and heel and direct pressure towards lower risk areas. The heel support is particularly well suited to patients with normal or slightly limited mobility. It is also ideal for prevention in agitated patients or those subject to constant movement of the lower limbs.

This comfort boot comes with a washable, breathable cover and is fitted using soft Hook and Loop straps, to tailor the perfect fit even to users with lower limb oedema.

KEY FEATURES & BENEFITS

- Heel area hollowed out to provide protection against friction effects
- Sides of the forefoot sections hollowed out to divert pressure towards the surrounding areas at lower risk
- Achilles tendon protected by side supports that direct pressure towards areas at lower risk

MATERIAL

Visco foam; moulded and of very high density (80kg/m³, medical grade) memory foam.

Able to mould to the body to provide better reduction of pressure on areas at high risk of ulcers.

1 year warranty



Product Information

Product Code	Description	Dimensions
P901T1HW	Integral Heel Support Viscoelastic Foam with Memory Foam	550(L) x 250 (W) x 75(H)mm

HEEL PAD



The Heel Pad is made of soft viscoelastic foam with memory effect and stronger viscoelastic foam for support. The multibearing surface concept shifts pressure from areas at high risk to areas at low risk and the end can be tucked under and secured if full pressure relief for the heels is required.

Available in single or double width, these clever positioning aids are excellent for protecting the heels and feet in those that are less mobile.

Well suited for people with moderate loss of mobility in lower limbs at medium to very high risk of pressure ulcers.

KEY FEATURES & BENEFITS

- Raised sides to avoid limbs slipping off
- Sloped to decrease pressure either partly or fully on the heel
- Long to support the knee and avoid hyperextension
- Width enables freedom of movement
- Available in half size to enable different positions for each leg

MATERIAL

Visco foam; moulded and of very high density (80kg/m³, medical grade) memory foam.

Able to mould to the body to provide better reduction of pressures on areas at high risk of ulcers.

1 year warranty



Product Information

Product Code	Description	Dimensions
P902T1HW	Heel Pad	730(L)x640(W)x105(H)mm
P904T1HW	Half Heel Pad	730(L)x320(W)x105(H)mm

HEEL SUPPORT BOOT



The Heel Support Boot is made of viscoelastic foam with memory effect. Well suited for people with low mobility in lower limbs to assist in the prevention and treatment of pressure ulcers. This zero heel pressure boot comes with a washable, breathable cover.

The high density visco foam is moulded for comfort and to prevent rotation and heel/tendon pressure, protecting all aspects of the foot.

KEY FEATURES & BENEFITS

- Heel area hollowed out for zero pressure
- Hollow area beneath the achilles tendon to reduce pressure
- Curved outline to distribute pressure evenly
- Prevents equinus distortion and rotation

MATERIAL

Visco foam; moulded and of very high density (80kg/m³, medical grade) memory foam.

Able to mould to the body to provide better reduction of pressures on areas at high risk of ulcers.

1 year warranty



Product Information

Product Code	Description	Dimensions
P903T1HW	Half Support Boot	635(L)x220(W)x310(H)mm

SEMI-FOWLER RETRACTION WEDGE



harvesthealthcare
a prismhealthcare company

The Semi-Fowler Retraction Wedge is made of high resiliency foam enabling positions to be maintained over time for effective treatment. Recommended for the prevention and treatment of pressure ulcers on the sacrum, this wedge can be used in conjunction with a range of heel-specific aids for more tailored care.

KEY FEATURES & BENEFITS


- Heel area hollowed out for zero pressure
- Hollow area beneath the achilles tendon to reduce pressure
- Curved outline to distribute pressure evenly
- Prevents equinus distortion and rotation

MATERIAL

Visco foam; moulded and of very high density (80kg/m³, medical grade) memory foam.

Able to mould to the body to provide better reduction of pressure on areas at high risk of ulcers.

1 year warranty

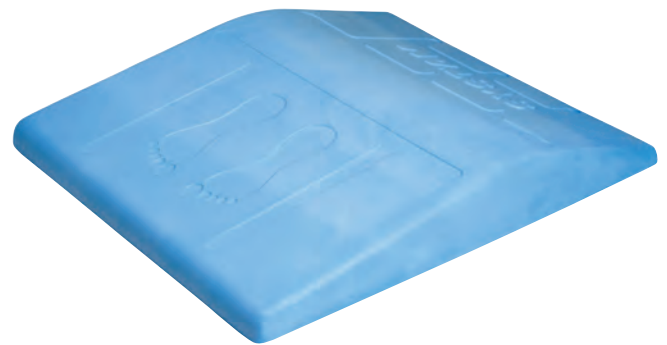

Product Information

Product Code	Description	Dimensions
P905T1HW	Semi-Fowler Retraction	620(L)x700(W)x75(H)mm

30° POSITIONING WEDGES



The Positioning Wedge is made of memory foam with a stabilisation section made of high resiliency foam, designed to keep residents comfortable and stable in a 30-degree position, whilst also maintaining pressure relief. This comfortable positioning aid reduces the need for pillows or other bedding being used to support turns, gives carers peace of mind that their resident is in the optimum position, and also helps the resident maintain position when unmonitored.

The stabilisation insert enables the 30° position to be maintained over time by acting as a support structure. Available full sized or in half for the left or right side.

KEY FEATURES & BENEFITS

- Supports lower limb at the thigh level
- Hollowed out in sacrum area to relieve pressure and reduce painful contact areas
- Anatomical shape moulds to the body
- Reduced width of wedge allows the patient to be centred on the bed, avoids accidental contact with side rails and stops it from becoming lodged between the side rails

MATERIAL

- Visco foam*
- High resiliency foam**

1 year warranty

Visco Foam: Moulded and of very high density (80kg/m³, medical grade) memory foam.

Able to mould to the body to provide better reduction of pressures on areas at high risk of ulcers.*

High Resiliency Foam: Its mechanical properties (elasticity, flexibility) and high density (40kg/m³) enable positions to be maintained over time, enabling effective treatment.**



Product Information

Product Code	Description	Dimensions
P910L1HW	30 Degree Lateral Positioning Wedge with Cover	990(L)x550(W)x150(H)mm
P910LD1HW	30 Degree Lateral Positioning Half Wedge RHS Tilt	990(L)x280(W)x150(H)mm
P910LG1HW	30 Degree Lateral Positioning Half Wedge LHS Tilt	990(L)x280(W)x150(H)mm

POSITIONING WEDGES FOR UPPER LIMBS



The Wedge is made of viscoelastic foam with memory effect. Designed to prevent and treat oedema. Available to support the whole arm or just the lower part.

Anatomically shaped to help comfort, can be used on both the left and right side, extension of the elbow in lying position avoids restricting vein perfusion in the forearm.

KEY FEATURES & BENEFITS

- Anatomically shaped to help comfort and maintain the position of the upper limbs over time
- Allows access for a drip to be used
- Can be used on both the left and right side
- Extension of the elbow in lying position avoids restricting vein perfusion in the forearm

MATERIAL

Visco foam; moulded and of very high density (80kg/m³, medical grade) memory foam.

Able to mould to the body to provide better reduction of pressure on areas at high risk of ulcers.

1 year warranty



Product Information

Product Code	Description	Dimensions
P912L1HW	Positioning wedge for Upper Limbs	650(L)x370(W)x160(H)mm
P915L1HW	Positioning wedge for Upper Limbs	490(L)x190(W)x160(H)mm

ABDUCTION WEDGE



The Abduction Wedge is made of high resiliency foam. Recommended for the prevention of postoperative dislocation and to help stabilise the legs in a position.



KEY FEATURES & BENEFITS

- Anatomically shaped to fit the natural angle between the lower limbs to ensure optimal distribution of pressure areas
- Fine curved tip to help carers position the wedge easily
- High enough for use in a chair and low enough to fit under a table

MATERIAL

High resiliency foam; its mechanical properties (elasticity, flexibility) and high density (40kg/m³) enable positions to be maintained over time, enabling effective treatment.

1 year warranty

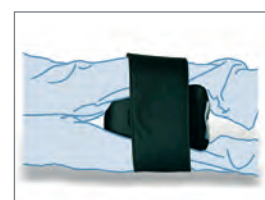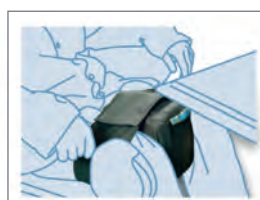

Product Information

Product Code	Description	Dimensions
P911L1HW	Abduction Wedge	240(L)x160(W)x180(H)mm

OCCIPITAL SYSTEM



The Occipital System provides prevention and positioning aid for the occipital area.

MATERIAL

Visco foam; moulded and of very high density (80kg/m³, medical grade) memory foam.

Able to mould to the body to provide better reduction of pressures on areas at high risk of ulcers.

Viscoelastic foam with shape memory;

Defined body moulding and increased body surface area in contact with the foam

More effective distribution of pressure; reduction of transcutaneous pressure on high risk zones, facilitated blood flow

Increased comfort and more effective maintenance of patient position

1 year warranty



Product Information

Product Code	Description	Dimensions
P913L1HW	Occipital Support	500(L)x430(W)x95(H)mm

UNIVERSAL CUSHION







The Universal Cushion provides postural aid and prevention for pressure ulcers. Helps to support the user and relieves pressure on the heels, protects the internal lateral knee areas and ankles. Can be used as postural aid in chairs and wheelchairs.

MATERIAL

Visco foam; moulded and of very high density (80kg/m³, medical grade) memory foam.

Able to mould to the body to provide better reduction of pressures on areas at high risk of ulcers.

1 year warranty

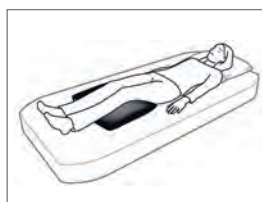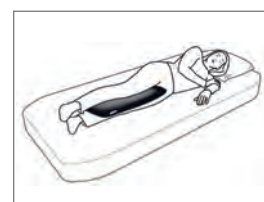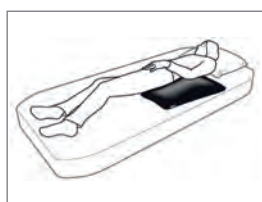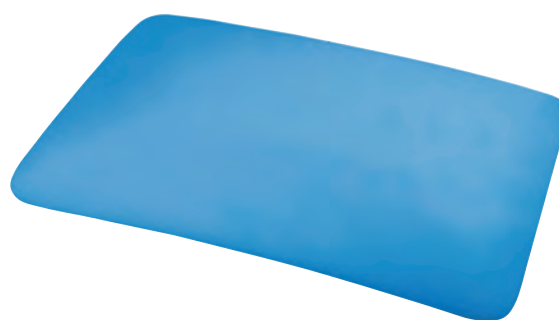

Product Information

Product Code	Description	Dimensions
P914L	Universal Cushion Large	600(L)x400(W)x130(H)mm

ABDUCTION BLOCK FOR KNEES CYLINDRICAL CUSHION



The **Abduction Block for Knees** is made of microbeads. The simultaneous use of the abduction block and cylindrical cushion achieves a comfortable Semi-Fowler position for the user. It also supports the dorsal position by firmly holding the pelvis and abducted areas in place thanks to the semi-circle cushion. Provides inner protection for the inner sides of the knees.

The **Cylindrical Cushion** is made of microbeads. Can be placed under the knees to maintain lower limbs in a Semi-Fowler position.

MATERIAL

Expanded polystyrene microbeads

Very light material: expanded polystyrene microbeads with small granulometry for enhanced comfort and adaptability of support to the patient's body.

Synthetic polyester fibre

Strong and resistant to stretching and shrinking. Retains air volume providing a soft sensation and increased comfort.

1 year warranty



Abduction Block for Knees



Cylindrical Cushion



Product Information

Product Code	Description	Dimensions
P9703B1HW	Abduction Block for Knees	620(L)x320(W)
P9702B1HW	Cylinder Cushion	600(L)x200(Diameter)mm

UNIVERSAL CUSHION

RING CUSHION



The **Universal Cushion** is made of microbeads. It helps support the patient, relieves the pressure on the heels, protects the internal lateral knee areas, ankles and can be used as a pillow. It also comes in two sizes, small and standard.

The **Ring Cushion** is also made of microbeads. It comfortably supports the head both in a sitting or semi-reclined position and has an easy-to-use, self gripping strip.

MATERIAL

Expanded polystyrene microbeads

Very light material: expanded polystyrene microbeads with small granulometry for enhanced comfort and adaptability of support to the patient's body.

Synthetic polyester fibre

Strong and resistant to stretching and shrinking. Retains air volume providing a soft sensation and increased comfort.

1 year warranty



Universal Small and Standard Cushion



Ring Cushion



Product Information

Product Code	Description	Dimensions
P9700B1HW	Universal Cushion	350(L)x260(W)mm Small / 560(L)x400(W)mm Standard
P9704B1HW	Ring Cushion	500(Diameter)mm

HALF-MOON CUSHION

LATERAL POSITIONING CUSHION



The **Half-Moon Cushion** is made of microbeads. It holds the pelvis and upper body firmly in place when the user is in the Semi-Fowler or dorsal position and relieves pressure on elbows. Fitted with a strap for extra support to keep the head in position, it can be used on a wheelchair.

The **Lateral Positioning Cushion** is made of microbeads. Designed with an integrated pillow for optimal comfort, it ensures comfortable positioning of the knees, pelvis and upper body when in the lateral position. It is compatible with all patient sizes.

MATERIAL

Expanded polystyrene microbeads
Very light material: expanded polystyrene microbeads with small granulometry for enhanced comfort and adaptability of support to the patient's body.

Synthetic polyester fibre
Strong and resistant to stretching and shrinking. Retains air volume providing a soft sensation and increased comfort.

Both of them use an Easy Shape system to maintain a position that is adapted to the patient's body, thereby increasing comfort. By pulling on the straps, the curvature of the cushion increases and the shape is maintained by tying the straps.

1 year warranty



Half-Moon Cushion



Lateral Positioning Cushion



Product Information

Product Code	Description	Dimensions
P9705B1HW	Half-Moon Cushion	1350(L)x850(W)mm
P9707B1HW	Lateral Positioning Cushion	1800(L)x540(W)mm

CROSS PAD WEDGE WITH FRONT/BACK SLOPE



The **Cross Pad** is made of microbeads. It enables hip and knee positioning as well as full pressure relief of the heels and has an adjustable strap to fix on the mattress.

The **Wedge with Front/Back Slope** enables the positioning of upper or lower limbs. The front/back slope relieves the pressure points on the heels or the hand.

MATERIAL

Expanded polystyrene microbeads

Very light material: expanded polystyrene microbeads with small granulometry for enhanced comfort and adaptability of support to the patient's body.

Synthetic polyester fibre

Strong and resistant to stretching and shrinking. Retains air volume providing a soft sensation and increased comfort.


1 year warranty



Cross Pad



Wedge with Front/Back Slope

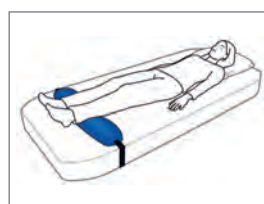

Product Information

Product Code	Description	Dimensions
P9709B1HW	Cross Pad	850(L)x250(W)x150(H)mm
P9712B1HW	Wedge with Front/Back Slope	460(L)x300(W)x150(H)mm



HIP ABDUCTION WEDGE DELTA CUSHION



harvesthealthcare
a prismhealthcare company

The **Hip Abduction Wedge** enables hip abduction and the reduction of shear and friction forces on the lower limbs.

The **Delta Cushion** is a positioning cushion for the Semi-Fowler position.

MATERIAL

Expanded polystyrene microbeads

Very light material: expanded polystyrene microbeads with small granulometry for enhanced comfort and adaptability of support to the patient's body.

Synthetic polyester fibre

Strong and resistant to stretching and shrinking. Retains air volume providing a soft sensation and increased comfort.


1 year warranty



Hip Abduction Wedge



Delta Cushion

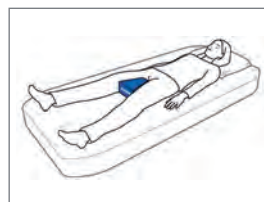

Product Information

Product Code	Description	Dimensions
P9711B1HW	Hip Abduction Wedge	300(L)x100-300(W)x150(H)mm
P9710B1HW	Delta Cushion	600(L)x400(W)x300(H)mm

4. NURSE CALL ALARM SYSTEM

The prevention and management of falls is an important issue in maintaining a patient's quality of life in a care home environment. Harvest Healthcare produce a variety of products to assist in the likelihood of falls and sound an alert on patient movement. Our advanced alert system can be used in conjunction with an existing nurse call system and integrated into the framework of our Woburn bed range.

If you're not sure which wired nurse call system suits your needs please see our fall reduction flowchart guidance (page 91).

NURSE CALL ALARM WIRED SYSTEM



Harvest Healthcare alarm systems are available for use with both seating and mattresses, and can also be connected to your current nurse call system. Be alerted immediately, or with a time delay, whenever one of your residents leaves their bed or chair.

Adjustable volume, optional alarm tones, and a silent mode function when connected to your nurse call system, all mean that you can set up your system to work best for you.

The alarm systems are supplied with a charger as standard, but can also be battery operated (batteries not included).

KEY FEATURES & BENEFITS


- Compact
- Easily connected
- Can be used with existing nurse call system
- Adjustable volume and silent mode
- Delayed alarm mode
- Can be battery operated




Nurse Call Floor Pressure Pad



Seat Pressure Pad

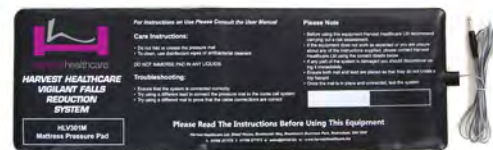

Bed Alarm System

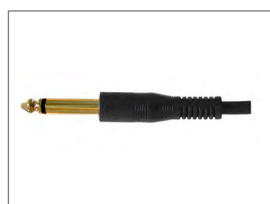

Mono Lead
HJM308



*Mono to Mono Lead
HJM309



**Y connector
HLV006


*Mono to Mono Lead, to be used in conjunction with seat pressure pad, mattress pressure pad and control box
** Y connector, to be used if the floor pad, mattress pad and seat pads are to be used together with 1 Nurse Call adapter in the room

Product Information

Product Code	Description
HJM301	Nurse Call Floor Pressure Pad
HJM302	Sensor Mat Carpet Cover (750 x 1000mm) (for the above floor pressure pad)
HJR075	Alarm System Control Box (to be used with the below mattress & seat pressure pads)
HLV301M	Alarm System Mattress Pressure Pad
HLV301S	Alarm System Seat Pressure Pad
HJM330.M	Nurse Call Handset with Mono Jack Plug (1.8m lead)
HJM330.S	Nurse Call Handset with Stereo Jack Plug (1.8m lead)

NURSE CALL ALARM FLOW CHART

If you're not sure which nurse call system suits your needs please see our fall reduction flowchart guidance.

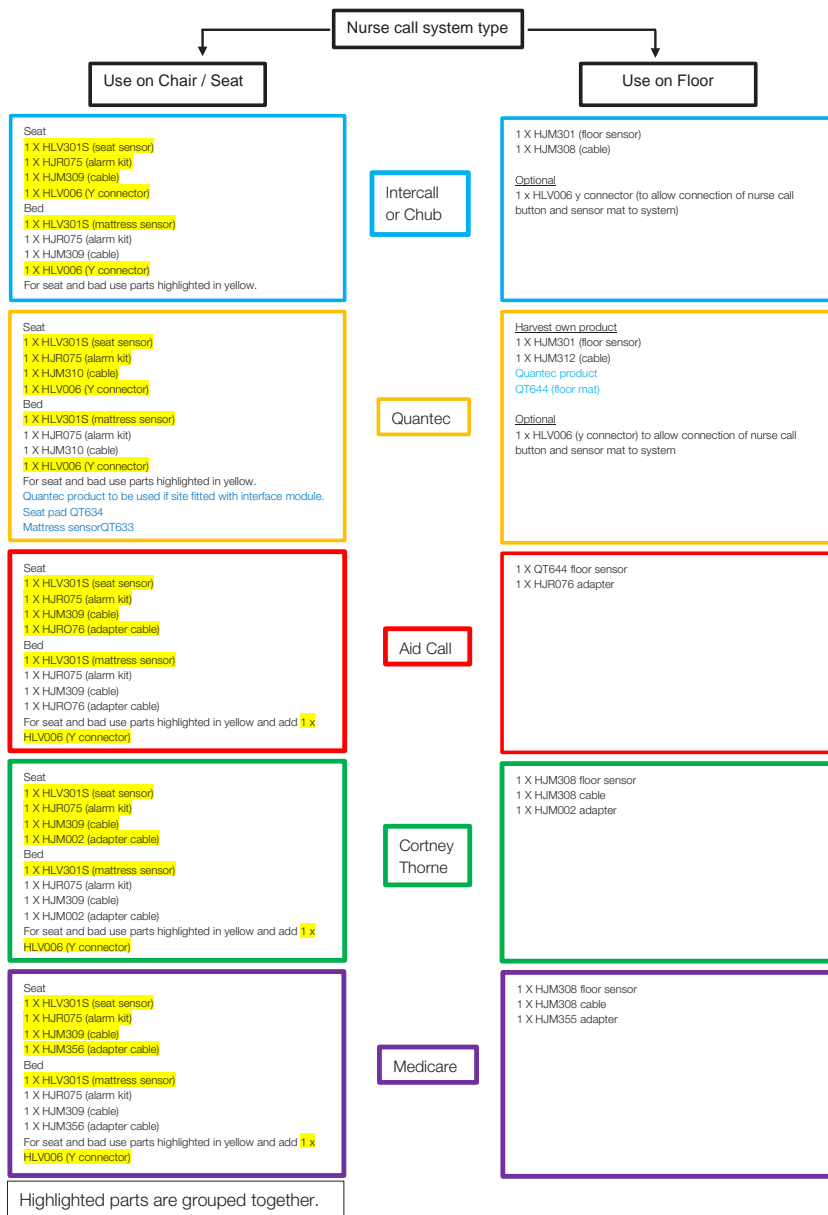

To find out more information about the different nurse call system connections please visit our website:

<https://harvesthealthcare.co.uk/product-category/fall-reduction/>

If your nurse call system is not listed above or if there is anything you're unsure about, please contact us on 01709 377172

Product Information

Product Code	Description
HJM002	Adapter Cable
HJM305	Infrared Bedside Motion Detector
HJM115	1/4" Mono Plug to Two Mono 1/4" Sockets Adapter
HJM310	Nurse Call Connection Lead 3m (Mono to Stereo 6.3mm Plug)
HJM312	Mat Lead 3m (6.3 Male Stereo to 3.5 Mono Male)

NURSE CALL ALARM WIRELESS SYSTEM



The Guardian is a modern wireless falls alarm system. Its fundamental purpose is to act as a pressure-sensitive warning system to alert nurses and carers to unexpected and potentially dangerous movements from users.

Triggered by a release or applied pressure to a surface, the Guardian can be configured to suit an extensive range of situations and requirements.

The package includes 1 Receiver, 1 transmitters and a power supply. The Guardian can also be used with Harvest wired sensor mats HLV301S and HLV301M . Additional components can be purchased separately.

KEY FEATURES & BENEFITS

- Safe and versatile wireless system
- Compatible with any nurse call system*
- Simple to set up and use
- 6 months plus battery life
- Pair up to 8 transmitters
- Can be paired with a variety of different pressure mats to suit various different situations
- Features adjustable alarm times to prevent accidental triggers



Guardian Receiver



Guardian Transmitters

*Different connections are required for different nurse call systems, please ask for more information.

Product Information

Product Code	Description
HJR079	Power Supply
HJM340	Receiver
HJM341	Transmitter
HJM343S	Seat Mat with Pocket
HJM343M	Mattress Mat with Pocket
HJM301	Floor Mat
HJM342	Adapter Cable to use HJM301
HJR078	Adapter Cable to use HLV301S and HLV301M

5. MOVING AND HANDLING

Harvest Healthcare offer a comprehensive range of moving and handling equipment which consists of hoists, stand aids and slings and bed management systems. Suited to a wide range of care environments the equipment provide excellent solutions for moving and handling problems. Promoting mobility and allowing daily handling routines effectively, the products are designed with the comfort, safety and ergonomics of both caregiver and user in mind.

In this section you will find:

Stand Aids, Hoists and Transfer Aids

Slings

Bed Management Systems

Wheelchairs

OLYMPIC PATIENT TURNER



Harvest Healthcare's Olympic Turner is designed with the comfort, safety and ergonomics of both caregiver and patient in mind.

Designed for patients that have some capacity to pull themselves up but need a helping hand, the Olympic Turner requires patients to use their own strength and muscles, helping them to stay active.

KEY FEATURES & BENEFITS

- Helps to assist a user with sit to stand transfers
- Lightweight and easy to manoeuvre
- Padded leg support for added comfort
- Multiple grip options for both user and carer
- Freeway Raiser Strap available to enable a carer to assist a user to their feet
- Compatible with the Freeway safety and ladder belts
- Independent locking rear brakes
- 360° turning radius

1 year warranty



160kg/25st
max user weight



2002ACCRAISLAD	Freeway Ladder Belt
2002ACCRAISBELTS	Freeway Safety Belt - Small
2002ACCRAISBELTM	Freeway Safety Belt - Medium
2002ACCRAISBELTL	Freeway Safety Belt - Large



Freeway Safety Belt



Freeway Ladder Belt
590(W) x 120(D)mm

Product Information

Product Code	Overall Height	Assembled Unit Weight	Front and Rear Castors/Wheels
HLA045	1201mm	19kg	35mm/125mm

OLYMPIC PATIENT MOVER



Simple and intuitive to use, the Olympic Patient Mover is an attractive, quick and universal resident transport solution. Using the Olympic Patient Mover allows resident transfers with minimal handling and with more efficiency than other transfer options. The Olympic Patient Mover supports residents to move comfortably from a seated to a standing position, and promotes early mobilisation.

KEY FEATURES & BENEFITS

- Provides easy and more efficient transportation to the bathroom and toilet
- Offers a suitable alternative to a wheelchair
- Allows for better eye contact between the user and carer
- Padded knee support for extra comfort
- Manual leg open & closing
- Compatible with the Mackworth back support belt


1 year warranty



185kg/29st
max user weight



1201TAB550C	Mackworth Back Support Belt S
1201TAB540C	Mackworth Back Support Belt M
1201TAB530C	Mackworth Back Support Belt L
1201TAB520C	Mackworth Back Support Belt XL

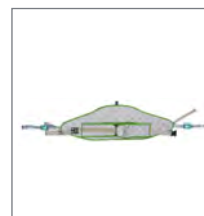

Mackworth
Back Support Belt

Product Information

Product Code	Length	Width	External Width (Legs Open)	Turning Radius
HLA046	1500mm	650mm	1060mm	1125mm

ASCENT 180 STAND AID



The Ascent 180 Stand Aid has an impressive range of ergonomic features to provide a safe and dignified lift.

It is uniquely designed to provide three functions; as a transfer hoist, as a sit to stand hoist and as an aid to rehabilitation by assisted standing.

Suitable for transferring patients into a wheelchair, onto toileting facilities, onto a bed, or for many other assisted transfers, the Ascent 180 Stand Aid is also excellent for standing exercises and promoting mobility.

KEY FEATURES & BENEFITS

- Ability to lift a wide range of people from 4ft – 6ft 5”
- Compact and easily manoeuvred
- Profiled push handles designed in conjunction with the REBA assessment tool allows for safe moving and positioning
- Powered leg opening allows for a stable and safer lift without restricting access around chairs and beds
- Front castors and rear castors for maximum control and positioning on all floor surfaces
- Multi-directional knee pad allows for individual adjustment to provide full support
- Rear foot bar helps to negotiate the stand aid over uneven floor surfaces and door thresholds

5 year warranty

180kg/28st
max user weight



OPTIONAL EXTRAS:

- Jumbo Charger Complete with Wall Mount Bracket & Lead: [ELE022](#)
- Jumbo Battery Pack: [ELE053A](#)



For Clip Stand Aid Slings
HLA039-Powered Leg



For Loop Stand Aid Slings
HLA042-Powered Leg

Product Information

Product Code	Length	Height (Lowest/Highest)	Leg Inside Width (Closed/Spread)	Lifting Range
HLA039-Clip with Powered Leg	1190mm	1145mm/1720mm	540mm/1010mm	885mm
HLA042-Loop with Powered Leg				

ASCENT PRO 200 PREMIUM STAND AID



harvesthealthcare
a prismhealthcare company

With a sleek, ergonomic design, the Ascent Pro 200 provides ideal assistance for transfers while maintaining a home from home aesthetic.

Suitable for transferring patients into a wheelchair, onto toileting facilities, onto a bed, or for many other assisted transfers, the Ascent Pro 200 is also excellent for standing exercises and promoting mobility.

KEY FEATURES & BENEFITS

- Premium, sleek, aesthetic design
- Battery operated height and leg adjustment
- Soft stop/start motion
- Handset operation with battery level display
- LCD display for service notifications, actuator cycles and overloads
- Ergonomic handles
- Manual emergency lowering system
- Flexible extra soft knee supports with safety straps
- Removable non-slip footplate
- Enclosed castors keep out dirt and are easy to maintain

5 year warranty



200kg/31.5st
max user weight



OPTIONAL EXTRAS:

- Spare Battery: [ELE053A](#)
- Jumbo Charger Complete with Wall Mount Bracket & Lead: [ELE022](#)

Product Information

Product Code	Length	External Width (Closed/Spread)	Adjustable Height
HLA006	1027mm	700mm/1044mm	1168mm-1781mm

OLYMPIC 180 HOIST



The Olympic 180 Hoist benefits from the rigorous ergonomic design procedures and quality checks. It has been constructed from high-quality powder coated steel that will provide years of reliable service. The hoist boasts an impressive 180kg lifting capacity, powered leg opening and market leading Linak® actuators and motors.

The Hoist Actuator Pad is also available as an optional extra. The Actuator Pad provides a protective layer to minimise any risk of issue surrounding body parts coming in to contact with hoist actuators, without losing any functionality to the hoist.

KEY FEATURES & BENEFITS

- Cost-effective solution without compromising on quality
- Robust powder coated steel frame
- Linak control system
- Electric leg opening as standard
- Red emergency lowering screw
- Easy to position with 4 large swivel castors and ergonomic handle
- Tool free assembly
- Single person installation
- Can lift a patient from the floor
- 75mm front castors and 100mm rear castors
- Available with a scales unit incorporating QRS feature which is a dignified, accurate and convenient method of determining body mass by weighing an individual whilst situated in a sling and suspended by the Olympic 180 hoist

5 year warranty



185kg/29st
max user weight



OPTIONAL EXTRAS:

- Jumbo Charger Complete with Wall Mount Bracket & Lead: [ELE022](#)
- Jumbo Battery Pack: [ELE053A](#)
- Hoist with QRS Feature: [HLA034QRS](#)
- Pivot Frame for Clip Slings Available
- Hoist Actuator Pad Available



Hoist Actuator Pad



QRS

Product Information

Product Code	Length	Height (Highest/Lowest)	Leg Inside Width (Closed/Open)	Lifting Range
HLA034 HLA034QRS	1265	1335mm/2010mm	490mm/1075mm	1755mm

OLYMPIC PRO 185 PREMIUM HOIST



The Olympic Pro 185 hoist provides a beautifully designed, premium solution for labour-intensive resident transfers, allowing powered, single carer use. The powered dynamic positioning system (PDPS) provides precise and stable positioning without the need for manual intervention.

The height, leg spreader and carrier frame are all electronically controlled, and the anti-rotation mechanism prevents unintentional rotational movement during lifting and transferring. The carrier frame can be lowered to the floor to allow picking up residents from the ground. Suited to a wide range of care environments.

The Hoist Actuator Pad is also available as an optional extra. The Actuator Pad provides a protective layer to minimise any risk of issue surrounding body parts coming in to contact with hoist actuators, without losing any functionality to the hoist.

KEY FEATURES & BENEFITS

- Aesthetic design, extra lightweight
- Emergency stop and emergency lowering
- Carrier frame and lifting column made of aluminium and steel
- Enclosed lockable double castors
- Electrically adjustable height, cradle and spreadable carrier frame
- Electrical Positioning System (EPS)
- Anti-rotation mechanism
- Battery compatible with Ascent Pro 200 stand aid
- Slow decline feature, in case of overloading
- LOLER test reminder light

5 year warranty



185kg/29st
max user weight



OPTIONAL EXTRAS:

- Spare Battery: [ELE053A](#)
- Jumbo Charger Complete with Wall Mount Bracket & Lead: [ELE022](#)
- Electronic Weighing Scale: [HLA010](#)
- Hoist with Weigh Scale: [HLA020](#)
- Two Point Yoke Frame for Loop Slings Available
- Hoist Actuator Pad Available



Hoist Actuator Pad



Hoist with the Actuator Pad

Product Information

Product Code	Chassis Length/Height	Chassis Width (Closed/Spread)	Maximum Height	Lifting Range
HLA002 HLA020 with Weigh Scale	1290mm/125mm	700mm/1180mm	2170mm	1400mm

OLYMPIC EXTRA 230 PLUS SIZE HOIST



The Olympic Extra 230 is a sling hoist specially made for lifting, moving and transferring plus size patients. Due to the special construction and the steel components, the Olympic Extra is extra light and easy to move but still offers the required stability for the care of plus size patients.

We provide compatible slings for our hoists.

The Hoist Actuator Pad is also available as an optional extra. The Actuator Pad provides a protective layer to minimise any risk of issue surrounding body parts coming in to contact with hoist actuators, without losing any functionality to the hoist.

KEY FEATURES & BENEFITS

- Aesthetic design
- Universal use, also appropriate for bathing
- Extra light and easy to manoeuvre
- Electrical positioning system (EPS) places the patient exactly without additional force
- Electrically adjustable height, cradle and spreadable carrier frame, for easier patient transfer
- Lifts from the floor
- No rotational movements during lifting and transferring due to anti-rotation mechanism
- Hand control with battery charge status display
- In-built electronic scales (optional extra)
- Wall charging station (optional extra)

5 year warranty



230kg/36st
max user weight



OPTIONAL EXTRAS:

- Jumbo Charger Complete with Wall Mount Bracket & Lead: [ELE022](#)
- Jumbo Battery Pack: [ELE053A](#)



Pivot frame
To be used with clip slings



2 point yoke frame
To be used with loop slings

Product Information

Product Code	Length	External Width (Legs Open)	Min Height / Max Height
HLA007	1400mm	1240mm	1340mm / 2170mm

FREEWAY SCALES UNIT INCORPORATING QRS

freeway

A dignified, accurate and convenient method of determining body mass by weighing an individual whilst situated in a sling and suspended by a hoist.

This product is compatible with the Olympic 180 hoist (product code: HLA034QRS) .


KEY FEATURES & BENEFITS

- Fits in-line between the hoist and carry bar
- Load of 10 to 350kg (functional maximum load is limited by the hoist or sling, i.e. 180kg when used with the Olympic 180 hoist).
- QRS (Quick Release System)
- 1 button operation
- Auto shut down after 60 seconds (programmable)
- 15mm LED screen
- Display resolution to 0.2kg
- Class III accuracy

5 year warranty



180kg/28st
max user weight

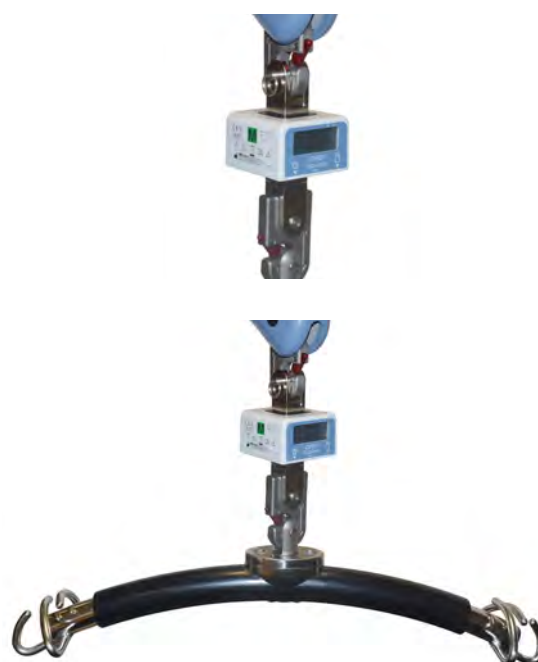

Product Information

Product Code	Description
610000	QRS & Scales Assembly
1102CBQRS3HK	3 Hooked Clipped Carry Bar

FREEWAY TRANSACTIVE XTRA CEILING TRACK HOIST

freeway


The Freeway Transactive Xtra is a product of world-class design and state of the art features, providing safe, reliable and effective hoisting for carers and users alike.

The Transactive Xtra benefits from a low profile design which is unobtrusive, yet provides a much greater lifting height than any other hoist option. Together with the extensive range of track configurations, the Transactive Xtra can provide the perfect solution for any hoisting requirements.

KEY FEATURES & BENEFITS

- Ergonomically designed handset
- State-of-the-art display & diagnostics
- Freeway Spreader Bar as standard
- Emergency stop
- Powered emergency lowering
- Manual lowering & raising facility

2 year warranty

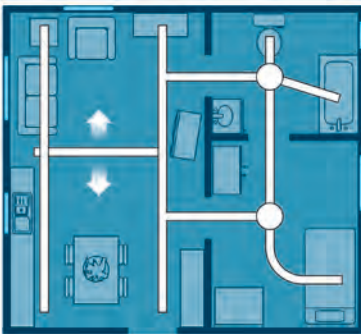

Installed Ceiling Track Systems – Multiple options, subject to site survey

Straight Track – Point to Point transfer

Room Covering – Pick up and transfer to multiple points within the room

Room to Room – Transfers between multiple linked rooms

Product Information

Product Code	Description
1001TX270D	Transactive Xtra Drive - 270kg
1001TX270	Transactive Xtra Manual - 270kg
1001TX200D	Transactive Xtra Drive - 200kg
1001TX200	Transactive Xtra Manual - 200kg
1001TX160D	Transactive Xtra Drive - 160kg
1001TX160	Freeway Transactive Xtra Manual - 160kg
1001TX130D	Freeway Transactive Xtra Drive - 130kg
1001TX130	Freeway Transactive Xtra Manual - 130kg

UNIVERSAL DELUXE SLING



This multi-purpose sling is designed to offer excellent support, and to be used for the majority of patients, including for bathing when using super soft material.

Available in polyester which is versatile, breathable and durable.

It is also available in super soft 3-dimensional fabric which is breathable, very durable and provides excellent pressure distribution. Super soft fabric is also fast drying and therefore suitable for bathing.

KEY FEATURES & BENEFITS

- Can be applied and removed while seated in a chair
- Full padded support of upper thigh reduces hip rotation and abduction
- Available in a range of sizes
- Available in Polyester and Super Soft material
- Pocket in the back of the sling to aid proper fitting
- Available with head support
- Also available with clip fastening
- Extra stitching to improve life and safe working load of the sling
- Information stitched into label to protect when washing




300kg/47st (super soft)
max user weight



400kg/62st (Polyester)
max user weight


1 year warranty

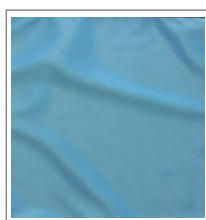

Available in
Polyester Fabric



Available in
Super Soft Fabric

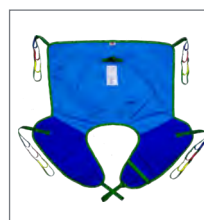

Universal Deluxe
Loop Sling
HSL100



Universal Deluxe
Clip Sling
HSL110

Product Information

Product Code	Sizes	Width Range (cm)	Back Length Range (cm)
HSL100/HSL110 (Please see page 107 & 108)	X Small Small Medium Large X Large	45-52 52-60 60-76 76-92 92-100	48-56 53-63 63-73 68-78 68-78

IN-CHAIR HAMMOCK SLING



harvesthealthcare
a prismhealthcare company

The In Chair Hammock Sling is designed to provide full support to the back, neck and shoulders.

Available in polyester which is versatile, breathable and durable.

It is also available in super soft 3-dimensional fabric which is breathable, very durable and provides excellent pressure distribution. Super soft fabric is also fast drying and therefore suitable for bathing.

KEY FEATURES & BENEFITS

- Hip Straps
- Pommel
- Pockets to tuck lifting straps in
- Available in a range of sizes
- Available in Polyester and Super Soft materials
- Also available with clip fastening
- Extra stitching to improve life and safe working load of the sling
- Information stitched into label to protect when washing



300kg/47st (super soft)
max user weight



400kg/62st (Polyester)
max user weight

1 year warranty



Available in
Polyester Fabric



Available in
Super Soft Fabric



In Chair Hammock
Loop Sling
HSL400

Product Information

Product Code	Sizes	Width Range (cm)	Back Length Range (cm)
HSL400/HSL410 (Please see page 108)	X Small	45-52	48-56
	Small	52-60	53-63
	Medium	60-76	57-67
	Large	76-92	62-72
	X Large	92-100	65-75

TOILET / DRESS SLING




The Toilet/Dress Sling has been designed to help facilitate the toileting process, with good access to allow for the removal of clothing. Available in polyester which makes it versatile, breathable and durable.

KEY FEATURES & BENEFITS

- Easy to fit
- Fully padded
- Adjustable Hook and Loop body belt
- Wide legs for additional comfort
- Slide sheet material on legs to reduce friction and ease fitting
- Available in a range of sizes
- Also available with clip fastening
- Extra stitching to improve life and safe working load of the sling
- Information stitched into label to protect when washing

1 year warranty

400kg/62st
max user weight

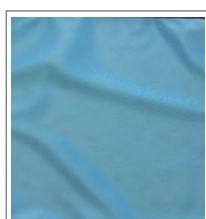

Available in
Polyester Fabric



Toilet/Dress
Sling
HSL200

Product Information

Product Code	Sizes	Width Range (cm)	Back Length Range (cm)
HSL200 (Please see page 109)	X Small	40-50	26-30
	Small	50-60	29-34
	Medium	62-68	34-40
	Large	68-78	36-42
	X Large	78-85	38-44

STAND AID SLING



harvesthealthcare
a prismhealthcare company

The Stand Aid Sling is designed to assist with the standing process using a hoist. It can be used with Harvest Healthcare Ascent Extra Stand Aid.

Available in polyester which makes it versatile, breathable and durable.


KEY FEATURES & BENEFITS

- Easy to fit
- Fully padded
- Adjustable Hook and Loop body belt
- Available in a range of sizes
- Also available with clip fastening
- Extra stitching to improve life and safe working load of the sling
- Information stitched into label to protect when washing


1 year warranty



400kg/62st
max user weight

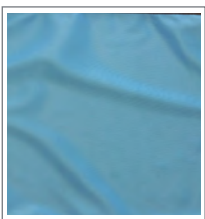

Available in
Polyester Fabric

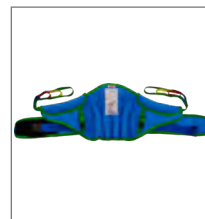

Stand Aid Loop
Sling
HSL500



Stand Aid Clip
Sling
HSL510

Product Information

Product Code	Sizes	Width Range (cm)	Back Length Range (cm)
HSL500/HSL510 (Please see page 109)	X Small	40-50	28-32
	Small	50-60	30-36
	Medium	60-70	33-39
	Large	70-80	37-42
	X Large	80-90	40-46

SLING CODES



Universal Deluxe Sling

Product Code	Description	Material	Size
HSL100.PO.S	Universal Deluxe Sling with Loop	Polyester	Small
HSL100.PO.M	Universal Deluxe Sling with Loop	Polyester	Medium
HSL100.PO.L	Universal Deluxe Sling with Loop	Polyester	Large
HSL100.PO.XL	Universal Deluxe Sling with Loop	Polyester	Extra Large
HSL100.SS.S	Universal Deluxe Sling with Loop	Super Soft	Small
HSL100.SS.M	Universal Deluxe Sling with Loop	Super Soft	Medium
HSL100.SS.L	Universal Deluxe Sling with Loop	Super Soft	Large
HSL100.SS.XL	Universal Deluxe Sling with Loop	Super Soft	Extra Large
HSL110.PO.S	Universal Deluxe Sling with Clip	Polyester	Small
HSL110.PO.M	Universal Deluxe Sling with Clip	Polyester	Medium
HSL110.PO.L	Universal Deluxe Sling with Clip	Polyester	Large
HSL110.PO.XL	Universal Deluxe Sling with Clip	Polyester	Extra Large
HSL110.SS.S	Universal Deluxe Sling with Clip	Super Soft	Small
HSL110.SS.M	Universal Deluxe Sling with Clip	Super Soft	Medium
HSL110.SS.L	Universal Deluxe Sling with Clip	Super Soft	Large
HSL110.SS.XL	Universal Deluxe Sling with Clip	Super Soft	Extra Large
HSL300.PO.S	Universal Deluxe Sling with Head Support & Loop	Polyester	Small
HSL300.PO.M	Universal Deluxe Sling with Head Support & Loop	Polyester	Medium
HSL300.PO.L	Universal Deluxe Sling with Head Support & Loop	Polyester	Large
HSL300.PO.XL	Universal Deluxe Sling with Head Support & Loop	Polyester	Extra Large
HSL300.SS.S	Universal Deluxe Sling with Head Support & Loop	Super Soft	Small
HSL300.SS.M	Universal Deluxe Sling with Head Support & Loop	Super Soft	Medium
HSL300.SS.L	Universal Deluxe Sling with Head Support & Loop	Super Soft	Large
HSL300.SS.XL	Universal Deluxe Sling with Head Support & Loop	Super Soft	Extra Large

SLING CODES

**Universal Deluxe Sling**

Product Code	Description	Material	Size
HSL310.PO.S	Universal Deluxe Sling with Head Support & Clip	Polyester	Small
HSL310.PO.M	Universal Deluxe Sling with Head Support & Clip	Polyester	Medium
HSL310.PO.L	Universal Deluxe Sling with Head Support & Clip	Polyester	Large
HSL310.PO.XL	Universal Deluxe Sling with Head Support & Clip	Polyester	Extra Large
HSL310.SS.S	Universal Deluxe Sling with Head Support & Clip	Super Soft	Small
HSL310.SS.M	Universal Deluxe Sling with Head Support & Clip	Super Soft	Medium
HSL310.SS.L	Universal Deluxe Sling with Head Support & Clip	Super Soft	Large
HSL310.SS.XL	Universal Deluxe Sling with Head Support & Clip	Super Soft	Extra Large

In Chair Hammock Sling

Product Code	Description	Material	Size
HSL400.PO.S	In Chair Hammock Sling with Loop	Polyester	Small
HSL400.PO.M	In Chair Hammock Sling with Loop	Polyester	Medium
HSL400.PO.L	In Chair Hammock Sling with Loop	Polyester	Large
HSL400.PO.XL	In Chair Hammock Sling with Loop	Polyester	Extra Large
HSL400.SS.S	In Chair Hammock Sling with Loop	Super Soft	Small
HSL400.SS.M	In Chair Hammock Sling with Loop	Super Soft	Medium
HSL400.SS.L	In Chair Hammock Sling with Loop	Super Soft	Large
HSL400.SS.XL	In Chair Hammock Sling with Loop	Super Soft	Extra Large
HSL410.PO.S	In Chair Hammock Sling with Clip	Polyester	Small
HSL410.PO.M	In Chair Hammock Sling with Clip	Polyester	Medium
HSL410.PO.L	In Chair Hammock Sling with Clip	Polyester	Large
HSL410.PO.XL	In Chair Hammock Sling with Clip	Polyester	Extra Large
HSL410.SS.S	In Chair Hammock Sling with Clip	Super Soft	Small
HSL410.SS.M	In Chair Hammock Sling with Clip	Super Soft	Medium
HSL410.SS.L	In Chair Hammock Sling with Clip	Super Soft	Large
HSL410.SS.XL	In Chair Hammock Sling with Clip	Super Soft	Extra Large

SLING CODES



Toilet/Dress Sling


Product Code	Description	Material	Size
HSL200.PO.S	Toilet/Dress Sling with Loop	Polyester	Small
HSL200.PO.M	Toilet/Dress Sling with Loop	Polyester	Medium
HSL200.PO.L	Toilet/Dress Sling with Loop	Polyester	Large
HSL200.PO.XL	Toilet/Dress Sling with Loop	Polyester	Extra Large

Stand Aid Sling

Product Code	Description	Material	Size
HSL500.PO.S	Stand Aid Sling with Loop	Polyester	Small
HSL500.PO.M	Stand Aid Sling with Loop	Polyester	Medium
HSL500.PO.L	Stand Aid Sling with Loop	Polyester	Large
HSL500.PO.XL	Stand Aid Sling with Loop	Polyester	Extra Large
HSL510.PO.S	Stand Aid Sling with Clip	Polyester	Small
HSL510.PO.M	Stand Aid Sling with Clip	Polyester	Medium
HSL510.PO.L	Stand Aid Sling with Clip	Polyester	Large
HSL510.PO.XL	Stand Aid Sling with Clip	Polyester	Extra Large

Sling Sizing Guide

1. Measure the width of your resident between the outer midpoints of each upper arm as per the diagram (W).
2. Compare your measurement with the 'User Width' range for the selected sling type in the table below.
3. Measure the back length of your resident as per the diagram (B1/B2).
4. Compare your measurement with the 'User Back Length' range for the selected sling type in the table below. Please note some ranges overlap, in these cases select the most appropriate size.
5. If these two measurements correlate, your resident fits a standard sling. If not, consider a bespoke sling or select the most appropriate standard sling after a risk assessment.

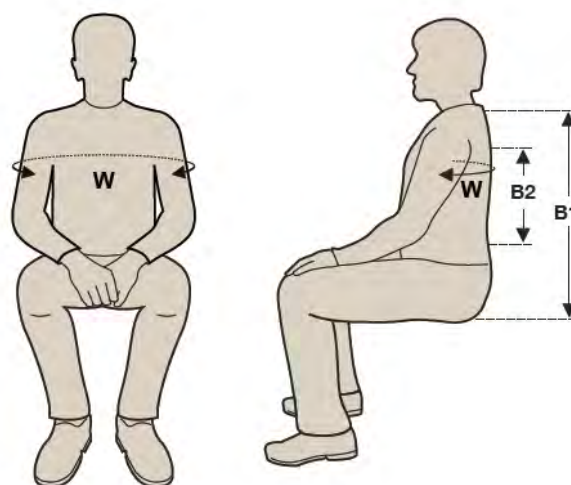

All sling selections should be accompanied by a risk assessment by a trained person. People of atypical anatomy may not be suited to a standard sling.

The Harvest Tape Measure (HSLTAPE) and Sling Assessment Guides are available on request.

FREEWAY FALLS RECOVERY SLING

freeway

The Freeway Falls Recovery Sling allows carers to safely transfer a person that needs assistance immediately after a fall. It has a versatile, one size fits all design, so is suitable for both adults and children.

The Freeway Falls Recovery Sling allows for the fallen person to be rolled onto the sling whilst in a lying position. Once moved onto the sling, the person can be easily slid to a more open area if required, for the suitable positioning of a hoist. The person can then be hoisted up off the floor and safely into a chair or bed.

KEY FEATURES & BENEFITS

- Versatile design – suitable for adults and children
- Strong, yet built for comfort
- Low friction material for comfortable transfer and sliding of a person (the yellow side)
- Easy application and removal
- Integral space saving storage and carry bag
- Machine washable



200kg/31st
max user weight



Product Information

Product Code	Description	Material	Size
12018E3-FRS	Falls Recovery Sling	Polyester /Slipfit	One Size

FREEWAY FALLS RECOVERY SLING

freeway



MACKWORTH ESSENTIALS SATIN SHEET

mackworth

The Mackworth Essentials Satin Sheet range is designed to reduce lifting associated strains and injuries to caregivers, whilst providing a solution for patients who require frequent repositioning.

The system consists of two components; a synthetic, breathable satin base sheet with a low friction central section and a soft top sheet with low friction backing. The Mackworth Essentials Satin Sheet range features handles on the Top Sheet to facilitate easier movement and use in conjunction with hoists.



KEY FEATURES & BENEFITS

Base Sheet

- Available in standard and bariatric mattress sizes
- Made of synthetic breathable satin and cotton

Top Sheet

- Maxi (full bed) and Midi (half bed) sizes
- Handles for easier positioning and use with hoists
- Made of synthetic breathable satin and cotton

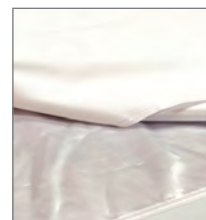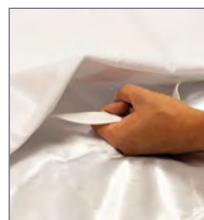

Product Information

Product Code	Description	Dimensions to fit mattress size
2002MAC9020	Mackworth Standard Base Sheet	2000(L)x900(W)mm
2002MAC1420	Mackworth Midi Top Sheet	2000(L)x1400(W)mm
2002MAC12020	Mackworth Satin Sheets Bariatric Base	2000(L)x1200(W)mm
2002MAC2020	Mackworth Satin Sheets Maxi Top	2000(L)x2000(W)mm

FREWAY BED MANAGEMENT SYSTEM

FREWAY BODY SUPPORT WEDGE



The **Freeway Bed Management System** (bariatric) provides an easy to use, low friction solution for patients who require frequent repositioning. It consists of two separate units – a nylon, water resistant under sheet with wide non-slip sides, and a fully breathable, water resistant, quilted polyester softshell top sheet. The system has been designed to provide comfort for patients, and can help to reduce incidences of pressure ulcers. It also provides ease of use by carers, and forms an ideal bed solution as part of a single handed care package.

Also available to work in conjunction with the Freeway Bed Management System, is the **Freeway Support Wedge**. This provides firm support for a patient in bed and helps them remain in a required position.

KEY FEATURES & BENEFITS

Bed Management System

- Provides comfort for patients
- Minimal effort from carers
- Ideal for use by a single carer
- Works easily with four section profiling beds
- Can be used on a pressure care mattress
- Can also be used with a hoist to assist turning

Body Support Wedge

- Provides a positive non-slip action when used against cotton bed linen
- Can be used to help reduce the incidence of pressure ulcers
- Made from 100% soft and flame retardant foam
- The polyurethane cover can be machine washed up to 80°C
- The nylon cover can be wiped clean with disinfectant



Product Information

Product Code	Product Description	Dimensions
2002ACCBMS-COM	Freeway Bed Management System - Complete - Single Bed Size	-
2002ACCBMS-BSS	Freeway Bed Management System - Base Sheet Only	2000(L)x900(W)mm
2002ACCBMS-QS	Freeway Bed Management System - Top Sheet Only	2000(L)x1400(W)mm
2002ACCBMS-BCOM	Freeway Bed Management System - Complete - Plus Size	-
2002ACCBMS-BSB	Freeway Bed Management System - Base Sheet Only - Plus Size	2000(L)x1200(W)mm
2002ACCBMS-QB	Freeway Bed Management System - Top Sheet only	2000(L)x1900(W)mm
200FW-WEB53	Freeway Positioning Wedge	53mm(W)
200FW-WEB60B	Freeway Positioning Wedge - Plus Size	60mm(W)

MULTIGLOVE LIMB LIFTER AND SLIDE SHEETS



The **Prism Multiglove** enables carers to facilitate a variety of moving and handling applications with ease due to its low friction. It is a simple but effective piece of equipment, and can be used to bring a person to a sitting position in a profiling bed, as the anti-shear material helps to smooth and release the tension on the skin.

Product Information

Product Code	Description
2002ACCPM-MGLOVE	Prism Multiglove (Pair)



The **Prism Limb Lifter** has been designed to support limbs, allowing a number of procedures to be completed, i.e. wound care, without the carer having to support the limb. The Prism Limb Lifter can be easily attached to a hoist, enabling the carer to perform procedures in the safest possible working environment, eliminating the possibility of strain or injury from supporting the limb themselves.

Coloured inset loops are featured to aid correct positioning of the Prism Limb Lifter when attaching to a hoist. Each Prism Limb Lifter is supplied with two removable washable covers. Additional washable covers can be purchased here.

The product is machine washable up to 72°C.

Product Information

Product Code	Description	Dimensions
2002ACCPLL150X1	Prism Limb Lifter	150mm(L)
2002ACCPLL150X5	Prism Limb Lifter	150mm(L) (box of 5)



The **Tubular and Flat Slide Sheets** are used as an aid for repositioning and moving clients. Various sizes available.

Product Information

Product Code	Description
HTS201	Tubular Slide Sheet (720 x 700mm)
HTS202	Tubular Slide Sheet (1100 x 800mm)
HTS204	Tubular Slide Sheet (1350 x 1000mm)
HTS267	Tubular Slide Sheet (1950 x 1400mm)
HTS207	Flat Slide Sheet (2000 x 1000mm)
HTS205	Flat Slide Sheet w/ Handles (2000 x 710mm)
HTS206	Flat Slide Sheet w/ Handles (2000 x 1000mm)



PRISM ONE WAY GLIDE PRISM HANDLING BELT



The **Prism One Way Glide Sheets** help a patient maintain a good sitting position and allow repositioning on a chair or bed.

The one-way glide will only slide in the direction of the arrows which are clearly labelled on the sheets. It can be used on all surfaces and helps prevent a patient from sliding forward in a wheelchair/day chair/bed and aids the positioning of a patient into a good sitting position.




Product Information

Product Code	Description
2002ACCPM-OWGS14	Prism One Way Glide - Single 1400 x 430
2002ACCPM-OWGS70	Prism One Way Glide - Single 700 x 430
2002ACCPM-OWGT60	Prism One Way Glide - Tube 600 x 430
2002ACCPM-OWGT73	Prism One Way Glide - Tube 730 x 600

The **Prism Comfi Belt** helps carers perform assisted manoeuvres and transfers safely whilst minimizing the risk of injury to themselves or the user. It has been designed to provide high levels of comfort to the user whilst providing the best possible grip and handling positions for the carer. Available in three separate sizes, each of which is also adjustable to ensure the perfect fit for waist sizes from 22inch (558mm) to 55inch (1395mm).

The product is machine washable up to 72°C.

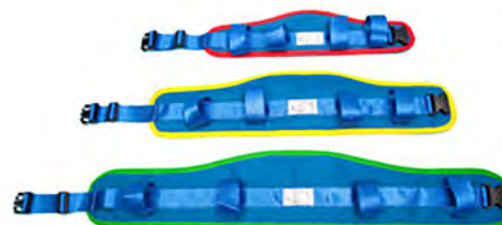

Product Information

Product Code	Description
2002ACCPM-CBLARG	Prism Comfi Belt - Large (1010 - 1410 mm)
2002ACCPM-CBMED	Prism Comfi Belt - Medium (710 - 1150 mm)
2002ACCPM-CBSMAL	Prism Comfi Belt - Small (560 - 810 mm)

SELF PROPEL WHEELCHAIR



The Self Propel Wheelchair is a cost-effective wheelchair, which has been designed around comfort and practicality.

Built with a strong steel framework, the chair can support a maximum user weight of 110kg.

KEY FEATURES & BENEFITS

- Strong steel framework
- 18 inch wide seat
- Black fire retardant seat material
- PU tyres (no more punctures)
- Rear wheel brakes
- Full length padded (detachable) armrests
- Armrests designed to enable the user to sit at tables
- Removable swing away footrests
- Folding backrest
- Heel support strap
- Optional Anti-wheelie Bars for the rear
- Supplied with instruction booklet and tools

5 year warranty



110kg/17st
max user weight



Product Information

Product Code	Dimensions	Seat	Front Wheel	Rear Wheel	Net Weight
HMS300B	1030 x 650mm	460(Wide)x410(Deep)mm	203mm (8") Diameter	610mm (24") Diameter	18kg

CAR TRANSIT WHEELCHAIR



The Car Transit Wheelchair is a fantastic economy wheelchair, designed around comfort and practicality.

Built with a strong steel framework, the chair is readily available to assist patients and carers when on the move.

KEY FEATURES & BENEFITS

- Strong steel framework
- 18 inch wide seat
- Black fire retardant seat material
- PU tyres (no more punctures)
- Rear wheel brakes
- Full length padded (detachable) armrests
- Armrests designed to enable the user to sit at tables
- Removable swing away footrests
- Folding backrest
- Heel support strap
- Supplied with instruction booklet and tools

5 year warranty



110kg/17st
max user weight



Product Information

Product Code	Dimensions	Seat	Front Wheel	Rear Wheel	Net Weight
HMA200B	920x600mm	460(Wide)x410(Deep)mm	203mm (8") Diameter	305mm (12") Diameter	15kg

6. SHOWERING AND TOILET AIDS

Harvest Healthcare stocks a comprehensive selection of Showering and Toilet Aids, suitable for a wide range of patients.

TILT IN SPACE SHOWER CHAIR



This shower chair provides comfort and dignity for patients during washing and toileting routines, while maintaining ease of use for the carer.

It provides a versatile bathing and toileting solution for a wide variety of patients.

KEY FEATURES & BENEFITS

- Tilt-in-space shower/commode chair
- Adjustable seat angle of 6° to 25° ensures suitability for a wide range of conditions and facilitates patient transfers
- Length adjustable, swing-away footrest for ease of access and for varying leg lengths
- Height adjustable headrest suitable for a range of patient heights
- Flip-up armrests manufactured in strong yet lightweight aluminium
- Durable and easy to transport; will not rust or corrode in wet environments
- Padded seat provides a comfortable and secure seating position for the patient
- 4 individually locking castors provides a stable seating position whilst aiding manoeuvrability
- Removable backrest can be laundered and replaced as necessary

1 year warranty



158kg/25st
max user weight



Dimensions

Overall Width	650mm	Seat Width	450mm	Back Height From Seat	1000mm
Overall Height	1100-1400mm	Seat Height	550mm	Armrest Height From Seat	250mm
Overall Length	1000mm	Seat Depth	420mm	Armrest Width x Armrest Length	50 x 240mm

Product Information

Product Code	Description
HMC119	Aluminium Tilt Transit Shower / Commode Chair

T20 MOBILE SHOWER CHAIR

freeway

The T20 benefits from all the attributes in terms of the toughness of the whole Freeway range and is available at a remarkably low price.

An antimicrobial and antibacterial coating protects the chair. It is resistant to the hospital superbug MRSA, E. coli, salmonella and listeria, providing a durable and hygienic option for toileting and bathing.

KEY FEATURES & BENEFITS

- Perfect for those who need a basic shower chair only
- Available either free-standing or with high-performance castors
- Frame constructed from heavy-duty, durable tubing to provide years of rust-free performance
- Coating of the shower chair frame includes an antibacterial and antimicrobial coating

5 year warranty



200kg/30st
max user weight



Dimensions/Technical Information

Seat Dimensions	440(W)mm x 400(D)mm
Height from floor to top of the seat:	460mm
Overall maximum width:	550mm
Overall maximum height:	900mm

Product Information

Product Code	Description
HMC100	Freeway T20 Mobile Shower Chair

FREEWAY ADJUSTABLE HEIGHT SHOWER TROLLEY

freeway

The Freeway Height Adjustable Shower Trolley benefits from optional cot sides, and its height can be adjusted by the hydraulic foot pump. Available in different sizes of 1400mm, 1500mm, 1700mm and 1900mm, the Freeway Height Adjustable Shower Trolley has a working height of between 600mm and 925mm.

The protective coating of all our shower trolleys and stretchers includes Freecote® - our unique antimicrobial and antibacterial coating, which is resistant against the hospital superbug MRSA, e-coli, salmonella and listeria.

KEY FEATURES & BENEFITS

- Available in 3 sizes – 1500 mm, 1700 mm and 1900 mm
- Has a working height of 600 mm to 925 mm.
- Practical ergonomic design
- Durable, manoeuvrable and extraordinary value for money.
- Top quality stainless steel bearing castors
- Special custom build available
- Includes FREECOTE® – our unique antimicrobial and antibacterial coating

1 year warranty



200kg/30st
max user weight

Product Information

Product Code	Product Description
1801AHT1900	Adjustable Height Trolley - 1900mm
1801AHT1700	Adjustable Height Trolley - 1700mm
1801AHT1500	Adjustable Height Trolley - 1500mm
1801AHT1900B	Adjustable Height Trolley with Adjustable Backrest - 1900mm
1801AHT1700B	Adjustable Height Trolley with Adjustable Backrest - 1700mm
1801AHT1500B	Adjustable Height Trolley with Adjustable Backrest - 1500mm
1801AHTPC1900	Adjustable Height Padded Changing Trolley - 1900mm
1801AHTPC1700	Adjustable Height Padded Changing Trolley - 1700mm
1801AHTPC1500	Adjustable Height Padded Changing Trolley - 1500mm
1801AHT1900L	Adjustable Height Trolley with Liner - 1900mm
1801AHT1700L	Adjustable Height Trolley with Liner - 1700mm
1801AHT1500L	Adjustable Height Trolley with Liner - 1500mm



With Adjustable Backrest



Padded Changing Trolley



With Liner

OPTIONAL EXTRAS:

- Low Cotside in 3 different heights 1300mm; **1802ACCWCU13**, 1500mm; **1802ACCWCU15** and 1700mm; **1802ACCWCU17**
- Padding for Cotside; **1802ACCWCP**

FREEWAY ASSISTED OR SELF-PROPELLED SHOWER, TOILET AND COMMUNE CHAIRS

freeway

Manufactured in the UK from only the highest quality materials, the comprehensive range of Freeway assisted or self-propelled shower, toilet and commode chairs are pioneering and firmly established as having a reputation second to none.

High-quality, robust and durable, they are designed to provide the most effective solution to a client's toileting and showering needs.

Freeway Shower Chairs are of modular design, which means they can be configured to a choice of interchangeable seats, backrests, armrests and footrests.

A multitude of accessories are available to meet individual requirements and ensure the most bespoke solution.

All adult models within the range are specifically designed to fit over a standard toilet or commode, with three widths available:

A – 540 mm (21")

B – 490 mm (19")



C – 440 mm (17")

Paediatric, plus size and models to fit over automatic toilets are also included in the range.

The interchangeable makeup of Freeway Shower Chairs makes the range the most adaptable solution on the market as accessories and parts can be added to suit the changing needs of a client.

Chairs can be manufactured with frames to be lowered and raised; rear wheels set back/forward and the frame depth extended and shortened.

Prism Medical UK has a team of experienced Assessors who can carry out a comprehensive shower chair assessment to determine the needs of the individual. To ensure the most effective and supportive solution to meet your client's needs, we highly recommend a shower chair assessment is made prior to purchase.

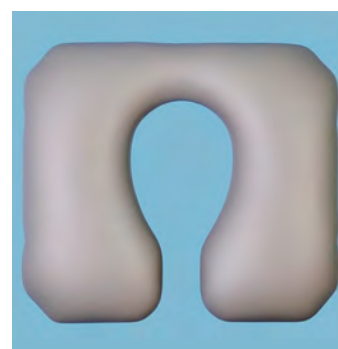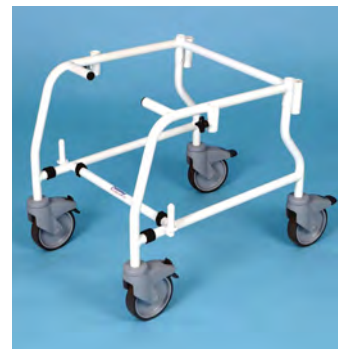

FREEWAY T80 RECLINING SHOWER CHAIR

freeway

The Freeway T80 Shower Chair is highly versatile and the reclined position has significant advantages for increasing comfort and relieving pressure for a wide range of conditions, including:

- Multiple Sclerosis
- Spinal injuries
- Motor Neurone Disease
- Stroke & Other neurological conditions
- Cerebral Palsy

It has two variations; Freeway T80 COM and Freeway T80 AUTO.

Freeway T80 COM: A variation of the T80 reclining shower chair, designed to work seamlessly with the Clos-o-Mat Samoa toilet.

T80 AUTO: A version of the standard T80 which has been modified for use with automatic toilets including the Geberit 8000 range and Clos-o-Mat Palma Vita.

KEY FEATURES & BENEFITS

- This chair reclines through an infinite number of positions up to 45 degrees and is easy to operate.
- The T80 is available in self-propelled versions with a choice of recliners and comes with a wide range of positioning accessories.
- Featuring two high-quality stainless steel gas struts, the tilt action is easy and light to operate.
- Available as assistant-propelled or self-propelled

1 year warranty



Freeway T80 COM



**Freeway T80
AUTO**

Product Information

Product Code	Description
1601T80C	T80 Tilt-in-Space Shower Chair – 440mm wide
1601T80B	T80 Tilt-in-Space Shower Chair – 490mm wide
1601T80A	T80 Tilt-in-Space Shower Chair – 540mm wide

FREEWAY T70 SHOWER CHAIR

freeway

The Freeway T70 Shower Chair is available as a shower, toilet or commode chair. One of the major benefits of the Freeway T70 is that it is fitted with large-braked rear wheels and push rims making this chair extremely easy to self-propel.

As with all Freeway Shower Chairs, only the finest materials have been used during manufacturing to ensure the shower chair is durable and long-lasting.

The Freeway T70 offers the ultimate in hygiene, as the durable seat, backrest and armrests can be simply wiped down after use. The protective coating of the shower chair includes FREECOTE® – our unique antimicrobial and antibacterial coating which is resistant to the hospital superbug MRSA, e-coli, salmonella and listeria.

KEY FEATURES & BENEFITS

- Available as a shower, toilet or commode chair
- Versions available to fit over Clos-o-Mat toilets, automatic bidets and for those who require easy side access
- Perfect for those who require a chair that is extremely easy to self propel
- Puncture proof rear braked wheels with removable push rims
- Large diameter castors to the rear allow easy manoeuvring over uneven surfaces
- Frame constructed from heavy duty, durable tubing in order to provide years of rust free performance
- Pliable backrest moulds to the user while providing excellent support
- Durable seat, backrest and armrests can be easily wiped clean to ensure excellent levels of hygiene
- Coating of the shower chair frame includes FREECOTE®, our unique antibacterial and antimicrobial coating

1 year warranty



200kg/30st
max user weight



Product Information

Product Code	Description
1601T70C	Self propelled Shower/Toilet/Commode Chair – 440mm wide
1601T70B	Self propelled Shower/Toilet/Commode Chair – 490mm wide
1601T70A	Self propelled Shower/Toilet/Commode Chair – 540mm wide

FREEWAY T40 ASSISTANT PROPELLED SHOWER CHAIR

freeway

The Freeway T40 is a highly versatile assistant propelled shower, toilet or commode chair that can be raised/lowered, extended forward, rear castors setback for amputee use.

KEY FEATURES & BENEFITS

- A highly versatile shower chair which can be made bespoke for a variety of applications
- Versions available to fit over Clos-o-Mat toilets and other automatic bidets
- Heavy duty total stop castors and footrests fitted as standard
- Frame constructed from heavy duty, durable tubing in order to provide years of rust free performance
- Durable seat, backrest and armrests can be easily wiped clean to ensure excellent levels of hygiene
- Coating of the shower chair frame includes FREECOTE®, our unique antibacterial and antimicrobial coating.

1 year warranty



200kg/30st
max user weight



Product Information

Product Code	Description
1601T40C	T40 Assistant Propelled Shower Chair – 440mm wide
1601T40B	T40 Assistant Propelled Shower Chair – 490mm wide
1601T40A	T40 Assistant Propelled Shower Chair – 540mm wide

ELEXO SHOWER CHAIR



The Elexo is a height-adjustable shower-toilet chair which provides comfort and security for both carers and patients. Built for safety, comfort and cleanliness, the Elexo is an exceptional tool for carers looking to give residents a secure and calm showering or toileting experience.

KEY FEATURES & BENEFITS

- Simple electric height adjustment
- Hand control with battery indicator
- Sleep mode setting after 60 minutes
- Tilt adjustment to a maximum of 17 degrees
- Easy clean, hygienic seat with toilet opening
- Inlay for the toilet opening
- Foldable upholstered armrests with front lock to enable sideways transfer
- Ergonomically formed seat and backrests
- Fold-out footrests with calf protection
- Fold-away and removable leg supports
- Rail for toilet bucket or bedpan
- Durable castors with swivel & lock features
- Extremely stable due to the fitting of two lifting motors
- Simple to recharge, with charger included

3 year warranty



180kg/28st
max user weight



**Ergonomically
formed seat for
toilet use**

Product Information

Product Code	Seat Width	Seat Depth	Adjustable Height
HLA013	480mm	460mm	50m-105mm

EVE SHOWER AND BATH LIFTER



The sleek Eve Shower and Bath Lifter allows users to be able to safely and comfortably shower or bathe. Highly manoeuvrable, and fitting in most modern care bathtubs, the Eve is suitable for residents who need extra support getting in and out of the bath and is also fully functional as a shower chair.

KEY FEATURES & BENEFITS

- Made entirely of aluminium and stainless steel
- 100% functional as a shower chair
- Simple operation by only one carer
- Unrestricted transfer from the Ascent Pro 200 to the seat lifter
- Ergonomically designed safety steering handles ensure accurate manoeuvring
- Easy-to-maneuver double castors
- Handheld control unit with visual indication of battery level
- Integrated service module reports technical problems with the lift as well as necessary maintenance.

3 year warranty



185kg/29st
max user weight



OPTIONAL EXTRAS:

- Digital scale
- Replacement battery
- External wall charging station

Product Information

Product Code	Length	Width	Height	Lifting Range	Free Space Under Carrier Frame
HLA008	950mm	680mm	1240mm	495mm-1065mm	30mm

7. EVACUATION PRODUCTS

The range of Evacuation Products from Harvest Healthcare are designed to be used in the event of an emergency. Patients will be kept secure and comfortable during an evacuation, as long as the appropriate equipment is used correctly. For more information, please contact a member of the sales team.

EVACUATION SHEET

EVACUATION SLEDGE



EVACUATION SHEET

The Fire Evacuation Sheet is designed to be used in conjunction with standard foam mattresses. During evacuation, the straps are secured around the mattress and the patient is then evacuated to safety. This sheet features a silicone coated side to aid swift and easy movement of the device across all surface types.

The red polyester fabric is Fire Retardant and extremely strong and has been vigorously tested.



Product Information

Product Code	Description
HTE001	Evacuation Sheet

EVACUATION SLEDGE

The Evacuation Sledge fits underneath the existing active mattress and is designed to suit most mattresses.

It contains a high density foam base so that in the event of an evacuation, the air mattress can be deflated by opening the CPR valves and disconnecting from the pump.

The patient is then securely strapped in and evacuated to safety.

The mattress base has a section of low friction material to make it easier to pull on most surfaces.

Easily accessible emergency straps are hidden in pockets to keep them from fouling profiling bed mechanisms.



Product Information

Product Code	Description
HTE008	Evacuation Sledge for Active Mattresses

EVACUATION MATTRESSES



EVACUATION SHEET FOR ACTIVE MATTRESSES

This Evacuation Sheet is designed to take an active air mattress. In the event of an evacuation, the service user is secured to the mattress using the straps on the evacuation sheet and moved to safety along with the mattress. Suitable for both Replacement and Overlay mattresses.

The red polyester fabric is Fire Retardant and extremely strong and has been vigorously tested.



Product Information

Product Code	Description
HTE001.2	Evacuation Sheet for Active Mattresses

WALL MOUNTED EVACUATION MATTRESS WITH STORAGE BAG

This popular product includes a storage bag and is designed to be mounted on a wall until needed. The Wall Mounted Evacuation Mattress includes a 50mm foam base to support the patient during evacuation. During evacuation, the patient is placed on the mattress, secured using the straps provided and evacuated to safety.

The red polyester fabric is Fire Retardant and extremely strong and has been vigorously tested.



Product Information

Product Code	Description
HTE002	Wall Mounted Evacuation Mattress with Storage Bag

STAIRWAY EVACUATION



EMERGENCY STAIRWAY EVACUATION DEVICES

Regulations require that your emergency plan should be appropriate for your premises, and arrangements should be in place for the safe evacuation of people identified as being especially at risk. With these unique chairs, evacuations can be carried out by only one person, giving you a better chance of getting your patients out in time. These chairs can fold down for easy storage and include wall mounting brackets.

EXITMASTER EGO

The Exitmaster Ego is easy to use and is ultra-slim for storage. It has two wheels with auto brakes when folded. It provides a comfortable and supportive seat position and includes a wall mounting bracket. The safe working load is 150kg / 24 Stone.

Product Information

Product Code	Description
HTE202	Exitmaster Ego



EXITMASTER VERSA

The Exitmaster Versa has a 4 wheel design meaning it travels well on landings, corridors and stairs. The seat offers excellent support and comfort and aids easy transfer with an advanced belt guide system. There is also a footrest for comfort. It has safety brakes on the rear wheels and includes a wall mounting bracket. The safe working load is 150kg / 24 Stone.

Product Information

Product Code	Description
HTE200	Exitmaster Versa



EXITMASTER ELITE

The Exitmaster Elite has a 4 wheeled design meaning that it travels well on landings, corridors and stairs. This Exitmaster provides extra comfort with a padded seat and arm rests and an advanced head rest with head band. It also has a footrest. There is an integral 4 point harness and a wall mounting bracket. This stairway evacuation device can double as a transit chair.

The safe working load is 150kg / 24 Stone.

Product Information

Product Code	Description
HTE201	Exitmaster Elite



8. REPOSE SPECIALIST SEATING

In January 2021, Repose Furniture Ltd was acquired by the Prism Group and became a sister company to Harvest Healthcare.

Repose are one of the market leaders in specialist seating solutions, designing and manufacturing a full range of specialist chairs ranging from entry-level riser-recliners to advanced bespoke healthcare chairs.

This acquisition allows Harvest to add the supply and service of specialist seating alongside our current range of beds, pressure care and moving & handling equipment.

All Repose chairs are manufactured in our own factory in the West Midlands. This ensures you receive the highest standards of quality as well as short lead times and bespoke manufacture to meet your needs.

In this section you will find:

Georgia Basic Recliner Chair


Lynton Riser Recliner Chair

Multi C-air Riser Recliner Chair

Multi Flex Chair Porter Chair

Madison Riser Porter Chair

Boston Porter Chair

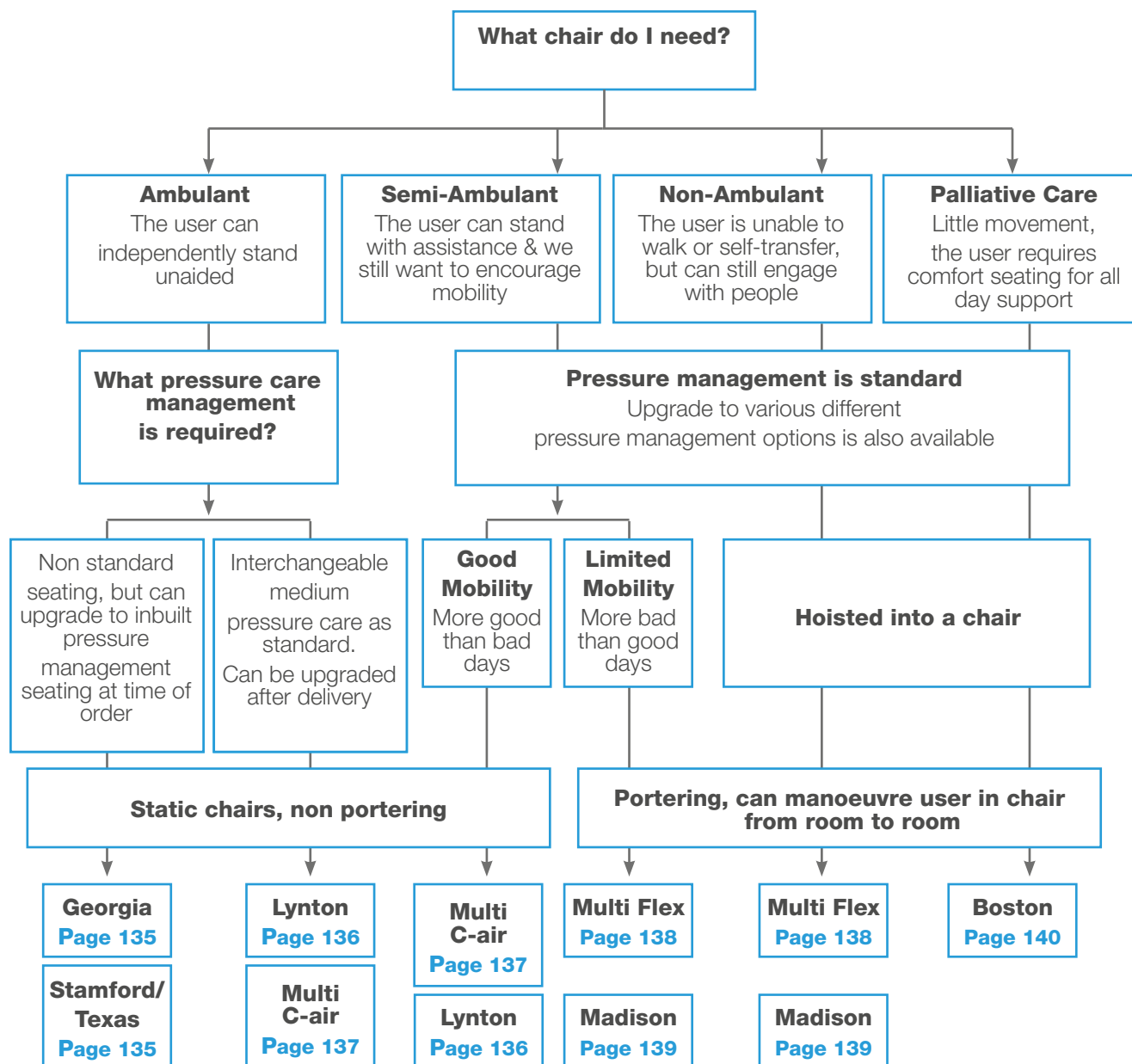

To find out more information about the different seating solutions please visit the Repose website:

reposefurniture.com

GEORGIA BASIC RECLINER CHAIR



Designed for general comfort, this chair has limited seams and simple lines to assist with infection control processes. It is ideal to be used as a bedside seating in care homes, hospital wards, day units and hospices. The Georgia is suitable for ambulant and semi-ambulant users who require small assistance using the rise feature on these chairs. However, it is not suitable for non-ambulant users.

Mechanism Movements:

- Manual recline only – users with good mobility and strength
- Electric recline only – users with reduced mobility but who are able to stand unassisted
- Electric recline and rise – users with reduced mobility who require assistance to stand

The chair is available in Vinyl fabrics as standard.

KEY FEATURES & BENEFITS

- Fixed back cushion style
- Crib 5 Healthcare fabric, easy to clean
- One standard size
- Various mechanism movements
- No pressure care seating
- Limited seams for infection control
- The maximum weight limit can be upgraded to 158kg / 25st

2 year warranty

127kg/20st
weight limit



Stamford Chair
(Wingless Georgia)



Texas Chair
(Georgia with Supportive Cushions)

Product Size Information

Seat Height	480mm (19")
Seat Width	510mm (20")
Seat Depth	510mm (20")
Back Height	710/740mm (28/28")
Arm Height	180mm (7")
Overall Height	1140mm (45")
Overall Width	760mm (30")

Optional Extras

- Magnetic handset storage to side of arm
- Medium to high pressure care seating inclusive of Dartex fabric to seat
- Full length arm cap protectors
- Drop down side table (not available on risers)
- Wall mounted transformer
- USB charger point to outside arm
- Supportive head cushions
- Supportive back cushions

LYNTON RISER RECLINER CHAIR



The Lynton has been specifically designed to be used in domiciliary homes, loan stores, discharge support environments as well as nursing and residential homes. With interchangeable pressure management and posture support, the chair is suitable for ambulant and semi-ambulant users who require a little assistance using the rise feature on these chairs. However, the Lynton is not suitable for non-ambulant users.

Mechanism Movements:

- Dual electric recline and rise – users with reduced mobility who require assistance to stand

Fabric Options:

- Waterproof, Anti-bacterial, stain resistant fabric outer
- Dartex inner

KEY FEATURES & BENEFITS

- Removable / interchangeable back posture support
- Interchangeable low pressure management
- Crib 5 Healthcare fabric, easy to clean
- 6 standard sizes from small to large
- Serviceable chair for easy maintenance in the community
- One person delivery item
- Accupak alarm wireless chair system
- Seat height adjustable
- High pressure care seating inclusive of Dartex fabric to seat (optional extra)
- Wall mounted transformer (optional extra)
- Supportive head cushions (optional extra)
- Fully replaceable parts i.e. wings, arms, seat and backs (optional extra)

2 year warranty

158kg/25st
weight limit



Product Size Information

	Small	Small Back Bracket Set 2	Medium	Medium Back Bracket Set 2	Large	Large Back Bracket Set 2
Seat Height	16"(406mm)	16"(406mm)	18"(457mm)	18"(457mm)	18"(457mm)	20"(508mm)
Seat Width	18"(457mm)	18"(457mm)	20"(508mm)	20"(508mm)	22"(559mm)	22"(559mm)
Seat Depth	20"(508mm)	18"(457mm)	20"(508mm)	18"(457mm)	20"(508mm)	18"(457mm)
Back Height	28"(711mm)	28"(711mm)	28"(711mm)	28"(711mm)	28"(711mm)	28"(711mm)
Arm Height	7"(178mm)	7"(178mm)	7"(178mm)	7"(178mm)	7"(178mm)	7"(178mm)

MULTI C-AIR RISER RECLINER CHAIR



The Multi C-Air has been designed to be used in domiciliary homes, loan stores, discharge support environments as well as nursing and residential homes. With interchangeable pressure management and posture support, the chair is suitable for ambulant and semi-ambulant users who require a little assistance using the rise feature on these chairs. It is not suitable for non-ambulant users.

Mechanism Movements:

- Single tilt in space recline and rise – users with reduced mobility who require assistant to stand
- Dual electric recline and rise – users with reduced mobility who require assistance to stand

Fabric Options:

- Standard fabric available is Vinyl, other fabrics are also available to suit the user's environment
- Dartex as standard

KEY FEATURES & BENEFITS

- Removable / interchangeable back posture support
- Interchangeable low pressure management
- Crib 5 Healthcare fabric, or any design to suit the setting, easy to clean
- 6 standard sizes from small to large
- Seat height adjustable, lockable castor optional
- Dartex seat and leg fabric as standard

2 year warranty

127kg/20st
weight limit



Product Size Information

	Small	Medium	Large
Seat Height	406mm (16")	457mm (18")	508mm (20")
Seat Width	457mm (18")	508mm (20")	559mm (22")
Seat Depth	457mm (18")	508mm (20")	508mm (20")
Seat Height <small>Fitting 50mm Leg adjuster</small>	457mm (18")	508mm (20")	559mm (22")
Depth Height <small>Fitting 50mm Seat Depth Pad</small>	406mm (16")	457mm (18")	457mm (18")
Back Height	711mm / 737mm (28" / 29")		
Arm Height	178mm (7")		

Optional Extras

- High pressure care seating inclusive of Dartex fabric to seat
- Wall mounted transformer
- Full length arm cap protectors
- Supportive head cushions
- Reading lamp
- Accupak alarm wireless chair system
- Various posture support back options

MULTI FLEX PORTER CHAIR



The Multi Flex chair is ideal for nursing and residential home, domiciliary and hospital ward environments. The chair is suitable to reduce manual handling for semi-ambulant and non-ambulant users who require daytime portering in a seated position with interchangeable pressure management and posture support. However, the Multi Flex is not suitable for those who need palliative care.

Mechanism Movements:

- 3 positions, ramped seat angle, to reduce risk of slipping, independent backrest recline

Fabric Options:

- Standard fabric available is Vinyl, other fabrics are also available to suit the user's environment
- Dartex as standard

KEY FEATURES & BENEFITS

- Removable / interchangeable back posture support
- Interchangeable seat low pressure management
- Adjustable seat angle for secure seating 10°-20°
- Manual back recliner
- Crib 5 Healthcare fabric, or any design to suit the setting, easy to clean
- 3 standard sizes from small to large
- Dartex seat and back fabric
- Sliding footplate, push handle and heavy duty lockable castors

2 year warranty

127kg/20st weight limit



Product Size Information

	Small	Medium	Large
Seat Height	410mm (16")	460mm (18")	510mm (20")
Seat Width	410mm (16")	460mm (18")	510mm (20")
Seat Depth	410-533mm (16"-21")	410-533mm (16"-21")	410-533mm (16"-21")
Overall Width	700mm (28")	750mm (29.5")	800mm (31.5")
Overall Depth	810mm (32")	810mm (32")	810mm (32")
Overall Height	960mm (38")	1010mm (40")	1070mm (42")

Optional Extras

- Leg support block
- Connecting straps for leg support block
- Supportive head cushions
- Medium to High pressure care seating inclusive of Dartex fabric to seat

MADISON RISER PORTER CHAIR



The Madison is ideal for nursing and residential home, domiciliary, hospital ward and day unit environments. The chair is suitable to reduce manual handling for semi-ambulant and non-ambulant users who require daytime portering in a seated position with interchangeable pressure management and posture support. It is ideal to encourage client independence when possible. However, the Madison is not suitable for those who need palliative care.

Mechanism Movements:

- Dual electric recline and rise – users with reduced mobility who require assistance to stand

Fabric Options:

- Standard fabric available is Vinyl, other fabrics are also available to suit the user's environment
- Dartex as standard

KEY FEATURES & BENEFITS

- Removable / interchangeable back posture support
- Interchangeable low pressure management
- Electrically operated with rise function
- Crib 5 Healthcare fabric, or any design to suit the setting, easy to clean
- 3 standard sizes from small to large
- Dartex seat and back as standard
- Push handle and heavy duty lockable castors
- Loop for handset storage

2 year warranty

127kg/20st
weight limit



Product Size Information

	Small	Medium
Seat Height	483mm (19")	483mm (19")
Seat Width	470mm (18.5")	521mm (20.5")
Seat Depth	508mm (20")	508mm (20")

Optional Extras

- Sliding footplate
- Supportive head cushions
- Accupak rechargeable wireless battery alarm system
- Removable wireless charging table
- Reading light
- USB charger
- Full length arm protector covers
- Magnet storage handset

BOSTON PORTER CHAIR



The Boston Porter chair is ideal for nursing and residential home, domiciliary, hospital ward, day unit and hospice environments. The chair is suitable to reduce manual handling for non-ambulant users and palliative care patients who require daytime portering in a seated position with interchangeable pressure management and posture support. It is also ideal for users with long term seating requirements as it can help them to have a change of environment and reduce bed use. However, the Boston is not suitable for ambulant and semi-ambulant users.

Mechanism Movements:

- Manual or electric tilt in space, leg rest and back recline.

Fabric Options:

- Standard fabric available is Vinyl, other fabrics are also available to suit the user's environment
- Dartex as standard

KEY FEATURES & BENEFITS

- Removable / interchangeable back posture support
- Interchangeable low pressure management
- Manual independent tilt in space, leg rest and back recliner
- Crib 5 Healthcare fabric, or any design to suit the setting, easy to clean
- D rings for use of optional extras, ie lap strap
- Dartex fabric as standard
- Push handle and heavy duty lockable castors
- Wing push handles
- Rear notes pocket

2 year warranty

127kg/20st
weight limit



Product Size Information

	Small	Medium	Large
Seat Height	635mm (25")	635mm (25")	635mm (25")
Seat Width	457mm (18")	508mm (20")	559mm (22")
Seat Depth	508mm (20")	508mm (20")	508mm (20")
Arm Height	178mm (7")	178mm (7")	178mm (7")
Overall Height	1270mm (50")	1270mm (50")	1270mm (50")
Overall Seat Width	699mm (27.5")	749mm (29.5")	800mm (31.5")
Overall Depth	927mm (36.5")	927mm (36.5")	927mm (36.5")

Optional Extras

- Height/angle adjustable footplate
- Electric option available
- Sliding footplate
- Supportive head cushions
- Full length arm protector covers
- Full length arm protector with width adjustment
- Back support posture control wedges
- Footplate pads

WHAT OUR SERVICE OFFERS YOU



- A tailored service to suit your requirements.
- A skilled, experienced and professional workforce.
- A cost-effective solution.
- Fully traceable work program and database with comprehensive reporting.
- A personal service with named contacts who will look after your needs.
- A complete parts service.

PLANNED/PREVENTATIVE MAINTENANCE

- Lifting Equipment – examination and servicing (LOLER).
- Beds – Inspection and service (PUWER).
- Mattresses – Inspection, decontamination and service/repairs.
- Wheelchairs – Inspection and service (PUWER).
- Suction machines/Nebulisers – Inspection and Service.
- Training packages for your maintenance staff.
- Production and maintenance of your asset register.

INSTALLATION & TRAINING

We have a team of highly experienced installation engineers who are available to install and demonstrate our products for you. We also provide moving and handling advice and training for you to be confident in using our products both efficiently and safely.

MAINTENANCE & REPAIRS

Our focus is to provide you with a first-time, quick and efficient fix, whenever possible. As a manufacturer of healthcare equipment, we have a comprehensively stocked parts department.



REACTIVE MAINTENANCE (ALL PRODUCTS)

- On or off-site repair options.
- Same day 24/7 response time.
- Mattress management – repairs, decontamination, exchange.
- Training packages for HH products.
- Telephone support for troubleshooting and spares.

Whether you want us to support you when you require, or to take on a complete asset management role, give us a call and let us know how you would like us to tailor a package for you.

MATTRESS DECONTAMINATION

Our air mattress repairs and decontamination division operates an experienced and dedicated service in repair, maintenance and cleansing of most makes of active pressure relief systems. Our on-site laundry uses OTEX technology which injects OZONE into the wash cycle destroying all bacteria assuring that your equipment is certified decontaminated with full traceability.

EQUIPMENT MANAGEMENT

From hospitals to walk-in centres and rehabilitation clinics, we have a dedicated department that specialises in the maintenance and servicing of equipment.



www.harvesthealthcare.co.uk

Sheaf House, Bradmarsh Business Park, Bradmarsh Way, Rotherham S60 1BW
T +44 (0)1709 377172 **F** +44 (0)1709 377173 **E** sales@harvesthealthcare.co.uk



The mattress management programme is a simple, easy and cost-effective way to manage your active mattresses.


Just ask for an active mattress when you need one and receive a fully decontaminated and serviced mattress the very next day*.

We will come and collect your old mattress, ready to be cleaned, serviced and stored for the next time it is needed.

All our decontamination is done to Public Health England standards with Coronavirus certified cleaner, and each mattress is given its own decontamination certificate.

This ensures safe transfer of mattresses between residents, and peace of mind for you.


Never store, wash, or repair active mattresses again, and keep your residents safe and well, with the Harvest Healthcare mattress management programme.

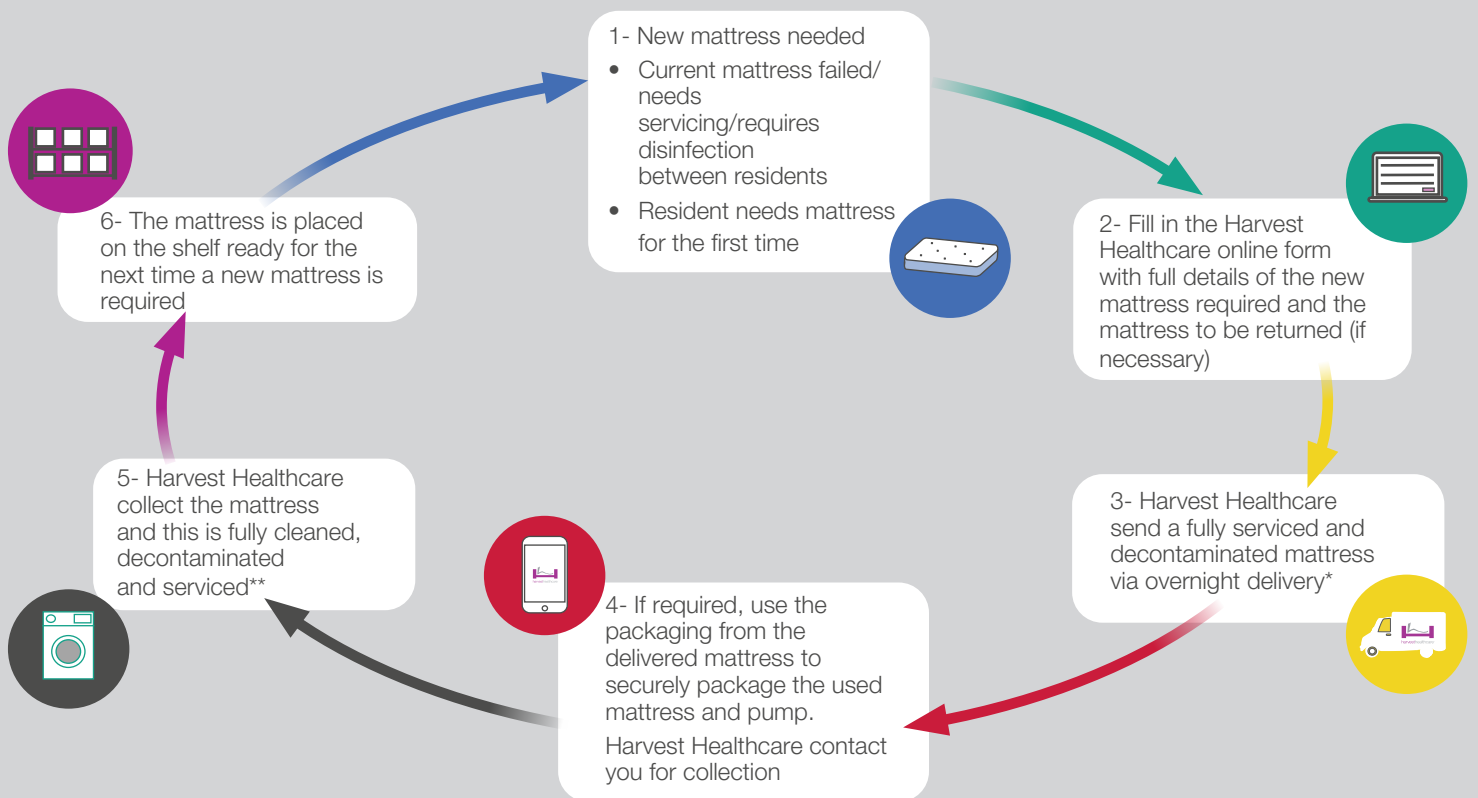

*Next day / overnight delivery available if ordered before 12pm.

**Decontaminated to HTM01-04 and serviced according to ISO13485. Mattress supplied will be similar specification to mattress collected (if applicable).

prismhealthcare

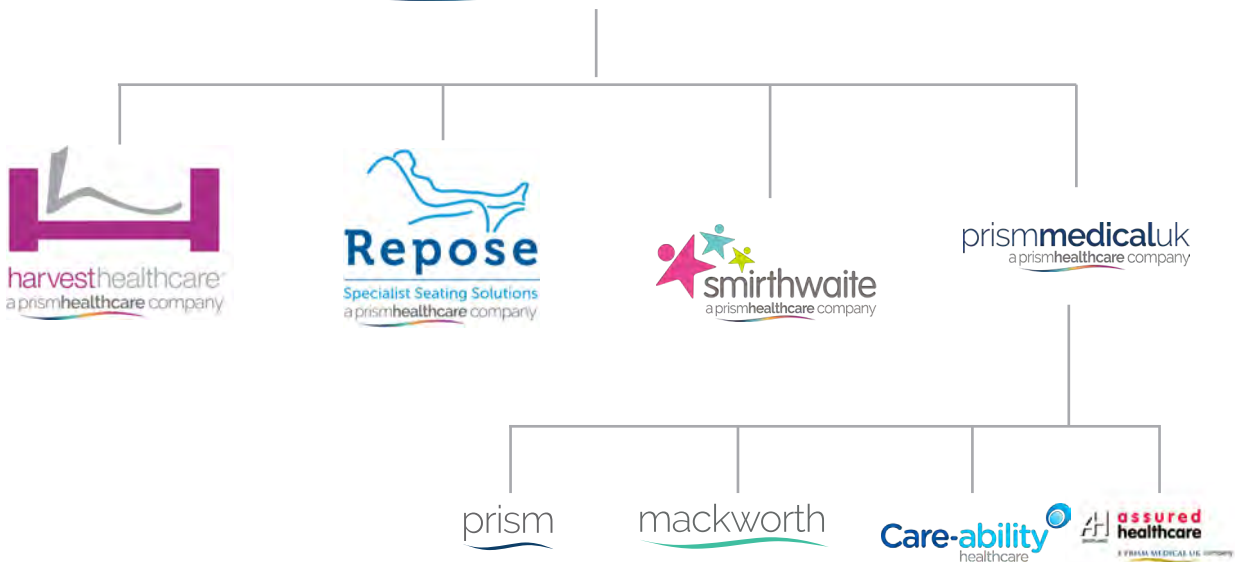

www.harvesthealthcare.co.uk

Email: sales@harvesthealthcare.co.uk

Tel: 01709 377172

Harvest Healthcare is a leading manufacturer and supplier of quality healthcare equipment. Specialised in manufacturing and distributing pressure care equipment, they continually develop their product range to deliver complete solutions for all modern care environments.



www.reposefurniture.com

Email: info@reposefurniture.com

Tel: 01384 567401

For over 20 years, Repose has designed and manufactured bespoke healthcare and riser recliner chairs for all care environments.



www.smirthwaite.co.uk

Email: info@smirthwaite.co.uk

Tel: 01626 835552

Smirthwaite have been designing, manufacturing and selling equipment for children and young people with special needs since 1986.



www.prismmedical.co.uk

Email: info@prismmedical.co.uk

Tel: 01924 840100

Prism Medical UK provides moving, handling and bathing products for a variety of client needs, with cost-effective solutions for tight budgets available from Mackworth, and premium options available through Freeway as required.



www.care-ability.co.uk

Email: info@care-ability.co.uk

Tel: 01914 159777

Care-ability Healthcare is a leading independent supplier and manufacturer of moving and handling products for the healthcare industry.



www.assured-hs.co.uk

Email: info@assured-hs.co.uk

Tel: 01786 478390

Assured Healthcare (Scotland) Ltd is an established and trusted supplier to the healthcare sector with over 20 years experience in the installation of patient lifting and bathing equipment.



harvesthealthcare®
a prismhealthcare company

Sheaf House, Bradmarsh Business Park,
Bradmarsh Way, Rotherham S60 1BW

T +44 (0)1709 377172
F +44 (0)1709 377173
E sales@harvesthealthcare.co.uk

www.harvesthealthcare.co.uk

Harvest Healthcare Brochure V3 - August 2023

© Copyright Harvest Healthcare

Harvest Healthcare Limited. Company No: 07210261

Our Full Terms & Conditions are available by request or can be found on our Website.

Harvest Healthcare reserves the right to alter or amend this document without notice.



harvest-healthcare-ltd



@HarvestHealthUK



Harvest Healthcare