


Visco MultiZone

Our new static mattress designed with enhanced features to provide exceptional pressure relief, comfort, and value.

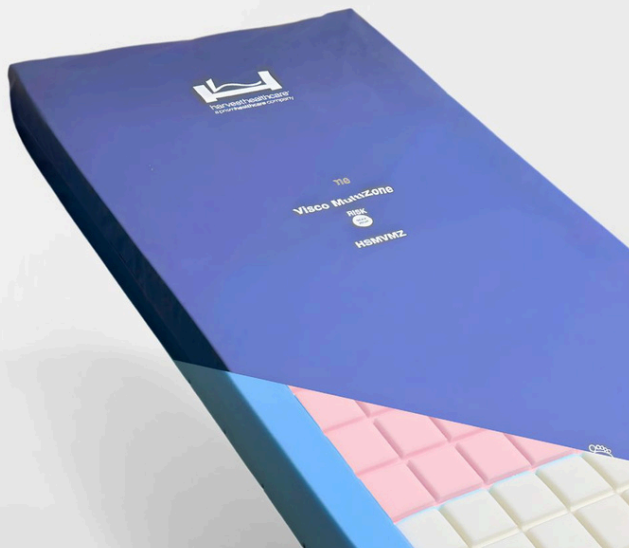

VISCO MULTIZONE

PRODUCT CODE: HSMVMZ

The Visco MultiZone Mattress is our enhanced static mattress.

Positioning specialist foam “Zones” across key areas of the mattress allows for more pressure relief around the most vulnerable areas of the body. The varying firmness in foam at zonal areas contributes to greater immersion and envelopment, resulting in greater pressure redistribution, as well as providing additional support in the areas where it is needed most.

HIGHLIGHTS

- **Reduces the need for dynamic mattresses.**
- **Very High-Risk**
- **Made in the UK, assembled by Harvest.**
- **Lighter weight compared to previous Viscoflex mattresses.**
- **Complete corners that reduce the risk of entrapment.**

CONTACT US

☎ +44 (0)1709 377172
✉ sales@harvesthealthcare.co.uk

FOLLOW US

🌐 [harvest-healthcare-ltd](https://www.harvest-healthcare-ltd.co.uk)
✕ @HarvestHealthUK

OUR LOCATION

📍 Sheaf House, Bradmarsh Business Park,
Bradmarsh Way, Rotherham S60 1BW

KEY PRODUCT FEATURES



Pressure Relief

Strategically placed foam zones and memory foam visco-elastic provide superior pressure redistribution.



Comfort and Stability

The high-resilience foam provides excellent support and stability, greater immersion and envelopment increases comfort.



Easy to Maintain

The mattress is supplied with a four-way stretch, fully welded, PU coated, polyester cover that is water-resistant and easy to clean.



High Maximum User Weight

The mattress has a maximum user weight of 177kg, making it suitable for a wide range of individuals.



Complete Corners

Complete corners reduce the risk of entrapment.

MORE INFO

• 3 Year Warranty


For More Information

Scan the QR Code below to visit the website.



TECHNICAL DATA

Product Code	HSMVMZ
Warranty	3 Years
Mattress Dimensions (mm)	2000mm x 900mm x 150mm
Max. User Weight	27 Stone / 177kg
Min. User Weight	N/A
Pressure Risk	Very High

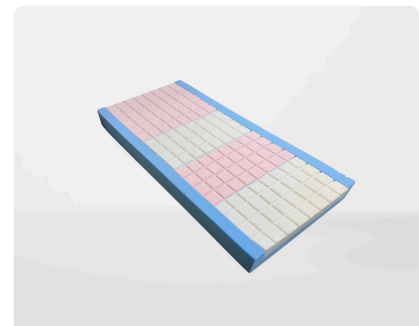

Important Safety Notice

Harvest Healthcare Ltd. supplies equipment and appliances to cover a range of medical conditions. Specific choice of equipment is the sole responsibility of the clinical practitioner or carer overseeing individual patient care. Harvest Healthcare Ltd supplies a range of products such as foam mattresses, active air mattresses, fall protection mats and bed rail protection bumpers. These products are usually covered in vapour-permeable PU fabric that could cause a risk of suffocation in the event of prolonged facial contact. Mattress risk categories are to be used as a guide only. Risk assessments specific to the user must be carried out prior to the use of all products. Products and data continually undergo development and may therefore deviate from the details given.

Version 2: Nov 2025